

Culture Agar Preparation

Overview

Prepare and make agar plates to cultivate bacteria.

A. Materials


- a. Agar power
- b. ddH₂O
- c. glycerol
- d. 500mL flask
- e. Magnetic Stirrer
- f. Antibiotics
 - i. Ampicillin
 - ii. Chloramphenicol
 - iii. Erythromycin

B. Procedure

1. Add 250mL ddH₂O, 9.25 LB agar powder, and 2mL glycerol to a 500mL flask
2. Add a stirring bar
3. Stir the mixture until all powder are dissolved
4. Cover the flask with aluminum foil and an autoclave indicator stucked on top
5. Autoclave for 30min
6. Put it in BSC while wait for it to cool down (can also stir before putting in BSC)
7. Add 250μL of antibiotics and gently mix it
8. Distribute 250μL of medium in plate

Culture Agar Preparation

Ingredients:

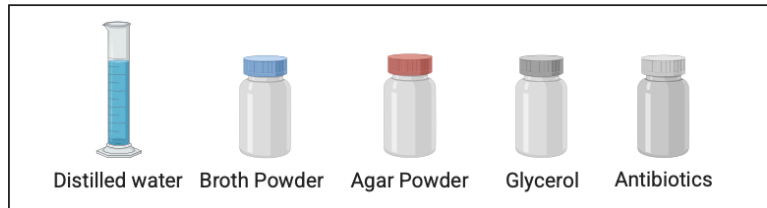


- ① Add 9.25g Broth Powder, 3.75g of agar powder, 2ml glycerol into the mixture




- ② Add a magnetic stirring bar and stir on the plate until the medium is dissolved completely



- ③ Transfer to flask and autoclave

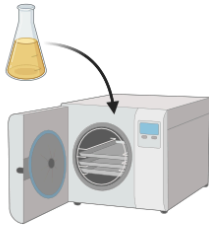

- ④ Add antibiotic once the solution reaches ~60°C

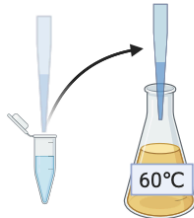

- ⑤ Pour medium in each petri dishes and let them solidify (operate in the BSC)

