

reference format

Font format: Arial, Front size:11

Date:

Title:

Purpose:

Procedure:

1. Heat the top of the glass vial with the Bunsen burner.
2. Add sterile water and crack it with tweezers.
3. Remove paper and cotton.
4. Add 500 μ L LB medium and resuspend the tablet.
5. Drop the supernatant and streak it with four-quadrant streak method.
6. Spread 100 μ L supernatant on LB plate with 100 μ L 25 mg/mL chloramphenicol.
7. Inoculate the remaining supernatant in 5 mL LB medium with 25 μ g/mL chloramphenicol.
8. Incubate in 37°C, started from 18:00

Done by Chars, recorded by Chars.

Results:

Common mistakes.

1. chloramphenicol
2. pNW33N-L2 (sfGFP)
3. Autoclave for 20 minutes at 127°C.
4. *Bacillus subtilis*

Make format consistent.

1. degrees Celsius->°C
2. ul-> μ L
3. 23444xg->23444 g
4. ddH₂O->ddH₂O
- 5.

July

2026/7/1 

Date: 20260701

Title: Colony PCR

Purpose: To verify bacterial colonies after a transformation step in molecular cloning

Procedure:

1. Prepare the following reaction (master mix) *21 and dispense solution in PCR tubes

	1 reaction	21 reactions
Forward Primer (T7 promoter)	1 μ L	21 μ L
Reverse Primer (SP6 promoter)	1 μ L	21 μ L
Qamp master mix	10 μ L	210 μ L
ddH ₂ O	8 μ L	168 μ L

2. Circle 3 colonies each from plates A to G

A: pGEM-Teasy-Csn: LEA
B: pGEM-Teasy-Csn: AmilCP
C: pGEM-Teasy-Csn: Accd
D: pGEM-Teasy-ccasm3
E: pGEM-Teasy-ccaR
F: pGEM-Teasy-Ho1: PcyA
G: pGEM-Teasy-pcpcG2-pveg
F

3. Pick up 3 colonies each agar plate with tip and streak the labelled grid
4. Dispense the tips with colonies into the PCR strips and wait for 5 mins
5. Take out the tips and close the lid
6. Place into PCR machine for the thermal cycle

Done by Naomi, Sarah

Recorded by Sophie H.

Date: 20260701

Title: 1.0% Agarese Gel preparation

Purpose: To prepare 1.0% Agarese Gel for further electrophoresis experiment

Procedure:

1. Prepare and assemble the mold and comb
2. Weight 0.4g of agarose powder and add 40ml of 1X TAE solution into flask
3. Cover the flask with plastic wrap, then use a tip to puncture holes on the plastic wrap
4. Microwave until the gel stops smoking but is still hot
5. Add 2 μ L of safety dye into the gel and agitate it
6. Pour the melted gel into the tray
7. Wait until the gel cools down
8. Repeat steps 1~8 to create second copy of gel

Done by Sophie L., Ethan L.

Supervised by Chars

Recorded by Sophie H.

Date: 20260701

Title: Gel Electrophoresis

Purpose: To test the concentration of respective solutions

Procedure:


1. Load the gel with the following layout(5 μ L each):

M: skb DNA ladder	M A1 A2 A3 B1 B2 B3 C1 C2 C3
A: pGEM-Teasy-Csn: LEA	D1 D2 D3 E1 E2 E3
B: pGEM-Teasy-Csn: AmilCP	
C: pGEM-Teasy-Csn: Accd	
D: pGEM-Teasy-ccasm3	
E: pGEM-Teasy-ccaR	
F: pGEM-Teasy-Ho1: PcyA	
G: pGEM-Teasy-pcpcG2-pv	

eg

2. Resolve at 100v for 30 min

Result:

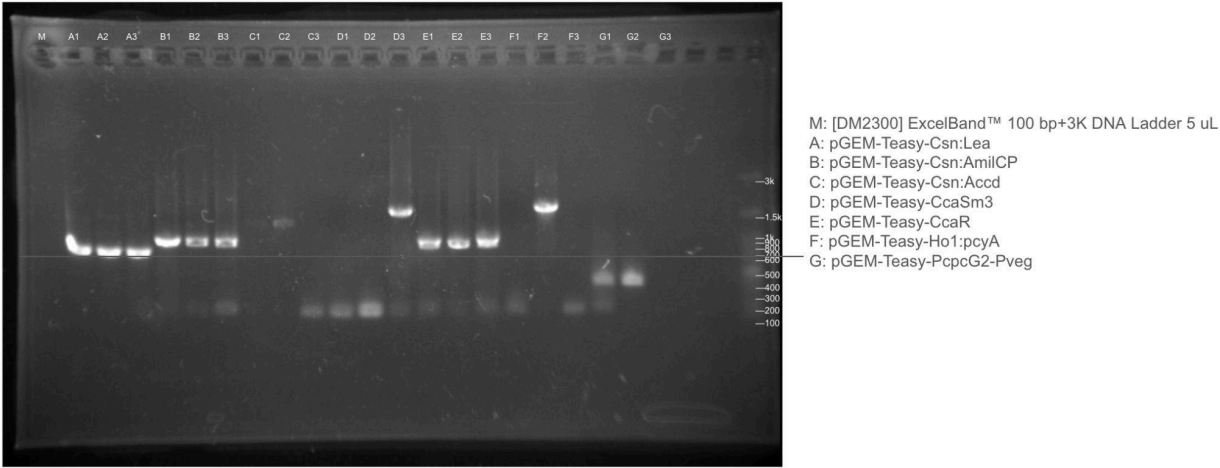

Done by Sophie L., Ethan L.
Supervised by Chars
Recorded by Sophie H.

Date: 20260701
Title: Colony PCR
Purpose: To verify bacterial colonies after a transformation step in molecular cloning
Procedure:

1. Prepare the following reaction and dispense them into PCR tubes:

	B (4rxn)	C(6rxn)	E(6rxn)	F(6rxn)
Forward Primer (T7 promoter)	5µL	7µL	7µL	7µL
Reverse Primer (SP6 promoter)	5µL	7µL	7µL	7µL
Qamp master mix	50µL	70µL	70µL	70µL
ddH ₂ O	40µL	56µL	56µL	56µL

2. Circle colonies without fluorescence form B.C.E>F

B: pSTK-1D-Pveg-EEBD-RBS0-Csn::ACCD-B0015
 C: pSTK-1C-Pveg-RBS0-CcaR-B0015
 E: pSTK-Exp2-Spacer x3-Pveg-EEBD-RBS0-Csn::amilCp-B0015
 F: pSTK-Exp1-Pveg-EEBD-RBS0-Spacer-GanA::Gaa-histag-B0015

3. Pick up the 22 colonies with a tip and streak the labelled grid
4. Dispense the tips with colonies into PCR strip and wait 5 min
5. Take out the tips and spindown
6. Place into PCR machine and run Q-amp thermal cycle

Done by: Chloe, Sophie L, Ethan

Supervised by: Chars

Date: 20260701

Title: Gel Electrophoresis

Purpose: verify the size of colony PCR products

Procedure:

1. Load the samples in the following layout:

<u>M</u>	<u>B1</u>	<u>B2</u>	<u>B3</u>	<u>B4</u>	<u>C1</u>	<u>C2</u>	<u>C3</u>	<u>C4</u>	<u>C5</u>	<u>C6</u>	<u>E1</u>	<u>E2</u>	<u>E3</u>	<u>E4</u>	<u>E5</u>	<u>E6</u>	<u>F1</u>	<u>F2</u>	<u>F3</u>	<u>F4</u>	<u>F5</u>	<u>F6</u>
5	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6

M: 3kb DNA ladder

B: pSTK-1C-Pveg-EEBD-RBS0-Csn:AmilCP-B0015


C: pSTK-Exp2-Pveg-RBS0-CcaR-B0015

E: pSTK-Exp2-Pveg-spacer x3 - Pveg-EEBD-RBS0-Csn:AmilCP-B0015

F: pSTK-Exp1-Pveg-EEBD-RBS0-spacer GanA-Gala-histag-B0015

2. Resolve at 100V for 30minutes

Results:

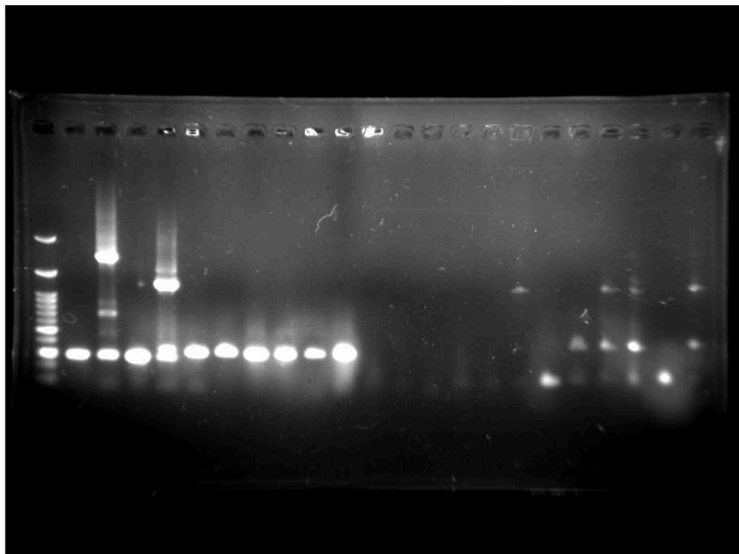

M: [DM2300] ExcelBand™ 100 bp+3K DNA Ladder 5 uL
 B: pSTK-EXP2-Pveg-EEBD-Csn:AmilCP-B0015
 C: pSTK-1C-Pveg-RBS0-CcaR-B0015
 E:pSTK-EXP2-Pveg-Spacerx3-EEBD-Csn:AmilCP-B0015
 F:pSTK-EXP1-Pveg-EEBD-RBS0-spacer-GanA-Gala-histag-B0015

Done by Sophie H.

Date: 20260701

Title: Transformation of TA cloning

Purpose: Transform and cultivate engineered bacteria on agar plate

Procedure:

1. Turn on water bath and set it to 42°C and prepare an ice bucket
2. Mix 100µL of competent cell with 10µL plasmid
3. Incubate on ice for 20min
4. Prepare reaction for blue-white screening → 1000µL of ddH₂O and 4µL of x-gal & 1µL of 1M IPTF, then spread on LBA plate
5. Heat shock at 42°C for 45 sec
6. Incubate on ice for 3 min
7. Add 1ml of SOC median to each microtube and incubate at 37°C for 1 hour
8. Centrifuge at 5000 g for 5 min and remove 1 ml of supernatant
9. Resuspend the cell pallet in the remaining supernatant & spread onto LBA plate

A: Csn:LEA

B: Csn:ACCD

C: PcpcG2-veg

D: Cpp:AmilCP

E: Ho1:pcyA

F: CcaR

G: Csn:AmilCP

H: Ccasm3

Done by Sophie L, Chloe, Ethan L

Supervised by Chars

2026/7/2



Date: 20260702

Title: Chlorophyll Extraction

Purpose: To isolate the green photosynthetic pigment from arabidopsis

Procedure:

1. Weight the fresh weight of arabidopsis (9 at a time)
2. Take the plant from the agar plate, careful not to hurt it by squeezing the tweezers.
3. Put the 9 plants on a tissue paper, lightly press to remove the agar striking to the plates
4. Pick up the plants and put them on a weighting plate to weight
5. Weight the plants
6. Pick up the plants and put them in a microtube, take care that the roots need to go in too
7. Repeat these steps for as many times as the sets of arabidopsis you are testing
8. Place 3 glass bits in each of the microtubes, place in -80°C for 10 min
9. Dilute acetone to 80%, add 1 mL to each tube (10mL in total)
10. Wait for the chlorophyll to extract, and can use a vortex to crush the plants
11. Add the liquid to the circus well, 200 µL in each well, 3 wells for each sample
12. ELISA reader insert the 96 well, run the result through the equation below:

*Total chlorophyll content equation:

$$(20.2 \cdot OD_{645} + 0.02 \cdot OD_{663}) / (1000 \cdot fw)$$

Done by Phoebe, Abby T.

Recorded by Sophie H.

Date: 20260702

Title: Colony PCR

Purpose: To verify bacterial colonies after a transformation step in molecular cloning

Procedure:

1. Prepare the following reactions and dispense them into PCR tubes

	Forward Primer	Reverse Primer	Q-Amp master mix	dd H ₂ O
C (4 rxn)	T7 promoter 5µL	Sp6 promoter 5µL	50µL	40µL
D (4 rxn)	T7 promoter 5µL	Sp6 promoter 5µL	50µL	40µL
F1 (4 rxn)	T7 promoter 5µL	Sp6 promoter 5µL	50µL	40µL
G (4 rxn)	T7 promoter 5µL	Sp6 promoter 5µL	50µL	40µL
E (6 rxn)	pSTK-Exp fw 7µL	pSTK-Exp rv: 7µL	70µL	56µL

F2-1 (6 rxn)	pSTK-Exp fw: 7µL	pSTK-Exp rv: 7µL	70µL	56µL
F2-2 (11 rxn)	pSTK-Exp fw: 7µL	pSTK-Exp rv: 12µL	100µL	96µL

2. Circle colonies without fluorescence from C, D, F, G, E, F2-1, F2-2

C: pGEM-Teasy-Csn: Accd

D: pGEM-Teasy-ccasm3

F1: pGEM-Teasy-Ho1:PcyA

G: pGEM-Teasy-pcpcG2-pveg

E: spacer*3-pveg-EEBD-RBS0-Csn:AmilCP-B0015 pExp 2

F2-1: 20260701pveg-EEBD-RBS0-spacer-Gana-Gala-histag-B0015 pExp1

F2-2: 20260630pveg-EEBD-RBS0-spacer-Gana-Gala-histag-B0015 pExp1

3. Pick up the colonies with a tip and streak the labeled grid
4. Dispense the tip with colonies into PCR strip and wait for 5 min
5. Take out the tips and spindown
6. Place into PCR machine and run Q-amp thermal cycle

Done by Chloe, Alex, Ryan Wei

Supervised by Ryan Yuan, Olivia

Record by Sophie H.

Date: 20260702

Title: Liquid Culture

Purpose: To grow bacteria inside liquid culture from agar plate

Procedure:

1. Label 9 liquid culture tubes
2. Turn on bunsen burner
3. Add 5µL of TB medium and 5µL of ampicillin into each tube
4. Dip a pipette tip into the colonies
5. Dispense tips in their respective tubes
6. Place the tubes in the incubator at a slanted angle

Done by Alex

Supervised by Chars

Record by Sophie H.

2026/7/3



Date: 20260703

Title: Prepare 0.8% Agarose Gel

Purpose: Make gel for electrophoresis

Procedure:

1. Prepare and assemble mold and comb
2. Weigh 0.16 agarose powder and add 80 mL of 1X TAE solution into flask
3. Cover the flask with plastic wrap, then use a tip to puncture holes on the plastic wrap.
4. Microwave until all particles are dissolved
5. Wait until the melted gel stops smoking but is still hot.
6. Add 4 μL of safety dye and agitate it
7. Pour the gel two trays with equal amount
8. Wait until the gel cools down

Done by Alex

Supervised by Chars

Date: 20260703

Title: Gel Electrophoresis

Purpose: To test the concentration of respective solutions


Procedure:

1. Load the gel with the following layout(5 μL each):

M C1 C2 C3 C4 D1 D2 D3 D4 F1(1) F1(2) F1(3) F1(4) G1 G2 G3 G4
M E1 E2 E3 E4 E5 E6 F2-1(1) F2-1(2) F2-1(3) F2-1(4) F2-1(5) F2-1(6)
F2-2(1) F2-2(2) F2-2(3) F2-2(4) F2-2(5) F2-2(6) F2-2(7) F2-2(8) F2-2(9) F2-2(10) F2-2(11)

2. Resolve at 100V for 30 min

Result:

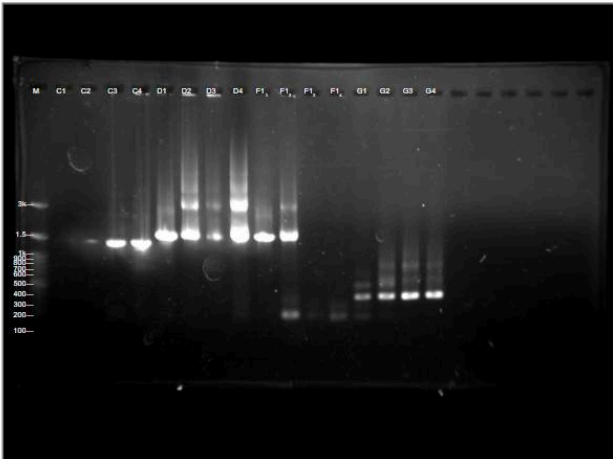

M: [DM2300] ExcelBand™ 100 bp+3K DNA Ladder 5 uL
C: pGEM-Teasy-Csn:Accd
D: pGEM-Teasy-CcaSm3
F1: pGEM-Teasy-Ho1:PcyA
G: pGEM-Teasy-PcpcG2-veg



M: [DM2300] ExcelBand™ 100 bp+3K DNA Ladder 5 uL
E: spacer * 3 - pveg-EEBD-RBS0-Csn:AmiCP-B0015 pExp2
F2-1: 20260701 pveg-EEBD-RBS0-spacer-Gana-Gala-histag-B0015 pExp1
F2-2: 20260630 pveg-EEBD-RBS0-spacer-Gana-Gala-histag-B0015 pExp1

Done by Alex
Supervised by Chars

Date: 20260703

Title: Miniprep

Purpose: Plasmid purification

Procedure:

A1: pGEM-Teasy-Csn:LEA

B1: pGEM-Teasy-Csn:AmilCP

E1: pGEM-Teasy-CcaR

A2: pGEM-Teasy-Csn:LEA

B2: pGEM-Teasy-Csn:AmilCP

E2: pGEM-Teasy-CcaR

A3: pGEM-Teasy-Csn:LEA

B3: pGEM-Teasy-Csn:AmilCP

E3: pGEM-Teasy-CcaR

1. Transfer 1.8 mL of medium to 2 mL microtube
2. Centrifuge at 13444 g 30 seconds
3. Discard the supernatant
4. Add 200µL of MX1 buffer into each microtube
5. Resuspend the pellet through pipetting
6. Add 250 µL of MX2 buffer
7. Gently mix by inverting microtube 8-10 times
8. Incubate at room temp for 2 mins
9. Add 350µL of MX3 buffer
10. Gently mix by inverting microtube 4-6 times
11. centrifuge for 5 minutes at 16000 g
12. Transfer supernatant (~800µL) into assembled column
13. Centrifuge for 30 seconds at 7000 g and discard flow through
14. Wash the column with 500µL WN buffer and centrifuge for 30 seconds at 7000 g and discard the ft
15. Wash the column with 700µL WS buffer and centrifuge for 30 seconds at 7000 g and discard the ft
16. Place the column into new 1.5 mL microtube
17. Add 50µL of ddH₂O to the center of membrane and let it stand for 2-3 min and centrifuge for 3 min at 10,000 g

Result:

Biodrop Results	Concentration ng/µL	A260/230	A260/280
A1	77	1.167	1.925
A2	60	1.500	1.935
A3	65	1.083	1.857
B1	55	1.375	1.897
B2	67	1.241	1.861
B3	49	1.195	1.885
E1	71	1.511	1.972
E2	77	1.481	2.026
E3	66	1.245	1.886

Done by: Abigail and Chloe

Date: 20260703

Title: Liquid culture

Purpose: Let bacteria grow in liquid medium

Procedure:

1. Label 7 liquid culture tubes
2. Turn on bunsen burner
3. Add 5mL of TB medium and 5 μ L of ampicillin into 6 tubes
4. Add 5mL of TB medium and 5 μ L of chloramphenicol into 1 tubes (E3)
5. Dip a pipette tip into the colonies
6. Dispense tips into their respective tubes
7. Place tubes in the incubator at a slanted angle

Done by: Alex & Sophie C

Supervised by: Chars

2026/7/4



Date: 07/04

Title: Miniprep

Purpose: extract plasmid DNA from liquid culture of colony PCR products

Procedure:

0701 E3: pSTK-Exp2- spacer ×3 - pveg - EEBD - RBS0 - Csn:AmilCP - B0015

0701 F2: pGEM -Teasy-Ho1:pcyA

0701 F4: pSTK-Exp1-Pveg-EEBD-RBS0-spacer-GanA-Gala-histag-B0015

0702 C3: pGEM-Teasy-Csn:ACCD

0702 D1: pGEM-Teasy-Ccasm3

0702 F1: pGEM-Teasy-Ho1:pcyA

0702 G3: pGEM-Teasy-PcpcG2-veg

1. Transfer 1.8 µL TB medium to 2 mL microtube.
2. Centrifuge at 16,000 g for 1 min.
3. Discard the supernatant.
4. Repeat steps 2-4 until collected all medium.
5. Add 200 µL MX1 buffer.
6. Resuspend the pellet by vortexing or pipetting.
7. The cells should completely resuspend.
8. Add 250 µL MX2 buffer.
9. Gently mix the buffer well.
10. Do not vortex. Invert microtube 8-10 times.
11. Incubate the mixture at room temperature for 2 min.
12. Add 350 µL MX3 buffer.
13. Immediately and gently mix the solution. (Do not vortex.)
14. Centrifuge at 16,000 g for 5 min.
15. While centrifuging, assemble the column and collection tube.
16. Transfer the supernatant, around 800 µL, to the column.
17. Centrifuge at 7,000 g for 30 sec.
18. Discard the flow through.
19. Add 500 µL WN buffer.
20. Centrifuge at 7,000 g for 30 sec.
21. Discard the flow through.
22. Add 500 µL WS buffer.
23. Centrifuge at 7,000 g for 30 sec.
24. Discard the flow-through.
25. Centrifuge the empty column at 10,000 g for 3 min to remove residual ethanol.
26. Place the column into a new 1.5 mL microtube.
27. Add 50 µL sterile distilled water on the center of the membrane.
28. Incubate the column at room temperature for 3 min.
29. Centrifuge at 10,000 g for 3 min.

	Concentration	A260/230	A260/280
E3	130	0.802	1.831
F2	77	0.987	1.878
F4	93	0.959	1.860
C3	142	0/940	1.844
D1	128	0.921	1.910
F1	136	0.845	1.838
G3	129	0.816	1.817

Done by Sarah, Sophie H

Date: 07/04

Title: Liquid Culture

Purpose: Bacteria culture

Procedure:

1. Label 4 liquid culture tubes
2. Turn on bunsen burner
3. Add 5mL of TB medium and 5µL of ampicillin into each tubes
4. Dip a pipette tip into the colonies
5. Dispense tips into their respective tubes
6. Place tubes in the incubator at a slanted angle

Done by Chloe, Jacquelyn

Supervised by Alex

Date: 07/04

Title: Making salinity agar plates

Purpose: prepare salinity agar plates for plant model Experimental Set #5 0,1,5,10 mM trehalose vs. 150, 250mM NaCl

Procedure:

The following protocol is for 600ml of agar

- 12 agar plates
 - control(no trehalose) x control(no NaCl)
 - control(no trehalose) x 150mM NaCl
 - control(no trehalose) x 250mM NaCl
 - 1mM trehalose x control(no NaCl)
 - 1mM trehalose x 150mM NaCl
 - 1mM trehalose x 250mM NaCl
 - 5mM trehalose x control(no NaCl)
 - 5mM trehalose x 150mM NaCl

- 5mM trehalose x 250mM NaCl
- 10mM trehalose x control(no NaCl)
- 10mM trehalose x 150mM NaCl
- 10mM trehalose x 250mM NaCl

1. Prepare the following:

½ MS medium powder(w/o sucrose)	1.38g
MES buffer	0.3g
Distilled water	600ml
Sucrose(1%)	6g
Agar powder(1%)	6g

2. Add the contents into a 1000ml flask
3. Use a hot plate to stir (with a magnetic stirring rod) until no sediments are left in the solution
4. Use the pH probe and measure the pH of the solution, make sure the pH is 5.7-5.8
 - a. If the pH is too low, add 0.1M KOH slowly to increase it
5. Pour 200ml of mixture out into another flask to separate into three flasks, label with [0mM], [150mM] and [250mM]
6. Add 2.192 g of NaCl into the flask labelled [150mM NaCl]
7. Add 3.652 g of NaCl into the flask labelled [250mM NaCl]
8. Prepare 4 flasks for different trehalose concentrations(0, 1, 5, 10mM)
9. Prepare 12 flasks and pour the respective trehalose and NaCl concentration in each of the flasks
 - a. making sure there is 50ml in each of the flasks(12)
 - b. label each flask clearly according to their trehalose concentrations
10. To their respective trehalose concentrations, add

Trehalose concentration	Amount of trehalose (for 50ml of medium)
0mM	0µL
1mM	500µL
5mM	2500µL
10mM	5000µL

11. Wrap the necks of the flasks with aluminum foil and stick a piece of indicator tape on top
12. Autoclave for 20 minutes at 127°C
 - a. Wait for the pressure to reach 0 before retrieving the flask
13. During the autoclave, label 8 of the petri dishes with a marker in the biosafety cabinet
 - a. right top: date and arrow pointing up to show direction for pictures
 - b. bottom lid: different trehalose concentrations and NaCl concentration or control

14. Transfer half of the flasks to BSC

- a. Store the other half of the flasks in the hot bath while pouring so that it does not solidify

15. Make sure its surface is dry by wiping the flasks with a paper towel before spraying 70% EtOH on them

16. Pour each of the medium mixtures into their square petri dishes (around 50ml)

17. Transfer the other half of the flasks into the BSC

18. Pour each of the medium mixtures into their square petri dishes

19. Let the petri dishes solidify at room temperature under the BSC

Done by Abby K, Abby T, Audrey H

2026/7/5



Date: 07/05

Title: Colony PCR

Purpose: Amplify part from 6/29 and screen bacterias for successful DNA inserts

Procedure:

1. Prepare the following reaction (master mix) x22 reactions and dispense solutions

Into PCR tubes	1rxn	22 rxn
Forward primer (Exp seq-fw 10mM)	1 µL	22µL
Reverse primer (Exp-rv 10mM)	1µL	22µL
Q amp master mix	10µL	220µL
ddH ₂ O	8µL	176µL

2. Circle 21 colonies each from plate C
(Exp-spacer-spacer-spacer-Pveg-EEBD-RBS0-Csn:AmilCP-B0015)
3. Pick up 21 colonies from each agar plate with tip and streak the labelled grid
4. Dispense the tips with colonies into the PCR strips and wait 5 minutes
5. Take out the tips and close the lid
6. Place into PCR machine for thermal cycle

95°C	2 min	1 cycle
95°C	30 secs	30 cycles
65°C	4 secs	30 cycles
72°C	2 min (bp 1800)	30 cycles
72°C	1 min	1 cycle

Done by Chloe, Abigail, Abby K

Supervised by Gabriel

Date: 07/05

Title: Miniprep

Purpose: plasmid purification

Procedure:

A: Pveg-EEBD-RBS0-Csn:AmilCP-B0015

B: PcpG2-veg-EEBD-RBS0-Csn:AmilCP-B0015

C: Pveg-RBS0-CcaR-B0015 (B2)

1. Transfer 1.8mL of medium to 2mL microtube
2. Centrifuge for 30 sec at 23444 g
3. Discard the supernatant
4. Repeat steps 2-3 until all medium is collected. Remember to empty the medium in the pipette tip from the medium tube
5. Add 200 microliters of Mx1 buffer
6. Resuspend by vortexing or pipetting
7. Add 250 microliters of Mx2 buffer
8. Invert slowly 8-10 times. Incubate at room temp for 2 minutes
9. Add 350 microliters of Mx3 buffer
10. Invert slowly 4-5 times
11. Centrifuge for 5 min at 16000 g
12. Transfer around 800 microliters of supernatant into the assembled column
13. Centrifuge at 7000 g for 30 seconds and discard the flow through
14. Wash the column with 500 microliters of WN buffer. Centrifuge at 7000 g for 30 seconds and discard the flow through
15. Wash the column with 700 microliters of WS buffer. Centrifuge at 7000 g for 30 seconds and discard the flow through
16. Centrifuge for 8 minutes at 10,000 g
17. Place column into new 1.5mL microtube
18. Add 50 microliters ddH₂O to the center of membrane
19. Incubate for 3 min at room temp
20. Centrifuge at 10,000 g for 3 min

Results:

Samples	Concentration(ng/microliter)	A260/A230	A260/A280
A	331.0	1.623	2.122
B	473.0	1.708	2.131
C	147.0	1.157	1.986

Done by Chloe, Jacquely, Noah

Date: 07/05

Title: Make 1.0% agarose gel

Purpose: making gel for electrophoresis

Procedure:

1. Assemble the tray and tray and mold and prepare the comb
2. Weigh 0.4g agarose powder and 40mL 1x TAE solution
3. Add 1x TAE and powder to the 100mL conical flask
4. Cover the flask with plastic wrap using a tip to puncture holes on it
5. Microwave for 30sec then take out and agitate it

6. Microwave for 10sec and agitate it. Repeat this step until there is no particles left in the gel
7. Wait until the melted gel stops smoking but is still hot
8. Add 2 microliters of safety dye and agitate it
9. Pour the melted gel into the tray and assemble the comb
10. Leave until solidify in room temp

Done by Abby K, Abigail

Date: 07/05

Title: Gel electrophoresis

Purpose: verify colony PCR product size and miniprep product size

Procedure:

1. Add the following samples in this layout(Add 1 μ l dye to A, B, C)

<u>M</u>	<u>E1</u>	<u>E2</u>	<u>E3</u>	<u>E4</u>	<u>E5</u>	<u>E6</u>	<u>E7</u>	<u>E8</u>	<u>E9</u>	<u>E10</u>	<u>E11</u>	<u>E12</u>	<u>E13</u>	<u>E14</u>	<u>E15</u>	<u>E16</u>	<u>E17</u>	<u>E18</u>	<u>E19</u>	<u>E20</u>	
3	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	
<u>E21</u>	<u>A</u>	<u>B</u>	<u>C</u>																		
5	6	6	6																		

(volume: μ l)

M: 1KB DNA ladder

E1-E21: EXP2 - spacer - spacer - spacer - Pveg - EEBD - RBS0 - Csn::AmilCP - B0015


A: Pveg - EEBD - RBS0 - csn:AmilCP - B0015

B: PcpG2 - veg - EEBD - RBS) - csn::AmilCP - B0015 -Amp

C: Pveg - RBS0 - ccaR - B0015

2. Resolve DNA at 100V for 40 minutes

Results:

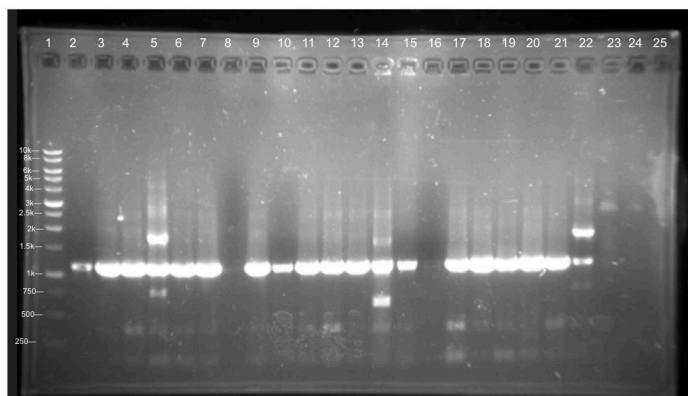

1. 1kb DNA ladder
- 2-22. EXP2-spacer-spacer-spacer-Pveg-EEBD-RBS0-csn::AmilCP-B0015
23. Pveg-EEBD-RBS0-csn:AmilCP-B0015
24. PcpG2-veg-EEBD-EBS0-csn::AmilCP-B0015-Amp
25. Pveg-RBS0-CcaR-B0015

Done by Abby K, Sophie H

Supervised by Gabriel

Date: 07/05

Title: Liquid Culture

Purpose: bacteria culture

Procedure:

1. Label 7 liquid culture tubes
2. Turn on bunsen burner
3. Add 3mL of TB medium and 3 μ L of ampicillin into 4 tubes; 3 μ L of chloramphenicol into 3 tubes
4. Dip a pipette tip into the colonies
5. Dispense tips into their respective tubes
6. Place tubes in the incubator at a slanted angle

Done by Jacquelyn and Chloe

Supervised by Chars

2026/7/6



Date: 20260706

Title: Colony pcr(redone)

Purpose: amplify parts from 6/29 and screen bacterias for successful DNA inserts

Procedure:

1. Prepare the following reaction(master mix) x 12 reactions and dispense solutions into PCR tubes

For A construct(Pveg-RBS0-Ho1:pcyA-B0015)

	1 reaction	13 reactions(for pipetting error)
Forward primer: pSTK-seq-fw 10mM	1 microliter	13 microliters
Reverse primer: pSTK-seq-rv 10mM	1 microliter	13 microliters
Q amp master mix	10 microliters	130 microliters
ddH ₂ O	8 microliters	104 microliters

For F construct(Pveg-EEBD-RBS0-GanA-Gala-B0015)

	1 reaction	13 reactions(for pipetting error)
Forward primer: pEXP-seq-fw 10mM	1 microliter	13 microliters
Reverse primer: pEXP-seq-rv 10mM	1 microliter	13 microliters
Q amp master mix	10 microliters	130 microliters
ddH ₂ O	8 microliters	104 microliters

2. Circle 12 colonies each from plate A and 3 colonies(F3, F4, F7) from plate F(chose the 3 colonies from F because others are empty)
 - a. A: Pveg-RBS0-Ho1:pcyA-B0015
 - b. F: Pveg-EEBD-RBS0-GanA-Gala-B0015
3. Pick up 12 colonies from A and 3 from F with tip and dispense in the PCR tubes(take colonies from incubated plates)
4. Take out tips and close the lid
5. Place into PCR machine for thermal cycle

95°C	2 min	1 cycle
------	-------	---------

95°C	30 sec	30 cycles
65°C	40 sec	
72°C	2.5 min	
72°C	1 min	1 cycles

Done by Chloe, Abby K, Abigail

Date: 20260706

Title: Prepare 0.8% Agarose gel for DNA electrophoresis (x2; one redone)

Purpose: DNA electrophoresis for A:Pveg-RBS0-Ho1:pcyA-B0015 and F:Pveg-EEBD-RBS0-GanA-Gala-B0015

Procedure:

1. Assemble the tray and tray and mold and prepare the comb
2. Weigh 0.32g agarose powder and 40mL 1x TAE solution
3. Add 1x TAE and powder to the 100mL conical flask
4. Cover the flask with plastic wrap using a tip to puncture holes on it
5. Microwave for 30sec then take out and agitate it
6. Microwave for 10sec and agitate it. Repeat this step until there is no particles left in the gel
7. Wait until the melted gel stops smoking but is still hot
8. Add 2 microliters of safety dye and agitate it
9. Pour the melted gel into the tray and assemble the comb
10. Leave until solidify in room temp

Done by Abigail, Abby T

Date: 20260706

Title: Gel electrophoresis for 0705 colony pcr

Purpose: validate the bps of pcr results

Procedure:

1. Load the gel with the following layout
A1(5 microliters) A2(5 microliters) A3(5 microliters) A4(5 microliters) A5(5 microliters) A6(5 microliters) A7(5 microliters) A8(5 microliters) A9(5 microliters) A10(5 microliters) A11(5 microliters) A12(5 microliters) 1kb DNA ladder(3 microliters) F1(5 microliters) F2(5 microliters) F3(5 microliters) F4(5 microliters) F5(5 microliters) F6(5 microliters) F7(5 microliters) F8(5 microliters) F9(5 microliters) F10(5 microliters) F11(5 microliters) F12(5 microliters)


A1-A12: Pveg-RBS0-Ho1:pcyA-B0015

F1-F12: Pveg-EEBD-RBS0-GanA-Gala-B0015

2. Resolve at 100v for 40min

Results:

1.

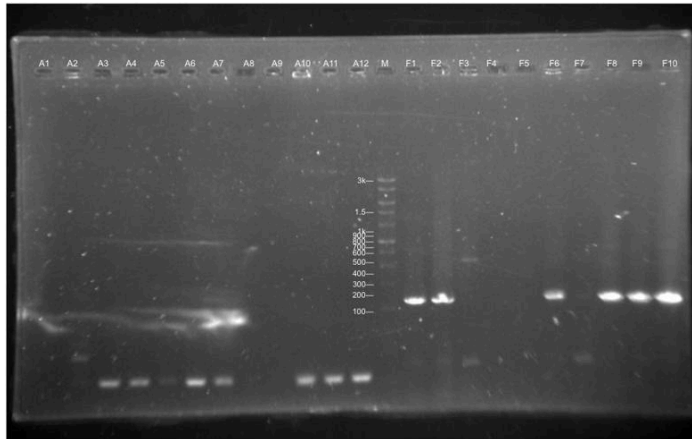


1-12. Pveg-RBS0-Ho1:pcyA-B0015

13. 1kb DNA ladder

14-25. Pveg-EEBD-RBS0-GanA-Gala:B0015

2. Redid the previous gel electrophoresis

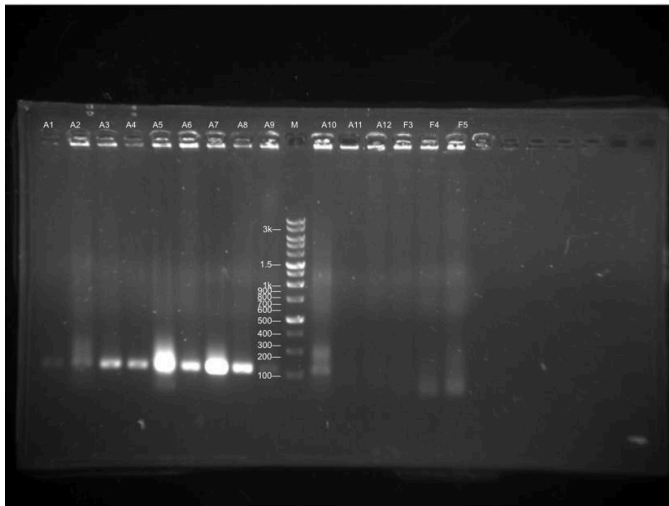

M: [DM2300] ExcelBand™ 100 bp+3K DNA Ladder 5 uL

A1-12: Pveg-RBS0-Ho1:pcyA-B0015

F3-F5: Pveg-EEBD-RBS0-GanA-Gala-B0015

Done by Abigail, Abby T

Date: 20260706

Title: Miniprep

Purpose:

Procedure:

C4: Exp2-spacer3x Pveg-EEBD-RBS0-csn:AmilCP-B0015

C21: Exp2-spacer3x Pveg-EEBD-RBS0-Csn:AmilCP-B0015

C13: Exp2-spacer3x Pveg-EEBD-RBS0-Csn:AmilCP-B0015

B2: Pveg-RBS0-CcaR-B0015

G3: PcpcG2-veg-EEBD-RBS0-Csn:AmilCP-B0015

H10: Pveg-EEBD-RBS0-Csn:AmilCP-B0015

1. Transfer 1.8mL of medium to 2mL microtube
2. Centrifuge for 30 sec at 23444 g

3. Discard the supernatant
4. Repeat steps 2-3 until all medium is collected. Remember to empty the medium in the pipette tip from the medium tube
5. Add 200 microliters of Mx1 buffer
6. Resuspend by vortexing or pipetting
7. Add 250 microliters of Mx2 buffer
8. Invert slowly 8-10 times. Incubate at room temp for 2 minutes
9. Add 350 microliters of Mx3 buffer
10. Invert slowly 4-5 times
11. Centrifuge for 5 min at 16000 g
12. Transfer around 800 microliters of supernatant into the assembled column
13. Centrifuge at 7000 g for 30 seconds and discard the flow through
14. Wash the column with 500 microliters of WN buffer. Centrifuge at 7000 g for 30 seconds and discard the flow through
15. Wash the column with 700 microliters of WS buffer. Centrifuge at 7000 g for 30 seconds and discard the flow through
16. Centrifuge for 8 minutes at 10,000 g
17. Place column into new 1.5mL microtube
18. Add 50 microliters ddH₂O to the center of membrane
19. Incubate for 3 min at room temp
20. Centrifuge at 10,000 g for 3 min

The results are measured by Biodrop

Samples	Concentration(ng/microliter)	A260/A280
C4	88	1.833
C13	90	1.837
C21	58	1.871
B2	134	1.861
G3	168	1.846
H10	184	min

Done by Sophie C and Sophia Y

Date: 20260706

Title: Gel electrophoresis(0.8% agarose gel)

Purpose: Validate the miniprep results

Procedure:

1. Assemble tray and mold, prepare the comb
2. Weigh 0.16g agarose powder and 20mL of 1x TAE solution
3. Add 1x TAE and powder to 100mL conical flask

4. Cover flask with the plastic wrap and puncture holes on top with a tip
- 5.
6. Microwave for 30 seconds, take it out and agitate it
7. Microwave 10 seconds and agitate and repeat until no particles are left in the melted gel
8. Wait until gel stops smoking and add 1 microliters of safety dye and agitate it
9. Pour the melted gel in the tray and assemble the comb
10. Leave to solidify in room temperature

Done by Sophie C

Date: 20260706

Title: Golden gate assembly

Purpose: assemble level 1 genetic circuit

Procedure:

A: (STK) Module 1(1A): Pveg–RBS0-Ho1:pcyA-B0015

1. prepare the reaction as follow:

Part	Original Amount	Dilution
Vector (1A) 115 ng/μL, 3120 bp	0.176348 μL	1.8 μL (dilute 1/10)
Insert promoter (Pveg) 280 ng/μL, 2361 bp	0.109618 μL	1.1 μL (dilute 1/10)
Insert RBS(RBS0) 188 ng/μL, 2113 bp	0.1461117 μL	1.5 μL (dilute 1/10)
Insert CDS(Ho1::PcyA) 238 ng/μL, 1556 bp	0.0849916 μL	1.7 μL (dilute 1/20)
Insert Terminator (B0015) 107 ng/μL, 2213 bp	0.268869 μL	1.3 μL (dilute 1/5)
T4 DNA ligase buffer	2 μL	2 μL
Sterile water	15.2140617 μL	8.6 μL
T4 DNA ligase	1 μL	1 μL

Bsal-HFv2	1 μ L	1 μ L
Total	20 μ L	20 μ L

Use the thermocycler to run the program below:

Temperature	Time	Cycles
37 °C (Restriction enzyme digestion)	5 minutes	30 cycles
16 °C (Ligation)	5 minutes	
37 °C (Restriction enzyme digestion)	10 minutes	1 cycle
80 °C (Heat inactivation of Bsal restriction enzyme)	20 minutes	1 cycle

B: (STK) Module 4-Pveg - EEED- RBS0 - CSN::ACCD - B0015

Vector (Exp1) - concentration(ng/μL): 162 - bp: 6929	0.278015 μ L	1.4 μ L (dilute $\frac{1}{5}$)
Insert promoter(Pveg) - concentration(ng/μL): 280 - bp: 2361	0.109618 μ L	1.1 μ L (Dilute 1/10)
Insert RBS(EEED-RBS0) - concentration(ng/μL): 150 - bp: 2392	0.207306 μ L	1.0 μ L (Dilute 1/5)
Insert CDS(CSN::ACCD) concentration(ng/μL): 105	0.147952 μ L	1.5 μ L(dilute 1/10)

bp: 1195		
Insert Terminator (B0015) - concentration(ng/μL): 107 - bp: 2213	0.268869 μL	1.3 μL (dilute 1/5)
T4 DNA ligase buffer	2 μL	2 μL
Sterile water		8.6 μL
T4 DNA ligase	1 μL	1 μL
BsaI-HFv2	1 μL	1 μL
Total	20 μL	20 μL

Use the thermocycler to run the program below:

Temperature	Time	Cycles
37 °C (Restriction enzyme digestion)	5 minutes	30 cycles
16 °C (Ligation)	5 minutes	
37 °C (Restriction enzyme digestion)	10 minutes	1 cycle
80 °C (Heat inactivation of BsaI restriction enzyme)	20 minutes	1 cycle

C: (STK) Module 4-Pveg - EEED- RBS0 - CSN::LEA - B0015

Vector (Exp1)	0.278015 μL	1.4 μL (dilute 1/5)
---------------	-------------	---------------------

- concentration(ng/μL): 162 - bp: 6929		
Insert promoter(Pveg) - concentration(ng/μL): 280 - bp: 2361	0.109618 μL	1.1 μL (Dilute 1/10)
Insert RBS(EEBD-RBS0) - concentration(ng/μL): 150 - bp: 2392	0.207306 μL	1.0 μL (Dilute 1/5)
Insert CDS(CSN::LEA) 77 ng/μL 634 bp	0.05351948051 μL	1.1 μL (dilute 1/20)
Insert Terminator (B0015) - concentration(ng/μL): 107 - bp: 2213	0.268869 μL	1.3 μL (dilute 1/5)
T4 DNA ligase buffer	2 μL	2 μL
Sterile water		10.0 μL
T4 DNA ligase	1 μL	1 μL
BsaI-HFv2	1 μL	1 μL
Total	20 μL	20 μL

Use the thermocycler to run the program below:

Temperature	Time	Cycles
37 °C (Restriction enzyme digestion)	5 minutes	30 cycles

16 °C (Ligation)	5 minutes	
37 °C (Restriction enzyme digestion)	10 minutes	1 cycle
80 °C (Heat inactivation of BsaI restriction enzyme)	20 minutes	1 cycle

D: (STK-EXP1) Pveg-EEBD-RBS0-SamyQ::TAT-Ha::AmilCP-B0015

Part	Original Amount	Dilution
Vector (shuttle vector): pSTK-EXP-1-sfGFP - concentration(ng/μL): 162 - bp: 6929	0.278015 μL	1.4 μL (dilute 1/5)
Insert promoter(Pveg) - concentration(ng/μL): 280 - bp: 2361	0.109618 μL	1.1 μL (dilute 1/10)
Insert RBS(EEBD-RBS0) - concentration(ng/μL): 150 - bp: 2392	0.207306 μL	1.0 μL (Dilute 1/5)
Insert CDS(SamyQ::TAT-Ha::AmilCP) - concentration(ng/μL): 52 - bp: 899	0.22475 μL	1.1 μL (dilute 1/5)
Insert Terminator (B0015) - concentration(ng/μL): 107 - bp: 2213	0.268869 μL	1.3 μL (dilute 1/5)

T4 DNA ligase buffer	2 µL	2 µL
Sterile water	14.9114418	10.0 µL
T4 DNA ligase	1 µL	1 µL
BsaI-HFv2	1 µL	1 µL
Total	20 µL	20 µL

Use the thermocycler to run the program below:

Temperature	Time	Cycles
37 °C (Restriction enzyme digestion)	5 minutes	30 cycles
16 °C (Ligation)	5 minutes	
37 °C (Restriction enzyme digestion)	10 minutes	1 cycle
80 °C (Heat inactivation of BsaI restriction enzyme)	20 minutes	1 cycle

Done by: Abigail, Chloe, Abby T

2026/7/7



Date: 20260707

Topic: Arabidopsis transformation

Purpose: To verify and transform the healthy Arabidopsis from the rest

Procedure:

1. Place the cap holder in 2000ml of 10% bleach for 10 min (200ml bleach and 1800ml of water)
2. After bleached for 10 min, take the cap holder out and cleanse it under ddH₂O tap (slightly tilt the holder to let the water flow through every corner of it) until no bleach smell left (at least 3 min)
3. Placed the cleansed holder in biosafety cabinet along with the plants you want to transform and tip box (blue)
4. Assemble the cap holder with tip box, than pour MS medium inside the tip holders to almost full, but never pour over the front side of the Arabidopsis holder
5. Select and place the desired plants on the cap holder

	1	2	3	4	5	6	7	8	9	10	11	12
A	DT	SPA	PA	PA	SPA	SPA	PA	PA	PA	PA	PA	PA
B	DT	PA	SPA	DT	PA	PA	DT	PA	SPA	SPA	DT	SPA
C	DT	DT	SPA	DT	PA	DT	PA	PA	DT	DT	DT	SPA
D	PA	PA	DT	DT	DT	DT	DT	SPA	DT	DT	DT	DT

*DT: Didn't Transform *PA(17caps): Prioritize Arabatopsis *SPA(10caps): Second Prioritize Arabatopsis

Done by Sarah, Sophie H, Sophie C, Alex, Naomi

Supervised by Alex, Gabriel

Recorded by Sophie H

Date: 20260707

Topic: Making ½ MS medium

Purpose: Make a total of 1 liter of hydroponic medium for hydroponic box

1. Measure and add to a flask:

½ MS medium (without sucrose)	2.3g
MES buffer	0.5g
ddH ₂ O	1L

2. Use the hot plate to stir (with a stirring rod) until no sediments are left in the solution
3. Determine pH of the solution and adjust the pH using a pH probe to around 5.6-5.7 using 1M KOH
 - a. Record original and final pH
 - i. original: 4.0
 - ii. final: 5.6
4. Wrap the neck with aluminum foil and place the indicator tape on top
5. Autoclave at 127°C for 20 minutes
6. Wait for the pressure to reach 0 before opening the autoclave
7. Pour 300ml of medium into each of the hydroponic boxes

Done by Abby K/T

Date: 20260707

Topic: Transfer arabidopsis to the salinity plates

Purpose: transfer the germinated seedlings from neutral plates to salinity plates to conduct salinity stress tests (experimental set 5).

Procedure:

1. Draw a standardized baseline on each of the salinity plates to signal where the roots will be lined up after transferring
2. Place all equipments into the BSC after sanitizing
3. Get the neutral plate inside the BSC one at a time
4. Lay a template sheet of paper underneath the salinity plate so that it is clear where to place the seedlings
5. Use the tweezer to slowly lift the arabidopsis plant from the top of the roots
6. Place the seedling softly onto the salinity plate, with the bottom of the root precisely onto the baseline ; 9 seedlings per row, per plate.
 - a. make sure to pick seedlings that have around the same root length
7. Repeat step 3-6 until all seedlings transferred into all of the salinity plates
8. Place into the incubator

Done by Abby K/T

Date: 20260707

Topic: Remake the incorrect concentrations salinity plates

Purpose: prepare salinity plates for experimental set 5

1. Prepare the following in a 500ml flask:

½ MS medium powder(w/o sucrose)	0.69g
MES buffer	0.15g
Distilled water	300 ml
Sucrose(1%)	3 g

2. Use a hot plate to stir (with a magnetic stirring rod) until no sediments are left in the solution
3. Use the pH probe and measure the pH of the solution, make sure the pH is 5.7-5.8
 - a. If the pH is too low, add 0.1M KOH slowly to increase it
4. Pour 50ml of mixture into the 6 flasks
5. Prepare 6 flasks and pour the respective trehalose and NaCl concentration in each of the flasks
 - a. making sure there is 50ml in each of the flasks (6)
 - b. label each flask clearly according to their trehalose concentrations
6. Add:

150mM NaCl, 5mM trehalose

trehalose (µL)	NaCl (g)	Agar powder (g)
2500	0.548	0.5

150mM NaCl, 10 trehalose

trehalose (µL)	NaCl (g)	Agar powder (g)
5000	0.548	0.5

250mM NaCl, 5mM trehalose

trehalose (µL)	NaCl (g)	Agar powder (g)
2500	0.913	0.5

250mM NaCl, 0mM trehalose

trehalose (µL)	NaCl (g)	Agar powder (g)
0	0.913	0.5

0mM NaCl, 0mM trehalose

trehalose (μL)	NaCl (g)	Agar powder (g)
0	0	0.5

0mM NaCl, 1mM trehalose

trehalose (μL)	NaCl (g)	Agar powder (g)
500	0	0.5

7. Wrap the necks of the flasks with aluminum foil and stick a piece of indicator tape on top
8. Autoclave for 20 minutes at 127°C
 - a. Wait for the pressure to reach 0 before retrieving the flask

During the autoclave, label 6 of the petri dishes with a marker in the biosafety cabinet

- b. right top: date and arrow pointing up to show direction for pictures
 - c. bottom lid: different trehalose concentrations and NaCl concentration or control
2. Transfer half of the flasks to BSC
 - a. Store the other half of the flasks in the hot bath while pouring so that it does not solidify
3. Make sure its surface is dry by wiping the flasks with a paper towel before spraying 70% EtOH on them
4. Pour each of the medium mixtures into their square petri dishes (around 50ml)
5. Transfer the other half of the flasks into the BSC
6. Pour each of the medium mixtures into their square petri dishes
7. Let the petri dishes solidify at room temperature under the BSC

Done by Abby K/T

Date: 20260707

Topic: Transformation

Purpose: Transform and cultivate engineered bacteria on agar plate

Procedure:

A	Pveg-RBS0-CcaR-B0015 (3243
B	PcpcG2-veg-EEBD-RBS0-Csn:AmilCP-B0015 (3485
C	PliaG-RBS0-CcaSm3-B0015 (3792
D	Pveg-EEBD-RBS0-Csn:AmilCP-B0015(7524
E	pGEM-Teasy-Csn:LEA

F	pGEM-Teasy-ccaR
G	pGEM-Teasy-csn:amilCP

1. Mix 15µL competent cell with 1µL plasmid
2. Incubate in ice for 20 min
3. Heat shock at 42°C for 45 sec
4. Place on ice for 3min to recover
5. Spread the mixture on TBA plate
6. Place in incubator at 25°C overnight

Done by Audrey, Felix, Noah, Chloe
Supervised by Gabriel

Date: 20260707

Topic: Transformation

Purpose: Transform and cultivate engineered bacteria on agar plate

A: Pveg-RBS0-CcaR-B0015

B: PcpG2-veg-EEBD0RBS0-Csn:AmilCP-B0015

C: PliaG-RBS0-CcaSm3-B0015

D: Pveg-EEBD-RBS0-Csn:AmilCP-B0015

E: pGEM-Teasy-Csn:LEA

F: pGEM-Teasy-CcaR

G: pGEM-Teasy-Csn:AmilCP

1. Mix 15µL competent cell with 1µL plasmid
2. Incubate in ice for 20 min
3. Heat shock at 42°C for 45 sec
4. CDsPlace on ice for 3min to recover
5. Spread the mixture on TBA plate
6. Place in incubator at 25°C overnight

Done by Sophie H, Naomi, Alex, Sophie C\
Supervised by Chars

2026/7/8 

Date: 20260708

Topic: Colony PCR

Purpose: To verify bacterial colonies after a transformation step in molecular cloning

Procedure:

1. Prepare the following reaction (master mix) and dispense solutions into PCR tubes

For samples A, B, C, D, and A1

Component	1 reaction	14 reactions
Forward primer: T7 promoter	1 μ L	14 μ L
Reverse primer: SP6 promoter	1 μ L	14 μ L
Qamp master mix	10 μ L	140 μ L
ddH ₂ O	8 μ L	112 μ L
Total	20 μ L	280 μ L

For samples E, F, and G

Component	1 reaction	10 reactions
Forward primer: T7 promoter	1 μ L	10 μ L
Reverse primer: SP6 promoter	1 μ L	10 μ L
Qamp master mix	10 μ L	100 μ L
ddH ₂ O	8 μ L	80 μ L
Total	20 μ L	200 μ L

For samples *B, *C, *D

Component	1 reaction	8 reactions
Forward primer: T7 promoter	1 μ L	8 μ L
Reverse primer: SP6 promoter	1 μ L	8 μ L
Qamp master mix	10 μ L	80 μ L
ddH ₂ O	8 μ L	64 μ L

Total	20 μ L	160 μ L
-------	------------	-------------

2.

a. Circle 3 colonies from each of the 7 plates plus 1 from A1

A: Pveg – RBS0 – CCaR – B0015

B: PcpG2 – veg – EEBD – RBS0 – Csn::AmiCP – B0015

C: PliaG – RBS0 – CcaSm3 – B0015

D: Pveg – EEBD – RBS0 – Csn::AmiCP – B0015

E: pGEM-TEASY – Csn::LEA

F: pGEM-TEASY – CCaR

G: pGEM-TEASY – Csn:: AmiCP

A1: Pveg – RBS0 – Ho1-psyA-B0015

b. Circle 1 colony each on plate *B and *C and 5 on plate *D

*B: Pveg-EEBD-RBS0-Csn::ACCD-B0015

*C: Pveg-EEBD-RBS0-Csn::LEA-B0015

*D: Pveg-EEBD-RBS0-SamyQ::TAT-HA::AmiCP-B0015

3. Pick up colonies from each of the plates with a tip and dispense into PCR tubes (colonies from incubated plates).

4. Remove tips and close the lid

5. Place into PCR machine for thermal cycle

95°C	2 min	1 cycle
95°C	30 secs	30 cycles
55°C	4 secs	30 cycles
72°C	2 min (bp 1800)	30 cycles
72°C	1 min	1 cycle

Done by: Sophie C, Abigail, Naomi

Supervised by: Gabriel

Date: 20260708

Title: 1.0% Agarese Gel preparation

Purpose: To prepare 1.0% Agarese Gel for further electrophoresis experiment

Procedure:

1. Prepare and assemble the mold and comb
2. Weight 0.4g of agarose powder and add 40ml of 1X TAE solution into flask
3. Cover the flask with plastic wrap, then use a tip to puncture holes on the plastic wrap
4. Microwave until the gel stops smoking but is still hot
5. Add 2µL of safety dye into the gel and agitate it
6. Pour the melted gel into the tray
7. Wait until the gel cools down
8. Repeat steps 1~8 to create second copy of gel

Done by Sophie L, Abby T

Supervised by Gabriel

Date: 20260708

Topic: Gel Electrophoresis

Purpose: verify the sizes of the colony PCR products

Procedure:

1. Load the samples in the following layout:

Gel 1:

<u>M</u>	<u>A1</u>	<u>A2</u>	<u>A3</u>	<u>B1</u>	<u>B2</u>	<u>B3</u>	<u>C1</u>	<u>C2</u>	<u>C3</u>	<u>D1</u>	<u>D2</u>	<u>D3</u>	<u>E1</u>	<u>E2</u>	<u>E3</u>	<u>F1</u>	<u>F2</u>	<u>F3</u>	<u>G1</u>	<u>G2</u>	<u>G3</u>	<u>*A1</u>	<u>*A2</u>	<u>M</u>
5	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	5

M: 1kb DNA ladder

A1-3: Pveg – RBS0 – CCaR – B0015

B1-3: PcpcG2 – veg – EEED – RBS0 – Csn::AmiCP – B0015

C1-3: PliaG – RBS0 – CcaSm3 – B0015

D1-3: Pveg – EEED – RBS0 – Csn::AmiCP – B0015

E1-3: pGEM-TEASY – Csn::LEA

F1-3: pGEM-TEASY – CCaR

G1-3: pGEM-TEASY – Csn:: AmiCP

*A1-*A2: Pveg – RBS0 – Ho1-psyA-B0015

Results:

Date: 20260708

Topic: Liquid Culture

Purpose: bacteria culture from successful colony PCR product

Procedure:

1. Label 9 liquid culture tubes
2. Turn on bunsen burner
3. Add 3ml of TB medium and 3 μ L of chloramphenicol into C4 ;3ml of TB medium and 3 μ L ampicillin into rest of the tubes
4. Dip a pipette tip into the colonies
5. Dispense tips in their respective tubes
6. Place the tubes in the incubator at a slanted angle around 45-50°

A1: Pveg-RBS0-CCaR-B0015

B1: PcpG2-veg-EEBD-RBS0-Csn::AmiCP-B0015

C1: PliaG-RBS0-CcaSm3-B0015

D2: Pveg-EEBD-RBS0-Csn::AmiCP-B0015

E1: pGEM-TEASY-Csn::LEA

F1: pGEM-TEASY-CCaR

G1: pGEM-TEASY-Csn:: AmiCP

*D2: Pveg-EEBD-RBS0-SamyQ::TAT-HA::AmiCP-B0015

C4: Exp2-spacer-spacer-spacer-Pveg-EEBD-RBS0-csn::AmiCp-B0015

Done by: Abby T

Supervised by: Gabriel

2026/7/9



Date: 20260709

Title: Miniprep

Purpose: Extract plasmid DNA from the liquid culture of colony PCR products for amplification

Procedure:

1. Prepare Glycerol stock, formula as follows:

	Exp2-spacer-spacer-spacer-Pveg-EEBD-RBS0-csn::AmiCP-B0015
50% Glycerol Stock	500 µL
Bacteria Culture Medium	500 µL

Then store in -80 deg C.

A1: Pveg-RBS0-CCaR-B0015

B1: PcpcG2-veg-EEBD-RBS0-Csn::AmiCP-B0015

C1: PliaG-RBS0-CcaSm3-B0015

D2: Pveg-EEBD-RBS0-Csn::AmiCP-B0015

E1: pGEM-TEASY-Csn::LEA

F1: pGEM-TEASY-CCaR

G1: pGEM-TEASY-Csn:: AmiCP

*D2: Pveg-EEBD-RBS0-SamyQ::TAT-HA::AmiCP-B0015

C4: Exp2-spacer-spacer-spacer-Pveg-EEBD-RBS0-csn::AmiCp-B0015

Results:

	Concentration(ng/µL)	A260/230	A260/230
A. pSTK-1C-Pveg-RBS0-CcaR-B0015	126	1.145	2.000
B. pSTK - PcpcG2 - v eg - EEBD - RBS0 - csn::AmiCP - B0015	186	1.094	1.958
C. pSTK - PliaG - RBS0 - Ccasm3 - B0015	134	1.411	2.062
D. pSTK - Preg - EEBD - RBS0 - csn::AmiCP - B0015	11	1.181	1.375

E. pGEM-TEASY - csnn::LEA	150	0.848	1.948
F. pGEM-TEASY - CcaR	112	1.008	1.800
G. pGEM-TEASY - csnn::AmilCP	128	1.356	1.910
*D2. Pveg - EEED - RBS0 - SamyQ::TAT-HA::AmilCP - B0015	160	1.552	2.000
C4. Exp2 - spacer - spacer - spacer - Pveg - EEED - RBS0 - csnn::AmilCP - B0015	72	0.524	1.940

Done by Abby T/K, Naomi, Sarah C
Supervised by Chars

Date: 20260709

Title: Colony PCR

Purpose: To amplify bacterial sequences after a transformation step in molecular cloning

Procedure

D: Pveg-EEED-RBS0-SamyQ-TAT-HA-AmilCP-B0015(Exp1)

A: 1A-Pveg-RBS0-Ho1:pcyA-B0015

1. Prepare the following reaction (master mix) and dispense solutions into PCR tubes

A	1 rxn	13 rxn
Fw Primer (pSTK-seq-fw, 10mM)	1µL	13 µL
Rv Primer (STK-seq-rv, 10mM)	1µL	13µL
Qamp master mix	10µL	130µL
ddH ₂ O	8µL	104µL

D	1 rxn	13 rxn
Fw primer	1µL	13 µL

(STK-Exp1-seq-fw, 10mM)		
Rv Primer (STK-Exp1-seq-rv, 10mM)	1µL	13µL
Qamp master mix	10µL	130µL
ddH ₂ O	8µL	104µL

- Circle 12 colonies with no fluorescent signal under UV light from plates A and D
A: 1A- Pveg – RBS0 – Ho1:pcyA – B0015
D: Pveg – EEED – RBS0 – SamyQ – TAT – HA – AmilCP – B0015 (Exp1)
- Pick up colonies from each of the plates with a tip and dispense into PCR tubes (colonies from incubated plates).
- Remove tips and close the lid
- Place into PCR machine for thermal cycle

95°C	5 min	1 cycle
95°C	20 secs	35 cycles
55°C	20 secs	
72°C	2.5 min (bp 2200)	
72°C	5 min	1 cycle

Done by: Sophie C., Audrey C. , Sophie H.

Supervised by: Chars

Recorded by Sophie H.

Date: 20260709

Title: 1.0% Agarose Gel preparation (x2 gels)

Purpose: To prepare 1.0% Agarose Gel for electrophoresis of miniprep

Procedure:

- Prepare and assemble the mold and comb
- Weight 0.4g of agarose powder and add 40ml of 1X TAE solution into flask
- Cover the flask with plastic wrap, then use a tip to puncture holes on the plastic wrap
- Microwave until the gel stops smoking but is still hot
- Add 2µL of safety dye into the gel and agitate it
- Pour the melted gel into the tray

7. Wait until the gel cools down
8. Repeat steps 1~8 to create second copy of gel

Done by Abby T/K, Sophie C, Audrey C, Sophie H
Supervised by Chars

Date: 20260709

Title: Gel Electrophoresis

Purpose: Verify the sizes of the cloning PCR product of bacterial colony A and D

Procedure:

1. Load the samples in the following layout:

<u>M</u>	<u>D6</u>	<u>D7</u>	<u>D8</u>	<u>D9</u>	<u>D10</u>	<u>D11</u>	<u>D12</u>	<u>D13</u>	<u>D14</u>	<u>D15</u>	<u>D16</u>	<u>D17</u>	<u>A1</u>	<u>A2</u>	<u>A3</u>	<u>A4</u>	<u>A5</u>	<u>A6</u>	<u>A7</u>	<u>A8</u>	<u>A9</u>
5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5
												<u>A10</u>	<u>A11</u>	<u>A12</u>							
												5	5	5							


M: 3 kb DNA ladder

D6-D17: Pveg - EEED - RBS0 - SamyQ - TAT - HA - AmilCP - B0015(Exp1)

A1 - A12: 1A - Pveg - RBS0 - Ho1:PcyA - B0015

2. Resolve at 100V for 30 minutes

20260709 Colony PCR A D

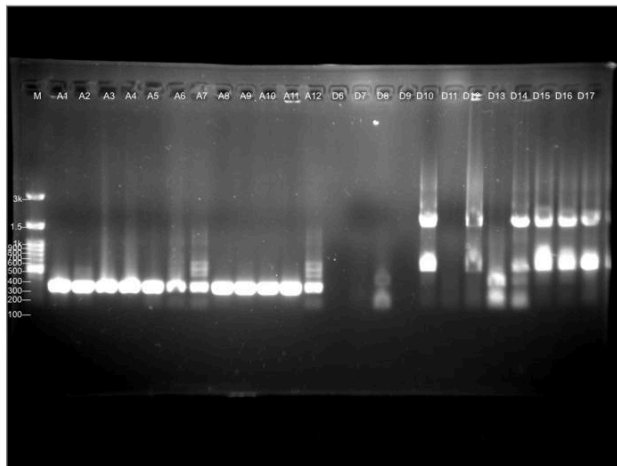

M: [DM2300] ExcelBand™ 100 bp+3K DNA Ladder 5 uL
A: 1A-Pveg-RBS0-H01:PcyA-B0015
D:Pveg-EEED-RBS0-SamyQ-TAT-HA-AmilCP-B0015(Exp1)

Done by Sophie C, Sarah, Sophie H, Naomi

Date: 20260709

Topic: Gel Electrophoresis

Purpose: verify the sizes of the miniprep product

Procedure:

1. Load the samples in the following layout:

Gel 1:

M	A	B	C	D	E	F	G	*D2	C4
3	6	6	6	6	6	6	6	6	6

A: pSTK-1C - Pveg - RBS0 - CcaR - B0015

B: pSTK - PcpcG2 - v eg - EEED - RBS0 - csn::AmilCP - B0015

C: pSTK - Pliag - RBS0 - Ccasm3 - B0015

D: pSTK - Preg - EEED - RBS0 - csn::AmilCP - B0015

E: pGEM-TEASY - csn::LEA

F: pGEM-TEASY - CcaR


G: pGEM-TEASY - csn::AmilCP

D2: Pveg - EEED - RBS0 - SamyQ::TAT-HA::AmilCP - B0015

C4 :Exp2 - spacer - spacer - spacer - Pveg - EEED - RBS0 - csn::AmilCP - B0015

2. Resolve the DNA at 100V for 30 minutes

20260709 miniprep

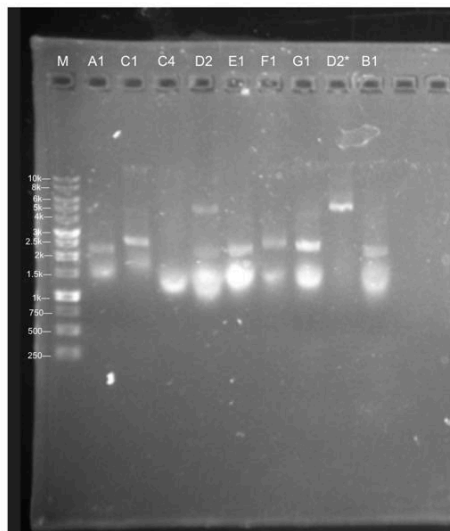

M: 1kb DNA ladder
A1: pSTK-1C-Pveg-RBS0-CcaR-B0015
C1: pSTK-PilaG-RBS0-CcaSm3-B0015
C4:
Exp2-spacer-spacer-spacer-Pveg-EEED-RBS0-csn::AmilC
P-B0015
D2: pSTK-Pveg-EEED-csn::AmilCP-B0015
E1: pGEM-TEASY-csn:LEA
F1: pGEM-TEASY-CcaR
G1: pGEM-TEASY-csn:AmilCP
D2*: Pveg-EEED-RBS0-SamyQ::TAT-HA::AmilCP-B0015
B1: pSTK-PcpcG2- veg-EEED-RBS0-csn:AmilCP-B0015

Done by Abby K, Abby T, Audrey Chu

Supervised by: Chars

Date: 20260709

Title: Golden gate assembly

Purpose:

Procedure:

A: pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015

1. Prepare the reaction as follow:

Part	Original Amount	Dilution
Insert vector(pSTK-1A) 115ng/ µl (bp:3120)	0.176348 µl	1.76 µl (1/10)
Insert Promoter(Pveg)	0.109618 µl	1.1 µl (1/10)
Insert RBS(RBS0)	0.1461117 µl	1.46 µl (1/10)
Insert CDS(Ho1::pcyA)	0.436455824 µl	1.31 µl ($\frac{1}{3}$)
Insert terminator(B0015)	0.268869 µl	1.34 µl ($\frac{1}{5}$)
T4 DNA ligase buffer	2 µl	2 µl
Sterile water	15.2040617 µl	9.03 µl
T4 DNA ligase	1 µl	1 µl
Bsal-HFv2	1 µl	1 µl
Total	20 µl	20 µl

2. Use the thermal cycler to run the program below:

Temperature	Time	Cycles
37 °C (Restriction enzyme digestion)	5 minutes	30 cycles
16 °C (Ligation)	5 minutes	
37 °C (Restriction enzyme digestion)	10 minutes	1 cycle
80 °C (Heat inactivation of Bsal restriction enzyme)	20 minutes	1 cycle

Done by Abby T/K, Audrey Chu, Sophie H, Sophie C, Naomi, Sarah
Supervised by Chars

Date: 20260709

Title: Liquid Culture

Purpose: bacteria culture

Procedure:

1. Label 4 liquid culture tubes
2. Turn on Bunsen Burner
3. Add 3 mL of TB medium into each tube
4. Add 3 μ l of chloramphenicol into C4 tube
5. Add 3 μ l of ampicillin into the rest of the 3 tubes(D10, D12, D14)
6. Dip a pipette tip into the colonies
7. Dispense tips in their respective tubes
8. Place the tubes in the incubator at a slanted angle

C4: spacer*3 - pveg - EEBD - RBS0 - csn:AmilCP - B0015 Exp2

D10, D12, D14: Pveg - EEBD - RBS0 - SamyQ - TAT - HA - AmilCP - B0015 Exp1

Done by Alex

Supervised by Chars

2026/7/10  (Only 1 task)

Date: 20260710

Title: Miniprep

Purpose: Isolate plasmid from DH5x for sequencing

Plasmid: C4: pSTK-EXP2-Spacerx3-Pveg-EEBD-RBS0-Csn:AmilCP-B0015

D10: pSTK-EXP1-Pveg-EEBD-RBS0-SamyQ: TAT-HA:AmilCP-B0015

D12: pSTK-EXP1-Pveg-EEBD-RBS0-SamyQ: TAT-HA:AmilCP-B0015

D14: pSTK-EXP1-Pveg-EEBD-RBS0-SamyQ: TAT-HA:AmilCP-B0015

Procedure:

1. Collect liquid culture 1.8mL in 2mL microtubes
2. Centrifuge at 23444g for 30 seconds. Discard the supernatant
3. Repeat step 1 until collect all of the liquid culture
4. Resuspend pellet in 250 microliters of P1 buffer(Qiagene)
5. Add 250 microliters of P2 buffer.
6. Invert the microtubes for 10 times. Incubate at RT for 2 min
7. Add 350 microliters of N3 buffer. Invert the microtubes for 10 times
8. Centrifuge at 13000g for 10 mins
9. Transfer supernatant 800 μ L to column and discard the pellet
10. Centrifuge at 17900g for 1 min. Discard the flowthrough.
11. Add 500 μ L PB buffer. Centrifuge at 17900g for 1 min and discard the flowthrough
12. Add 750 μ L PE buffer. Centrifuge at 17900g for 1 min and discard the flowthrough
13. Centrifuge the empty column at 17900g for 3 min
14. Place the column on the 1.5 mL microtube
15. Add 50 μ L sterile deionized water at 65°C
16. Incubate at RT for 5min
17. Centrifuge at 17900g for 3 min. Discard the column
18. Measure the concentration using Biodrop.

Biodrop results:

	concentration(ng/ μ l)	A260/280	A260/230
C4	94.0	1.880	2.186
D10	277.0	1.847	2.347
D12	306.0	1.913	2.354
D14	354.0	1.914	2.376

Done by Chars

2026/7/12 

Date: 20260712

Topic: Prepare 1x ½ MS medium agar

Purpose: To prepare agar for seed vernalization

Procedure:

1. The following protocol is for 250ml of 1x agar.

½ MS medium powder (w/o sucrose)	0.575g
MES buffer	0.125g
Distilled water	250ml
Sucrose(1%)	2.5g
Agar powder(1%)	2.5g

2. Add the ½ MS medium powder, MES buffer, Distilled water, and Sucrose into a 500ml flask.
3. Use a hot plate to stir (with a magnetic stirring rod) until no sediments are left in the solution.
4. Use the pH probe and measure the pH of the solution, make sure the pH is 5.7-5.8.
 - a. If the pH is too low, add 0.1M KOH slowly to increase it.
5. Add Agar powder into the solution.
6. Wrap the neck of the flask with aluminum foil and stick a piece of indicator tape on top.
7. Autoclave for 20 minutes at 127°C.
 - a. Wait for the pressure to reach 0 before retrieving the flask.
8. Prepare 48 eppendorf microtube lids.
9. Add 250µL of agar into each lid.
10. Wait until it solidifies.

Done by: Alex

Supervised by: Gina

Date: 20260712

Topic: Prepare 2x ½ MS medium agar

Purpose: To prepare agar for seed vernalization for later use

Procedure:

1. The following protocol is for 250ml of 2x agar.

½ MS medium powder (w/o sucrose)	1.15g
MES buffer	0.25g
Distilled water	250ml
Sucrose(1%)	5.0g
Agar powder(1%)	5.0g

2. Add the ½ MS medium powder, MES buffer, Distilled water, and Sucrose into a 500ml flask.

3. Use a hot plate to stir (with a magnetic stirring rod) until no sediments are left in the solution.
4. Use the pH probe and measure the pH of the solution, make sure the pH is 5.7-5.8.
 - a. If the pH is too low, add 0.1M KOH slowly to increase it.
5. Add Agar powder into the solution.
6. Wrap the neck of the flask with aluminum foil and stick a piece of indicator tape on top.
7. Autoclave for 20 minutes at 127°C.
8. Wait for the pressure to reach 0 before retrieving the flask.
9. Aliquote 25mL agar into 50mL centrifuge tubes.
10. Storage at room temperature.

Done by: Alex

Supervised by: Gina

Date: 20260712

Title: Prepare ½ MS medium (hydroponic medium)

Purpose: To prepare medium for hydroponics testing

Procedure:

1. The following protocol is for 1000ml of ½ MS medium

½ MS medium powder (w/o sucrose)	2.3g
MES buffer	0.5g
Distilled water	1000ml

2. Add the contents into a 1000ml flask
3. Use a hot plate to stir (with a magnetic stirring rod) until no sediments are left in the solution
4. Use the pH probe and measure the pH of the solution, make sure the pH is 5.7-5.8
5. If the pH is too low, add 0.1M KOH slowly to increase it
6. Wrap the neck of the flask with aluminum foil and stick a piece of indicator tape on top
7. Autoclave for 20 minutes at 127°C
8. Wait for the pressure to reach 0 before retrieving the flask
9. Add 2.2mL of medium into each tube rack.

Done by Naomi

Supervised by: Gina

Date: 20260712

Title: Golden Gate Assembly

Purpose: Assemble multiple DNA fragments into a construct

Procedure:

A: (STK) Module 2-PlepA-RBS0-CcaSm3-B0015

1. Prepare the following:

Vector(1B)	1.3µL(dilute 1/10)
Insert promoter(PlepA)	1.9µL(dilute 1/10)
Insert RBS(RBS0)	1.5µL(dilute 1/10)
Insert CDS(CcaSm3)	1.5µL(dilute 1/10)

Insert Terminator(B0015)	1.3μL (dilute 1/5)
T4 DNA ligase Buffer	2μL
Sterile water	8.5μL
T4 DNA ligase	1μL
Bsal-HFv2	1μL

Total	20μL
-------	------

2. Use thermocycler to run
 37°C-5 min
 16°C-5 min 30 cycles
 37°C-10 min
 80°C-20 min

B: (STK) Module 4-Pveg-EEBD-RBS0-CSN::ACCD-B0015

1. Prepare the following:

Vector(Exp1)	1.4μL(dilute 1/5)
Insert promoter(Pveg)	1.0μL(dilute 1/5)
Insert RBS(EEBD-RBS0)	1.1μL(dilute 1/10)
Insert CDS(CSN::Accd)	0.6μL(dilute 1/10)
Insert Terminator(B0015)	1.3μL (dilute 1/5)
T4 DNA ligase Buffer	2μL
Sterile water	10.58μL
T4 DNA ligase	1μL
Bsal-HFv2	1μL

Total	20μL
-------	------

2. Use thermocycler to run

37°C-5 min
 16°C-5 min 30 cycles
 37°C-10 min
 80°C-20 min

C: (STK) Module 4-PcpcG2-veg-EEBD-RBS0-csN::ACCD-B0015

Vector(1D)	1.5μL(dilute 1/5)
Insert promoter(Pveg)	1.2μL(dilute 1/50)
Insert RBS(PcpcG2-veg)	1.0μL(dilute 1/5)
Insert CDS(Csn:Accd)	0.6μL(dilute 1/10)
Insert Terminator(B0015)	1.3μL (dilute 1/5)
T4 DNA ligase Buffer	2μL

Sterile water	10.4μL
T4 DNA ligase	1μL
Bsal-HFv2	1μL

Total	20μL
-------	------

2. Use thermocycler to run

37°C-5 min

16°C-5 min 30 cycles

37°C-10 min

80°C-20 min

D: (STK) Module 4-Pveg-EEBD-RBS0-CSN::LEA-B0015

Vector(Exp1)	1.4μL(dilute 1/5)
Insert promoter(Pveg)	1.1μL(dilute 1/10)
Insert RBS(EEBD-RBS0)	1.0μL(dilute 1/5)
Insert CDS(CSN::LEA)	1.1μL(dilute 1/10)
Insert Terminator(B0015)	1.3μL(dilute 1/5)
T4 DNA ligase Buffer	2μL
Sterile water	10.0μL
T4 DNA ligase	1μL
Bsal-HFv2	1μL

Total	20μL
-------	------

2. Use thermocycler to run

37°C-5 min

16°C-5 min 30 cycles

37°C-10 min

80°C-20 min

E: (STK) Module 4-PcpcG2-veg-EEBD-RBS0-CSN::LEA-B0015

Vector(1D)	1.5μL(dilute 1/10)
Insert promoter(PcpeG2-veg)	1.2μL(dilute 1/50)
Insert RBS(EEBD-RBS0)	1.0μL(dilute 1/5)
Insert CDS(Csn:LEA)	1.1μL(dilute 1/10)
Insert Terminator(B0015)	1.3μL(dilute 1/5)
T4 DNA ligase Buffer	2μL
Sterile water	10.05μL
T4 DNA ligase	1μL
Bsal-HFv2	1μL

Total 20μL

2. Use thermocycler to run

37°C-5 min

16°C-5 min 30 cycles

37°C-10 min

80°C-20 min

F: Pveg-EEBD-RBS0-GanA-Gala-B0015

Vector(Shuttle vector) (pSTK-Exp-1sfGFP)	1.5μL(dilute 1/5)
Insert promoter(Pveg)	1.2μL(dilute 1/50)
Insert RBS(EEBD-RBS0)	1.0μL(dilute 1/5)
Insert CDS()	0.6μL(dilute 1/10)
Insert Terminator(B0015)	1.3μL (dilute 1/5)
T4 DNA ligase Buffer	2μL
Sterile water	10.4μL
T4 DNA ligase	1μL
Bsal-HFv2	1μL

Total 20μL

2. Use thermocycler to run

37°C-5 min

16°C-5 min 30 cycles

37°C-10 min

80°C-20 min

Done by Abby K

Supervised by Chars

Date: 20260712

Title: Fresh weight of Arabidopsis

Purpose: To measure stress's impact to the plant's fresh weight

Procedure:

Control-1(contains 3 plants) 0.0154g

Control-2 // 0.0111g

Control-3 // 0.014g

1mM trehalose // -1 0.0095g

1mM trehalose // -2 0.0014g

1mM trehalose // -3 0.0087g

5mM trehalose // -1 0.0101g

5mM trehalose // -2 0.1050g

5mM trehalose // -3 0.0112g

10mM trehalose // -1 0.0114g

10mM trehalose // -2 0.0134g

10mM trehalose // -3 0.0024g

150mM NaCl // -1 0.0090g

150mM NaCl // -2 0.0013g

150mM NaCl // -3 0.0024g

150mM NaCl+1mM trehalose // -1 0.0020g

150mM NaCl+1mM trehalose // -2 0.0015g

150mM NaCl+1mM trehalose // -3 0.0025g

150mM NaCl+5mM trehalose // -1 0.0029g

150mM NaCl+5mM trehalose // -2 0.0023g

150mM NaCl+5mM trehalose // -3 0.0035g

150mM NaCl+10mM trehalose // -1 0.0040g

150mM NaCl+10mM trehalose // -2 0.0036g

150mM NaCl+10mM trehalose // -3 0.0022g

250mM NaCl (contains 7 plants) 0.0013g

250mM NaCl+1mM trehalose // -1 0.0008g

250mM NaCl+1mM trehalose // -2 0.0009g

250mM NaCl+1mM trehalose // -3 0.0014g

250mM NaCl+5mM trehalose // -1 0.0008g

250mM NaCl+5mM trehalose // -2 0.0009g

250mM NaCl+5mM trehalose // -3 0.0011g

250mM NaCl+10mM trehalose // -1 0.0008g

250mM NaCl+10mM trehalose // -2 0.0004g

250mM NaCl+10mM trehaslose // -3 0.0008g

Done by Sophie H, Sophie Yeh, Naomi
Supervised by Gina

2026/7/13 

Date: 20260713

Title: Chlorophyll Extraction

Purpose: To isolate the green photosynthetic pigment from arabidopsis

Procedure:

1. Take out the microtubes with arabidopsis from -80°C (34 microtubes)
2. Dilute acetone to 80%, add 200µL to each tube
3. Use a crusher the plants in the acetone, crush until no visible green plants/leaves
4. Add 800µL of acetone to the tubes
5. Vortex and store in a dark place
6. Centrifuge at 7000 g for 30sec
7. Add the liquid to the circus 96 well, 200µL in each well. 3 wells for each sample

	1	2	3	4	5	6	7	8	9	10	11	12
A	Con trol 1	1T 1	5T 1	10T 1	150 NaC 1	150 NaC 1T 1	150 NaC 5T 1	150 NaC 10T 1	250 NaC 1	250 NaC 1T 1	250 NaC 5T 1	250 NaC 10T 1
B	Con trol 2	1T 2	5T 2	10T 2	150 NaC 2	150 NaC 1T 2	150 NaC 5T 2	150 NaC 10T 2		250 NaC 1T 2	250 NaC 5T 2	250 NaC 10T 2
C	Con trol 3	1T 3	5T 3	10T 3	150 NaC 3	150 NaC 1T 3	150 NaC 5T 3	150 NaC 10T 3		250 NaC 1T 3	250 NaC 5T 3	250 NaC 10T 3

8. ELISA reader, insert 96 well, run the results through the equation below, make sure the select the correct wavelength (A663 and A645)

*Total chlorophyll content equation (mg/g): $(20.2 \cdot OD_{645} + 8.02 \cdot OD_{663}) / V$ (1000*fw),
V=1ml

*250 NaCl 1T 3 mly use 200µL acetone

Done by Abby T, Sophie C, Audrey H, Abby K.

Supervised by Gabriel

Date: 20260713

Title: Gel Electrophoresis

Purpose: To run the gel to verify the length of amplified DNA sequence

Procedure:

1. Load the samples in the following order:(A: pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015, w1 promega ligase; a: pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015, w1 new ligase)

Gel 1:

Done by Sophie H.

Date: 20260713

Title: Transformation

Purpose: Transform bacterial liquid on agar plate

Procedure:

B:(STK) Module 4-Pveg-EEBD-RBS0-csn:ACCD-B0015

D:(STK) Module 4-Pveg-EEBD-RBS)-csn:LEA- B0015

F:Pveg-EEBD-RBS0-GanA-Gala-B0015

Z: SamyQ-TAT-HA-AmilCP

a: pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015

1. Turn on the water bath and set it to 42°C
2. Prepare an ice bucket
3. Thaw the competent cells on ice until 50% thawed
4. Add 15 µl of plasmid to 100 µl of competent cell(do not pipette)
5. Place the mixture on ice for 20 minutes
6. Place the tubes in a 42°C at water bath for 45 seconds-1 minute
7. Place tubes back on ice for 3 minutes to recover
8. Add 1 mL of soc medium under a bunsen burner
9. Incubate at 37°C for 1 hour
10. Centrifuge for 5 minutes at 5000 g
11. Remove 1 mL of supernatant
12. Resuspend the cell pellet in the remaining supernatant
13. Spread the mixture on LBA plate
14. Incubate at 37°C for 16 hours

Done by Chloe, Abby K, Alex, Abby T, Sophie C, Olivia

Date: 20260713

Title: TA cloning

Purpose: Insert PCR-amplified DNA fragments into a plasmid vector without relying on restriction enzymes

Procedure:

1. **Prepare the following reaction:**

Z: samyQ-TAT-HA-AmilCP(277ng/µL; 892bp)

Reagents	Standard Rxn
2x rapid ligation buffer, T4 DNA ligase	5µL

pGEM Teasy vector	1μL
pSTK SamyQ-TAT-HA-AmilCP(Z)	1.6μL
T4 DNA ligase	1μL
Dd H ₂ O	1.4μLL
Total:	10μL

Done by: Chloe

2026/7/14 

Date: 20260714

Title: Colony PCR

Purpose: To amplify bacterial sequences after a transformation step in molecular cloning

Procedure: (A: pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015)

1. Prepare the following reaction (master mix) x24 reaction and dispense solution into PCR tubes

	1 reaction	26 reactions
Fw Premier (pSTK-seq-fw, 10mM)	1µL	26 µL
Rv Premier (STK-seq-rv, 10mM)	1µL	26 µL
Qamp master mix	10µL	260 µL
ddH ₂ O	8µL	208 µL
Total	20µL	520 µL

2. Circle 24 colonies from the plate
3. Pick up colonies from the plate with tip and dispense on PCR tube respectfully
4. Take out tips and close the lid after incubating in room temperature for 5 mins(can also shake the tubes gently)
5. Place into PCR machine for thermal cycle

Cycle:

95°C	5min
95°C	30sec
55°C	30sec
72°C	2.5min
72°C	5min

Done by Sarah, Sophie C, Sophie H.

Supervised by Sophie C.

Date: 20260714

Title: Colony PCR

Purpose: To amplify bacterial sequences after a transformation step in molecular cloning

Procedure:

1. Prepare the following reaction (master mix) x24 reactions and dispense solution into PCR tubes

	1 reaction	26 reactions
Fw Premier (Exp-seq-fw)	1 μ L	26 μ L
Rv Premier (Exp-seq-rv)	1 μ L	26 μ L
Qamp master mix	10 μ L	260 μ L
ddH ₂ O	8 μ L	208 μ L
Total	20 μ L	520 μ L

A: Pveg-EEBD-RBS0-Csn:AccD-B0015

B: Pveg-EEBD-RBS0-GanA-Gala-B0015

C: Pveg-EEBD-RBS0-Csn:LEA-B0015

2. Circle colonies from each plate (A:3 B:18 C: 13)
3. Pick up colonies from each plate with tip and dispense on PCR tubes respectfully
4. Take out plate and close the lid after incubating in room temperature for 5min (or gently shake the PCR tubes horizontally)
5. Place into PCR machine for thermal cycle

95°C	5min
95°C	30sec
55°C	30sec
72°C	2.5min
72°C	5min

Done by Alex, Ryan Wei Supervised by Chloe

Date: 20260714

Title: Prepare 0.8% Gel for Gel Electrophoresis (for 2)

Purpose: Make gel for electrophoresis

Procedure

1. Assemble tray and mold, prepare comb
2. Weigh 0.64g agarose powder and 80 mL 1XTAE buffer
3. Add 1XTAE and powder to flask
4. Cover flask with plastic wrap and puncture holes
5. Microwave for 30 seconds and agitate it. Repeat this step until all particles dissolved
6. Wait until the gel stops smoking, add 2 μ L safety dye and agitate it.
7. Pour the melted gel in the tray and assemble the comb with
8. Leave to solidify in room temperature

Done by: Sarah, Sophie C, Sophie H

Date: 20260714

Title: LB Plate (LBA and LBA with glycerol)

Purpose: prepare agar plate

Procedure

1. Formula as below:
 - a. Broth powder: 6.25g + (Glycerol 2mL)
 - b. Agar powder: 3.75g
 - c. ddH₂O: 250 mL
2. Add one magnetic stirring bar to the flask
3. Place the flask containing the mixture on the hot plate, and stir until no sediments
4. Remove the stirring bar from the flask
5. Cover the flask with aluminum foil and apply an indicator tape to ensure the mixture is successfully autoclaved
6. Autoclave for 20 minutes at 127°C.
7. Place the new petri dishes in BSC and turn on the UVC
8. Set up water bath at 55°C
9. Spray the working surface of the BSC and equipments with 70% ethanol
10. Add 500 μ L of antibiotics (ampicillin) into each corresponding flask
11. Use serological pipette to add 10-15 mL of melting LB agar into the petri dish, avoid overfilling (pipette first to mix the antibiotics)
12. Arrange the petri dishes in stacks within the cabinet
13. Let the agar plates dry for at least 1 hour in the BSC

Done by: Sophie C, Sarah, Naomi

Date: 20260714

Title: Gel electrophoresis

Purpose: Check base pair of colony pcr product

Procedure

1. Load the samples in the following:

a: pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015

A: Pveg-EEBD-RBS0-Csn:ACCD-B0015

B: Pveg-EEBD-RBS0-GanA-Gala-B0015

C: Pveg-EEBD-RBS0-Csn:LEA-B0015

Gel 1: a1 a2 a3 a4 a5 a6 a7 a8 a9 a10 a11 a12 M a13 a14 a15 a16 a17
5 5 5 5 5 5 5 5 5 5 5 5 5 5 5 5 5 5

a18 a19 a20 a21 a22 a23 a24
5 5 5 5 5 5 5


Gel 2: A1 A2 A3 A4 B1 B2 B3 B4 B5 B6 B7 C1 M C2 C3 C4 C5 C6 C7
5 5 5 5 5 5 5 5 5 5 5 5 5 5 5 5 5 5

C8 C9 C10 C11 C12 C13
5 5 5 5 5 5

2. Resolve at 100V for 30 minutes


Result:

20260714 Colony PCR A

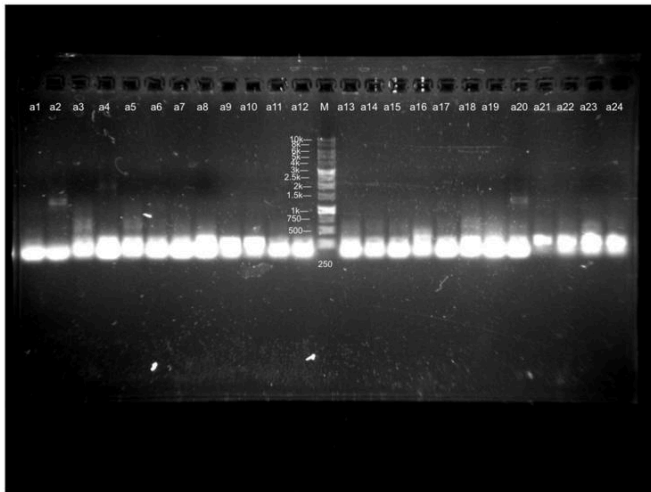

M: 1kb DNA ladder
a1-a24:
pSTK-1A-Pveg-RBS0-ho1::pcyA-B0015
Forward Primer (pSTK-seq-fw)
Reverse Primer (pSTK-seq-rw)
bp: 2231

20260714 Colony PCR B

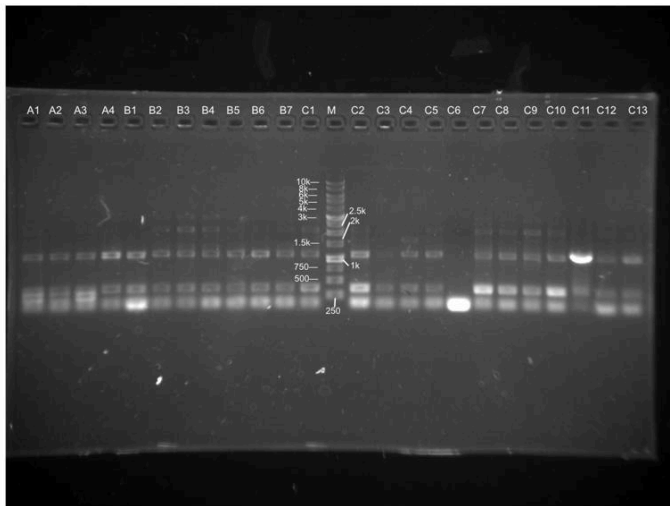

M: 1kb DNA ladder
A: Pveg-EEBD-RBS0-Csn:AccD-B0015
Forward Primer (Exp-seq-fw)
Reverse Primer (Exp-seq-rw)
bp: 1999
B: Pveg-EEBD-RBS0-GanA-Gala-B0015
Forward Primer (Exp-seq-fw)
Reverse Primer (Exp-seq-rw)
bp: 2993
C: Pveg-EEBD-RBS0-Csn:LEA-B0015
Forward Primer (Exp-seq-fw)
Reverse Primer (Exp-seq-rw)
bp: 1262

Done by: Sophie C, Alex, Chloe, Audrey C, Ryan W.

Date: 20260714

Title: Liquid Culture

Purpose: Transfer colony to medium

Procedure

1. Label 7 centrifuge tubes and 8 liquid culture tubes
2. Turn on Bunsen Burner
3. Add 3 μ L of ampicillin and 3 mL of TB medium into every tube
4. Dip a pipette tip into the colonies
5. Disperse tips in their respective tubes
6. Place the tubes in the incubator at a slanted angle

A1: pSTK-1A-Pveg-RBS0-Ho1:PcyA-B0015

B2, B3, B4: pSTK-Exp1-Pveg-EEBD-RBS0-GanA-Gala-B0015

C1, C2, C11: pSTK-Exp1-Pveg-EEBD-RBS0-Csn:LEA-B0015

D2: pSTK-Exp1-Pveg-EEBD-RBS0-SamyQ-TAT-HA-Csn:AmilCP-B0015

Done by: Chars

2026/7/15 

Date: 20260715

Title: Colony PCR

Purpose: Verify DNA insert of Ho1:pcyA for light circuit experiment

Procedure:

1. Prepare the following master mix x 22 reactions and dispense solutions into pcr tubes

	1 reaction	22 reactions
Forward primer(biobrick prefix)	1 μ l	22 μ l
Reverse primer(pcy A reverse primer: ATA TAT GAA GAC ATG AGG GGT CTC ATG TCT TAC TGG ATG ACA TCA AA)	1 μ l	22 μ l
Q-amp master mix	10 μ l	220 μ l
dd water	8 μ l	176 μ l

2. Circle 20 colonies from plate: a) pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015
3. Pick up 20 colonies from the agar plate with tip and streak the labelled grid
4. Dispense tips with colonies into PCR strips and wait 5 mins
5. Take out tips and close lids
6. Place in PCR machine for thermal cycle

95°C	2 min	1 cycle
95°C	30 sec	30 cycles
65°C	40 sec	
72°C	2 min(bp 1800)	
72°C	1 min	1 cycle

Done by Ryan Y, Joshua

Supervised by Chars

Date: 20260715

Title: Enzyme Digestion

Purpose: verify the effectiveness of Bsal-HFv2 RE

Procedure:

Date: 20260715

Title: Make 0.8% Gel(x2)

Purpose: make agarose gel for electrophoresis

Procedure:

1. Assemble tray and mold, prepare comb
2. Weigh 0.48g agarose powder and 60ml 1xTAE solution
3. Add 1xTAE and powder to 250ml conical flask
4. Cover flask with plastic wrap and puncture holes
5. Microwave for 30 seconds and agitate it. Repeat this step until all particles dissolved
6. Wait until the gel stops smoking, add 2 μ l safety dye and agitate it
7. Pour the melted gel in the tray and assemble the comb
8. Leave to solidify in room temperature

Done by Abby T/K

Date: 20260715

Title: Gel electrophoresis

Purpose: verify colony PCR products

Procedure:


1. Loading the samples in following layout:

<u>M</u>	<u>1</u>	<u>2</u>	<u>3</u>	<u>4</u>	<u>5</u>	<u>6</u>	<u>7</u>	<u>8</u>	<u>9</u>	<u>10</u>	<u>11</u>	<u>12</u>	<u>13</u>	<u>14</u>	<u>15</u>	<u>16</u>	<u>17</u>	<u>18</u>	<u>19</u>	<u>20</u>
5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5

M: 1 kb ladder

1-20: pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015

2. Resolve at 100V for 30 mins

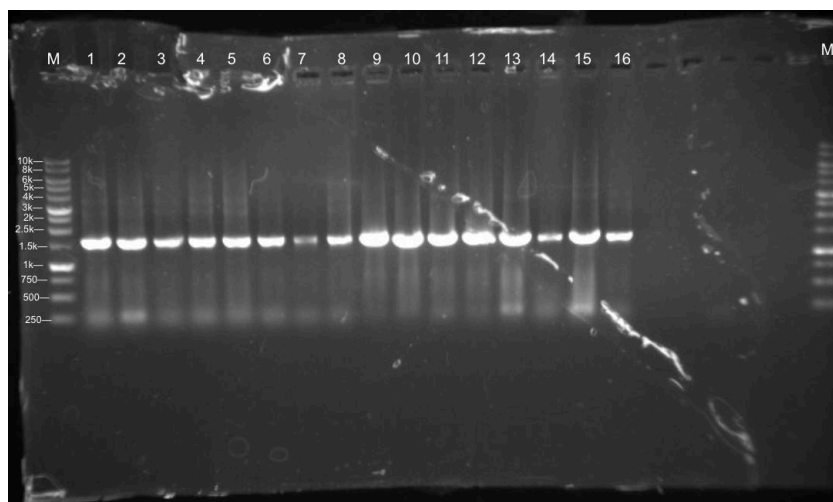

Done by: Abby K, Ryan Y

Supervised by: Chars

Date: 20260715

Title: Gel electrophoresis**Purpose: verify enzyme digestion products****Procedure:**

1. Loading the samples in following layout:

<u>M</u>	<u>B2</u>	<u>B3</u>	<u>B4</u>	<u>C1</u>	<u>C2</u>	<u>#4</u>	<u>SamyQ</u>
5	5	5	5	5	5	5	5

M: 1kb ladder


B: pSTK-EXP1 Pveg-EEBD-RBS0-GanA-Gala-B0015

C: pSTK-Exp1 Pveg-EEBD-RBS0-csn:LEA-B0015

#4: pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015

Samy Q: Pveg-EEBD-RBS0-SamyQ-TATHA-csn:AmilCP-B0015

2. Resolve at 100V for 30 mins

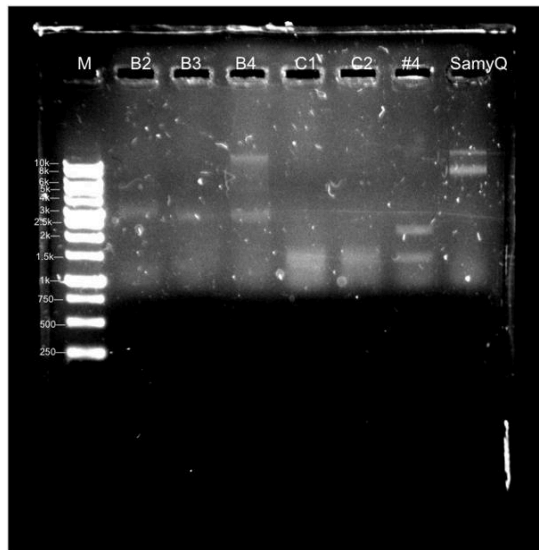Done by: Sophie C, Ryan Y**Supervised by: Chars****Results:**

M: 1kb DNA ladder

B: pSTK-EXP1-Pveg-EEBD-RBS0-GanA-Gala-B0015

C: pSTK-EXP2-Pveg-EEBD-RBS0-Csn::LEA-B0015

#4: pSTK-1A-Pveg-RBS0-ho1::pcyA-B0015

SamyQ: Pveg-EEBD-RBS0-SamyQ-TATHA-Csn::AmilCP-B0015

(failed Enzyme Digestion)

Date: 20260715**Title: Transformation****Purpose: Transfer bacterial liquid on agar plate****Procedure**

A: PlepA-RBS0-ccaSm3-B0015

C: pcpcG2-veg-EEBD-RBS0-csn:ACCD-B0015

E: pcpcG2-veg-EEBD-RBS0-csn:LEA-B0015

1. Turn on waterbath and set it to 42°C and prepare ice bucket
2. Get LBA plate

3. Thaw competent cells on ice until they're 50% thawed and add 15 μ L of plasmid
4. Incubate on ice for 20 min
5. Heat shock for 45 sec - 1 min
6. Recover on ice for 3 min
7. Add 1 mL soc medium & incubate at 37°C for 1 hour
8. Centrifuge at 5000 x g for 5 min and remove 1 mL supernatant
9. Resuspend pellet in remaining supernatant & spread on LBA plate
10. Incubate at 37°C for 16 hours

Done by: Abby K, Chloe

Date: 20260715

Title: Transformation

Purpose: amplification of arrived level 0 DNA parts (plasmids)

Procedure:

- A. pTwist-Amp-High-SamyQ-LEA-v2t2
 - B. pTwist-Amp-High-SamyQ-AmilCP-v2t1
 - C. pTwist-Amp-High-SamyQ-LEA-v2t1
 - D. pTwist-Amp-High-SamyQ-ACCD-v3t1
 - E. pTwist-Amp-High-SamyQ-ACCD-v2t1
 - F. pTwist-Amp-High-SamyQ-AmilCP-v3t1
1. Turn on water bath and set it to 42°C
 2. Get LBA plates and label corresponding DNA parts
 3. Add 15 μ L competent cells and 1 μ L plasmid into 1.5ml microtube
 4. Incubate on ice for 45 sec
 5. Recover on ice for 3 min

Done by Sophie H.

Date: 20260715

Title: Miniprep

Purpose: Collect concentrated bacterial liquid

Procedure:

B2, B3, B4: pSTK-Exp1 Pveg-EEBD-RBS0-GanA-Gala-B0015

C1, C2: pSTK-Exp1 Pveg-EEBD-RBS0-csn:LEA-B0015

#4: pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015

SamyQ: Pveg-EEBD-RBS0-SamyQ-TATHA-csn:AmilCP-B0015

1. Transfer 1.8mL of TB medium to each 2mL microtube
2. Centrifuge at 23444 g for 30 seconds. Discard the supernatant
3. Repeat steps 1-2 until all medium is collected. Remember to empty the medium in the pipette tips too.

4. Add 200 μ l of Mx1 buffer
5. Vortex or pipette to resuspend
6. Add 250 μ l of Mx2 buffer
7. Invert slowly 8-10 times. Incubate at RT for 2 mins
8. Add 350 μ l of Mx3 buffer
9. Invert slowly 4-5 times
10. Centrifuge for 5 mins at 16000 g
11. Transfer ~900 μ l of supernatant into assembled column(until all transferred) using a 200 μ l pipette
12. Centrifuge at 7000 g for 30 seconds and discard flow through
13. Wash the column with 500 μ l of WN buffer. Centrifuge at 7000 g for 30 seconds. Discard the flow through
14. Wash the column with 700 μ l of WS buffer. Centrifuge at 7000 g for 30 seconds. Discard the flow through
15. Centrifuge for 8 mins at 10000 g
16. Place column into new 1.5mL microtube
17. Add 50 μ l of EB buffer to the center of the membrane
18. Incubate for 3 mins at RT
19. Centrifuge at 10000 g for 3 mins

Biodrop results:

	concentration(ng/ μ l)	A260/280
B2	139	1.958
B3	312	2.026
B4	168	1.976
C1	164	1.843
C2	134	1.836
SamyQ	141	1.958
#4	174	2.000

Done by Abigail and Sophie C

2026/7/16 

Date: 20260716

Topic: Prepare LB plates

Purpose: Prepare plates for bacteria growing

Procedure:

1. Formula as below
 - a. LB broth powder: 6.25g
 - b. Agar Powder: 3.75g
 - c. ddH₂O: 250mL
2. Add one magnetic stirring bar to the flask
3. Place the flask containing the mixture on the hot plate and stir until no sediments
4. Remove the stirring bar from the flask
5. Cover the flask with aluminum foil and apply an indicator tape to make sure the mixture is successfully autoclaved
6. Autoclave for 20 minutes at 127°C.
7. Place the new petri dishes in BSC(adjust sash)
8. Set the water bath at 55°C

Done by Sophie C/H, Sarah, Naomi

Date: 20260716

Topic: Prepare ½ MS medium(hydroponics medium)

Purpose: prepare medium for hydroponic testing

Procedure:

1. The following protocol is for 1000mL of ½ MS medium
 - a. ½ MS medium powder(w/o) sucrose: 2.3g
 - b. MES buffer: 0.5g
 - c. Distilled water: 1000mL
2. Add the contents in the 2000mL beaker
3. Use a hot plate to stir(with a magnetic stir bar) until no sediments are left in the solution
4. Use the pH probe and measure the pH of the solution, ensure the pH is 5.7-5.8
 - a. pH before adjustment: 4.1
 - b. pH after adjustment: 5.7
5. If the pH is too low, add 0.1M KOH slowly to increase it
6. Wrap the beaker with aluminum foil and stick a piece of autoclave indicator tape on top
7. Autoclave for 20 minutes at 127°C

Done by Sophie C/H, Sarah, Naomi, Alex

Date: 20260716

Topic: Colony PCR

Purpose: screen for successful cloning from 0715 transformation

Procedure:

A: pSTK-1B-PlepA-RBS0-CcaSm3-B0015 (pcr product bp:1667)

E:pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:LEA-B0015(pcr product bp:1135)

C: pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:ACCD-B0015(pcr product bp:1596)

1. Prepare the following reaction (3 master mix for each plate) and dispense solutions into PCR tubes

A: pSTK-1D-PlepA-RBS0-CcaSm3-B0015 (pcr product bp:1667)

	1 reaction	12 reactions
Forward primer (Biobrick Prefix)	1µL	12µL
CDS primer (CcaSm3 R, ATA TAT GAA GAC ATG AGG GGT CTC ATG TCT TAA GCT CGA GGC AAA TGG T)	1µL	12µL
Q-amp master mix	10µL	120µL
ddH ₂ O	8µL	96 µL

E:pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:LEA-B0015(pcr product bp:1135)

	1 reaction	12 reactions
Forward primer (Biobrick Prefix)	1µL	12µL
CDS Reverse primer (LEA R, ATA TAT GAA GAC ATG AGG GGT CTC ATG TCT TAA TGA TGA TGA TGA TGA T)	1µL	12µL
Q-amp master mix	10µL	120µL
ddH ₂ O	8µL	96 µL

C: pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:ACCD-B0015(pcr product bp:1596)

	1 reaction	12 reactions
Forward primer (Biobrick Prefix)	1µL	12µL
CDS reverse primer (ACCD R, ATA TAT GAA GAC ATG AGG GGT CTC ATG TCT TAA TGA TGA TGA TGA T)	1µL	12µL

Q-amp master mix	10µL	120µL
ddH ₂ O	8µL	96 µL

2. Circle 10 colonies each from plate
3. Pick up 10 colonies from each agar plate with tip and streak the labelled grid
4. Dispense the tips with colonies into the PCR strips and wait 5 minutes
5. Take out the tips and close the lid
6. Place into PCR machine for thermal cycle

95°C	2 min	1 cycle
95°C	30 sec	30 cycles
65°C	40 sec	
72°C	2 min	
72°C	1 min	1 cycle

Done by Abby T/K, Abigail, Mia, Audrey C, Olivia D
Supervised by Chars

Date: 20260716

Topic: Prepare 0.8% Gel for electrophoresis

Purpose: prepare gel for electrophoresis for colony pcr

Procedure:

A: pSTK-1B-PlepA-RBS0-CcaSm3-B0015 (pcr product bp:1667)

E:pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:LEA-B0015(pcr product bp:1135)

C: pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:ACCD-B0015(pcr product bp:1596)

1. Assemble the tray, the mold and the comb(x2)
2. Weigh 0.32g of agarose powder and add to the first flask
3. Add 40 mL of 1x TAE buffer to the flask
4. Repeat step 2 and 3 for another flask
5. Cover the flasks with plastic wrap and puncture holes on the wrap
6. Microwave until all particles are dissolved
7. Wait until the gel stops smoking but is still hot
8. Add 2 µl of safety dye to each gel
9. Pour the gel to the assembled trays and use a tip to squeeze out the bubbles
10. Leave it to solidify

Done by Naomi, Sophie C/H, Sarah C

Date: 20260716

Topic: Gel electrophoresis of colony PCR products

Purpose: verify plasmid sizes

Procedure:

1. Load samples in the following layout:

Gel 1:

<u>M</u>	<u>A1</u>	<u>A2</u>	<u>A3</u>	<u>A4</u>	<u>A5</u>	<u>A6</u>	<u>A7</u>	<u>A8</u>	<u>A9</u>	<u>A10</u>	<u>C1</u>	<u>C2</u>	<u>C3</u>	<u>C4</u>	<u>C5</u>
5	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6

Gel 2:

<u>M</u>	<u>C6</u>	<u>C7</u>	<u>C8</u>	<u>C9</u>	<u>C10</u>	<u>E1</u>	<u>E2</u>	<u>E3</u>	<u>E4</u>	<u>E5</u>	<u>E6</u>	<u>E7</u>	<u>E8</u>	<u>E9</u>	<u>E10</u>
5	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6

2. Resolve at 100V for 30min


Results

Gel 1



M: 1kb DNA ladder
A1-A10:
pSTK-1B-PlepA-RBS0-Ccasm3-B0015
C9-C10:
pSTK-1D-PcpcG2-veg-EEBD-RBS0-csn:ACC
D-B0015
E8-E10:
pSTK-1D-PcpcG2-veg-EEBD-RBS0-csn:LEA
-B0015

Gel 2

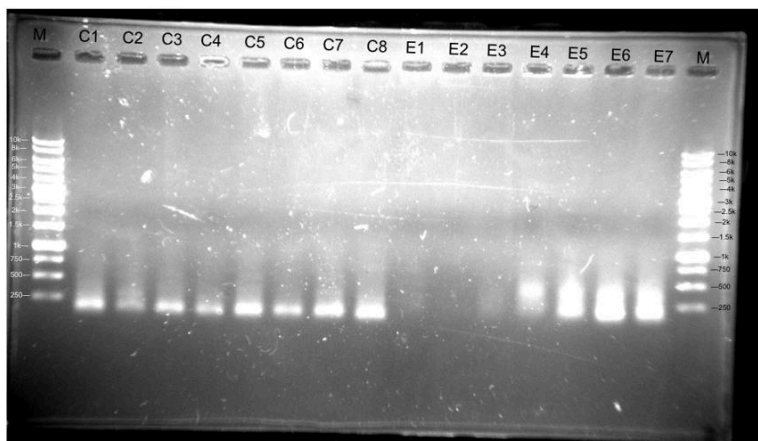

M: 1kb DNA ladder
C1-C8:
pSTK-1D-PcpcG2-veg-EEBD-RBS0-csn:ACC
D-B0015
E1-E7:
pSTK-1D-PcpcG2-veg-EEBD-RBS0-csn:LEA
-B0015

Done by Sarah C, Naomi, Sophie H, Audrey C

Date: 20260716

Topic: Golden Gate assembly

Purpose: redo of golden gate assembly of A, C, E constructs

Procedure:

C: Module 2-PlepA-RBS0-CcaSm3-B0015

I: Module 4- PcpG2-veg-EEBD-RBS0-csn:ACCD-B0015

G: Module 4-PcpG2-veg-EEBD-RBS0-csn:Lea-B0015

C: Module 2-PlepA-RBS0-CcaSm3-B0015

1. Prepare the reaction as follow

Vector (1B) • concentration(ng/μL): 151 • bp: 3120	1.3 μ l (dilute 1/10)
Insert promoter(PlepA) • concentration(ng/μL): 155 • bp: 2241	1.9 μ l (dilute 1/10)
Insert RBS(RBS0) • concentration(ng/μL): 188 • bp: 2113	1.5 μ l (dilute 1/10)
Insert CDS(CcaSm3) • concentration(ng/μL): 128 • bp: 1465	1.5 μ L(dilute 1/10)
Insert Terminator (B0015) • concentration(ng/μL): 107 • Bp: 2213	1.3 μ L (dilute 1/5)
T4 DNA ligase buffer	2 μ L
Sterile water	8.5 μ L
T4 DNA ligase	1 μ L

Bsal-HFv2	1 μ L
Total	20 μ L

2. Use the thermocycler to run the cycle below:

Temperature	Time	Cycles
37 °C (Restriction enzyme digestion)	5 minutes	45 cycles
16 °C (Ligation)	5 minutes	
37 °C (Restriction enzyme digestion)	10 minutes	1 cycle
60 °C (Heat inactivation of Bsal restriction enzyme)	20 minutes	1 cycle

I: Module 4-PcpcG2-veg - EEBD- RBS0 - CSN::ACCD - B0015

1. Prepare the reaction as follow:

Vector (1D) • concentration(ng/μL): 132 • bp: 3120	1.5 μ L (dilute 1/10)
Insert promoter(PcpcG2-veg) • concentration(ng/μL): 65 • bp: 235	1.2 (dilute 1/50)
Insert RBS(EEBD-RBS0) • concentration(ng/μL): 150 • bp: 2392	1.0 μ L (Dilute 1/5)
Insert CDS(CSN::ACCD) • concentration(ng/μL): 142 • bp: 1195	1.0 μ L (dilute 1/10)

Insert Terminator (B0015) - concentration(ng/μL): 107 - bp: 2213	1.3 μ L (dilute 1/5)
T4 DNA ligase buffer	2 μ L
Sterile water	10.4 μ L
T4 DNA ligase	1 μ L
Bsal-HFv2	1 μ L
Total	20 μ L

2. Use the thermocycler to run the cycle below:

Temperature	Time	Cycles
37 °C (Restriction enzyme digestion)	5 minutes	45 cycles
16 °C (Ligation)	5 minutes	
37 °C (Restriction enzyme digestion)	10 minutes	1 cycle
60 °C (Heat inactivation of Bsal restriction enzyme)	20 minutes	1 cycle

G: Module 4-PcpcG2-veg - EEED- RBS0 - CSN::LEA - B0015

1. Prepare the reaction as follow:

Vector (1D) - concentration(ng/μL): 132 - bp: 3120	1.5 μ L (dilute 1/10)
--	---------------------------

Insert promoter(PcpcG2-veg) • concentration(ng/μL): 65 • bp: 235	1.2 (dilute 1/50)
Insert RBS(EEBD-RBS0) • concentration(ng/μL): 150 • bp: 2392	1.0 μL (Dilute 1/5)
Insert CDS(CSN::LEA) 77 ng/μL 634 bp	1.1 μL (dilute 1/20)
Insert Terminator (B0015) • concentration(ng/μL): 107 • bp: 2213	1.3 μL (dilute 1/5)
T4 DNA ligase buffer	2 μL
Sterile water	10.1 μL
T4 DNA ligase	1 μL
BsaI-HFv2	1 μL
Total	20 μL

2. Use the thermocycler to run the cycle below:

Temperature	Time	Cycles
37 °C (Restriction enzyme digestion)	5 minutes	45 cycles
16 °C (Ligation)	5 minutes	
37 °C (Restriction enzyme digestion)	10 minutes	1 cycle

60 °C (Heat inactivation of Bsal restriction enzyme)	20 minutes	1 cycle
--	------------	---------

Done by Sophie C, Abby T/K

2026/7/17 

Date: 20260717

Title: Miniprep

Purpose: To isolate and purify small amounts of plasmid DNA from bacterial cultures done in 20260716

Procedure:

A: pTwist-SamyQ-LEA-v2t2

B: pTwist-SamyQ-AmilCP-v2t1

C: pTwist-SamyQ-LEA-v2t1

D: pTwist-SamyQ-AccD-v3t1

E: pTwist-SamyQ-AccD-v2t1

F: pTwist-SamyQ-AmilCP-v3t1

1. Transfer 1.8mL of TB medium to each 2mL microtube
2. Centrifuge at 23444 g for 30 seconds. Discard the supernatant
3. Repeat steps 1-2 until all medium is collected. Remember to empty the medium in the pipette tips too.
4. Add 200 μ l of Mx1 buffer
5. Vortex or pipette to resuspend
6. Add 250 μ l of Mx2 buffer
7. Invert slowly 8-10 times. Incubate at RT for 2 mins
8. Add 350 μ l of Mx3 buffer
9. Invert slowly 4-5 times
10. Centrifuge for 5 mins at 16000 g
11. Transfer ~900 μ l of supernatant into assembled column(until all transferred) using a 200 μ l pipette
12. Centrifuge at 7000 g for 30 seconds and discard flow through
13. Wash the column with 500 μ l of WN buffer. Centrifuge at 7000 g for 30 seconds. Discard the flow through
14. Wash the column with 700 μ l of WS buffer. Centrifuge at 7000 g for 30 seconds. Discard the flow through
15. Centrifuge for 8 mins at 10000 g
16. Place column into new 1.5mL microtube
17. Add 50 μ l of EB buffer to the center of the membrane
18. Incubate for 3 mins at RT
19. Centrifuge at 10000 g for 3 mins

Results:

BioDrop Results	Concentration	A260/230	A260/280
A			
B			
C			

D			
E			
F			

Date: 20260717

Title: Golden Gate Assembly

Purpose: assemble level 1 genetic circuit

Procedure:

1. Prepare the reaction as follows:

F: Pveg-EEBD-RBS0-GanA-Gala-B0015

Part	Original	Dilution
Vector (shuttle vector): pSTK-EXP-1-sfGFP • concentration(ng/μL): 162 • bp: 6929	0.278015 μL	1.39 μL (dilute 1/10)
Insert promoter(Pveg) • concentration(ng/μL): 280 • bp: 2361	0.109618 μL	1.1 μL (dilute 1/10)
Insert RBS(EEBD-RBS0) • concentration(ng/μL): 150 • bp: 2113	0.207306 μL	2.1 μL(dilute 1/10)
Insert CDS(GanA-Gala) • concentration(ng/μL): 177 • bp: 2080	0.1527683616 μL	1.53 (dilute 1/10) μL
Insert Terminator (B0015) • concentration(ng/μL): 107 • bp: 2213	0.268869 μL	1.34 μL (dilute 1/5)

Spacer (pSTK-OCX-spacer) - concentration(ng/μL): 127 - bp: 2095	0.10722440 μ L	1.07 μ L (dilute 1/10)
Histag (pSTK-OCZ-HisGX) - concentration(ng/μL): 135 - bp: 2402	0.1156185 μ L	1.16 μ L (dilute 1/10)
T4 DNA ligase buffer	2 μ L	2 μ L
Sterile water	μ L	8.54 μ L
T4 DNA ligase	1 μ L	1 μ L
BsaI-HFv2	1 μ L	1 μ L
Total	20 μ L	20 μ L

2. Run the thermal cycle as following:

Temperature	Time	Cycles
37 °C (Restriction enzyme digestion)	5 minutes	45 cycles
16 °C (Ligation)	5 minutes	
37 °C (Restriction enzyme digestion)	10 minutes	1 cycle
60 °C (Heat inactivation of BsaI restriction enzyme)	20 minutes	1 cycle

Done by Abby T, Abby K

Date: 20260717

Title: Transformation

Purpose: Transform golden gate assembly parts from 7/16 on agar plate with antibiotics to verify if the plasmids are successfully taken up by the bacteria through antibiotic selection

Procedure:

C: (STK) M2-PlepA-RBS0-CcaSm3-B0015

I: (STK) M4-PcpcG2-veg-EEBD-RBS0-Csn:ACCD-B0015

G: (STK) M4-PcpcG2-veg-EEBD-RBS0-Csn:Lea-B0015

1. Turn on the water bath, set it to 42°C and prepare the ice.
2. Prepare LBA plates without glycerol
1. Thaw competent cells on ice until they're 50% thawed and add 10 µL of plasmid into 100µL competent cells (don't pipette)
2. Incubate on ice for 20 minutes.
3. Heat shock at 42°C for 45 seconds.
4. Incubate on ice for 3 minutes.
5. Add 1 mL soc medium & incubate at 37°C for 1 hour
6. Centrifuge 5000 g for 5 minutes.
7. Take out 1000 ML and resuspend the cell.
8. Spread on the agar plate with Ampicillin.
9. Place in incubator for 37°C overnight.

Done by Chloe, Abigail

Supervised by Chars

Date: 20260717

Title: Transformation

Purpose: transform single construct and culture bacteria to verify if the plasmids are successfully taken up by the bacteria through antibiotic selection

Procedure:

pSTK-EXP1-Pveg-EEBD-RBS0-SamyQ-TAT-HA-AmilCP-B0015

1. Turn on the water bath, set it to 42°C and prepare the ice.
2. Prepare LBA plates with glycerol
3. Thaw competent cells on ice until they're 100% thawed
4. Take out 1 µL of plasmid into 15µL competent cells (don't pipette)
5. Incubate on ice for 20 minutes.
6. Heat shock at 42°C for 45 seconds.
7. Incubate on ice for 3 minutes.

Done by Abby K, Sophia Yeh

Supervised by Chars

2026/7/18 

Date: 20260718

Title: miniprep

Purpose: plasmid purification

Procedure:

A: pSTK-0A-Pveg

B: pSTK-0B-EEBD

C: pSTK-0D-B0015

1. Transfer 1.8mL of TB medium to each 2mL microtube
2. Centrifuge at 23444 g for 30 seconds. Discard the supernatant
3. Repeat steps 1-2 until all medium is collected. Remember to empty the medium in the pipette tips too.
4. Add 200 μ L of Mx1 buffer
5. Vortex or pipette to resuspend
6. Add 250 μ L of Mx2 buffer
7. Invert slowly 8-10 times. Incubate at RT for 2 mins
8. Add 350 μ L of Mx3 buffer
9. Invert slowly 4-5 times
10. Centrifuge for 5 mins at 16000 g
11. Transfer ~900 μ L of supernatant into assembled column(until all transferred) using a 200 μ L pipette
12. Centrifuge at 7000 g for 30 seconds and discard flow through
13. Wash the column with 500 μ L of WN buffer. Centrifuge at 7000 g for 30 seconds. Discard the flow through
14. Wash the column with 700 μ L of WS buffer. Centrifuge at 7000 g for 30 seconds. Discard the flow through
15. Centrifuge for 8 mins at 10000 g
16. Place column into new 1.5mL microtube
17. Add 50 μ L of EB buffer to the center of the membrane
18. Incubate for 3 mins at RT
19. Centrifuge at 10000 g for 3 mins

Biodrop Results	Concentration ng/ μ L	A260/230	A260/280
A	168	2.012	1.263
B	151	1.936	0.950
C	148	1.922	1.050

Done by Sarah, Naomi, Sophie C.

Date: 20260718

Title: Colony PCR

Purpose: screen for successful cloning from 0715 transformation

Procedure:

1. Prepare the following reaction (3 master mix for each plate) and dispense solutions into PCR tubes

C: pSTK-1B-PlepA-RBS0-CcaSm3-B0015 (pcr product bp:1667)

	1 reaction	10 reactions
Forward primer (Biobrick Prefix)	1µL	10 µL
CDS primer (CcaSm3 R, ATA TAT GAA GAC ATG AGG GGT CTC ATG TCT TAA GCT CGA GGC AAA TGG T)	1µL	10µL
Q-amp master mix	10µL	100µL
ddH ₂ O	8µL	80 µL

G: pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:LEA-B0015 (pcr product bp:1135)

	1 reaction	10 reactions
Forward primer (Biobrick Prefix)	1µL	10 µL
CDS Reverse primer (LEA R, ATA TAT GAA GAC ATG AGG GGT CTC ATG TCT TAA TGA TGA TGA TGA TGA T)	1µL	10µL
Q-amp master mix	10µL	100µL
ddH ₂ O	8µL	80 µL

I. pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:ACCD-B0015 (pcr product bp:1596)

	1 reaction	10 reactions
Forward primer (Biobrick Prefix)	1µL	10 µL
CDS reverse primer (ACCD R, ATA TAT GAA GAC ATG AGG GGT CTC ATG TCT TAA TGA TGA TGA TGA TGA T)	1µL	10µL
Q-amp master mix	10µL	100µL
ddH ₂ O	8µL	80 µL

2. Circle 8 colonies each from plate
3. Pick up 8 colonies from each agar plate with tip and streak the labelled grid
4. Dispense the tips with colonies into the PCR strips and wait 5 minutes
5. Take out the tips and close the lid
6. Place into PCR machine for thermal cycle

95°C	2 min	1 cycle
95°C	30 sec	30 cycles
65°C	40 sec	
72°C	1:45 min	
72°C	1 min	1 cycle

Done by Abby T/K, Audrey H

Date: 20260718

Title: prepare 2 plates of 0.8% agarose gel(25 holes each)

Purpose: for electrophoresis running for 20260717 colony PCR

Procedure:

1. Weight 0.32 g of agarose powder and 40ml of TAE buffer into the flask and mix by shaking (x2)
2. Cover the flask with plastic wrap and poke holes on it by using tip
3. Put it into the microwave and heat until the mixture become clear and with no
4. Remove the plastic wrap and wait until the mixture becomes less hot
5. Add 4µL of safety dye and mix by shaking
6. Pour the gel evenly into the two plates and wait for them to solidify

Done by Sophie C Sophie H

ELN by Sophie H

Date: 20260718

Title: prepare 2 plates of 0.8% agarose gel(25 holes each)

Purpose: for electrophoresis running for 0718 colony PCR

Procedure:

7. Weight 0.32 g of agarose powder and 40 ml of TAE buffer into the flask and mix by shaking (x2)
8. Cover the flask with plastic wrap and poke holes on it by using tip
9. Put it into the microwave and heat until the mixture become clear and with no
10. Remove the plastic wrap and wait until the mixture becomes less hot
11. Add 4µL of safety dye and mix by shaking
12. Pour the gel evenly into the two plates and wait for them to solidify

Done by Audrey H

Supervised by Abby T

Date: 20260718

Title: Gel electrophoresis

Purpose: verify the size of the 0718 colony PCR product

Procedure:

1. Load samples in the following layout

Gel 1:

<u>M</u>	<u>C1</u>	<u>C2</u>	<u>C3</u>	<u>C4</u>	<u>C5</u>	<u>C6</u>	<u>C7</u>	<u>C8</u>	<u>I1</u>	<u>I2</u>	<u>I3</u>	<u>I4</u>	<u>I5</u>	<u>I6</u>	<u>I7</u>	<u>I8</u>
5	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6	6

Gel 2:

<u>M</u>	<u>G1</u>	<u>G2</u>	<u>G3</u>	<u>G4</u>	<u>G5</u>	<u>G6</u>	<u>G7</u>	<u>G8</u>
5	6	6	6	6	6	6	6	6

M: 1kb DNA ladder


C1-C8: pSTK-1B-PlepA-RBS0-CcaSm3-B0015 (pcr product bp:1667)

I1-I8: pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:ACCD-B0015 (pcr product bp:1135)

G1-G8: pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:LEA-B0015 (pcr product bp:1596)

Results:

0718 colony pcr gel 1

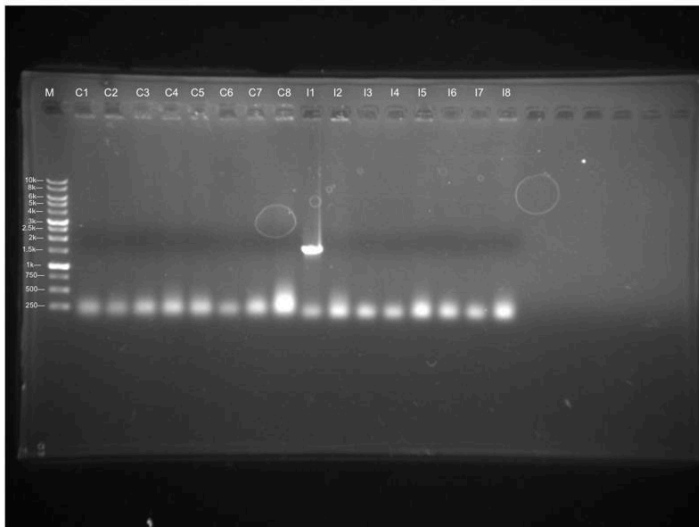

M: 1kb DNA ladder

C1-C8:pSTK-1B-PlepA-Ccasm3-B0015

I1-I8:pSTK-1D-PcpcG2-veg-csn:ACCD-B0015

0718 colony pcr gel 2



M: 1kb DNA ladder
G1-G8:pSTK-1D-PcpcG2-veg-csn:LEA-B0015

Done by Sophie H.

Date: 20260718

Title: Refilling LBA plates

Purpose: prepare agar plates for culturing bacteria

Procedure:

1. Formula as below:
 - a. Broth powder: 6.25g
 - b. Agar powder: 3.75g
 - c. ddH₂O: 250 mL
2. Add one magnetic stirring bar to the flask
3. Place the flask containing the mixture on the hot plate, and stir until no sediments
4. Remove the stirring bar from the flask
5. Cover the flask with aluminum foil and apply an indicator tape to ensure the mixture is successfully autoclaved
6. Autoclave for 20 minutes at 127°C.
7. Place the new petri dishes in BSC and turn on the UVC
8. Set up water bath at 55°C
9. Spray the working surface of the BSC and equipments with 70% ethanol
10. Add 500µL of antibiotics (ampicillin) into each corresponding flask
11. Use serological pipette to add 10-15 mL of melting LB agar into the petri dish, avoid overfilling (pipette first to mix the antibiotics)
12. Arrange the petri dishes in stacks within the cabinet
13. Let the agar plates dry for at least 1 hour in the BSC

Done by Sophie C, Naomi, Sophie H

2026/7/19



Date: 20260719

Title: Colony PCR (redo 0718 colony PCR)

Purpose: Transform and cultivate engineered bacteria on agar plate

Procedure:

1. Prepare the following reaction (3 master mix for each plate) and dispense solutions into PCR tubes

C: pSTK-1B-PlepA-RBS0-CcaSm3-B0015 (pcr product bp:1667)

	1 reaction	10 reactions
Forward primer (Biobrick Prefix)	1µL	10 µL
CDS primer (CcaSm3 R, ATA TAT GAA GAC ATG AGG GGT CTC ATG TCT TAA GCT CGA GGC AAA TGG T)	1µL	10µL
Q-amp master mix	10µL	100µL
ddH ₂ O	8µL	80 µL

G: pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:LEA-B0015 (pcr product bp:1135)

	1 reaction	10 reactions
Forward primer (Biobrick Prefix)	1µL	10 µL
CDS Reverse primer (LEA R, ATA TAT GAA GAC ATG AGG GGT CTC ATG TCT TAA TGA TGA TGA TGA TGA T)	1µL	10µL
Q-amp master mix	10µL	100µL
ddH ₂ O	8µL	80 µL

I. pSTK-1D-PcpcG2-veg-Pveg-EEBD-RBS0-Csn:ACCD-B0015 (pcr product bp:1596)

	1 reaction	10 reactions
Forward primer (Biobrick Prefix)	1µL	10 µL
CDS reverse primer (ACCD R, ATA TAT GAA GAC ATG AGG GGT CTC ATG TCT	1µL	10µL

TAA TGA TGA TGA TGA TGA T)		
Q-amp master mix	10µL	100µL
ddH ₂ O	8µL	80 µL

1. Circle 8 colonies each from plate
2. Pick up 8 colonies from each agar plate with tip and streak the labelled grid
3. Dispense the tips with colonies into the PCR strips and wait 5 minutes
4. Take out the tips and close the lid
5. Place into PCR machine for thermal cycle

95°C	2 min	1 cycle
95°C	30 sec	30 cycles
65°C	40 sec	
72°C	1:45 min	
72°C	1 min	1 cycle

Done by Sophie H.

Date: 20260719

Title: Gel electrophoresis

Purpose: verify the sizes of the colony PCR product

Procedure:

1. Load samples in the following layout

Gel 1:

M C1 C2 C3 C4 C5 C6 C7 C8 I1 I2 I3 I4 I5 I6 I7 I8
5 6 6 6 6 6 6 6 6 6 6 6 6 6 6 6 6

Gel 2:

M G1 G2 G3 G4 G5 G6 G7 G8
5 6 6 6 6 6 6 6 6

M: 1kb DNA ladder


C1-C8: pSTK-1B-PlepA-RBS0-CcaSm3-B0015 (pcr product bp:1667)

I1-I8: pSTK-1D-PcpcG2-veg-Pveg-EEBD-RBS0-Csn:ACCD-B0015 (pcr product bp:1135)

G1-G8: pSTK-1D-PcpcG2-veg-EEBD-RBS0-Csn:LEA-B0015 (pcr product bp:1596)


Results:

20260719 colony pcr C and I

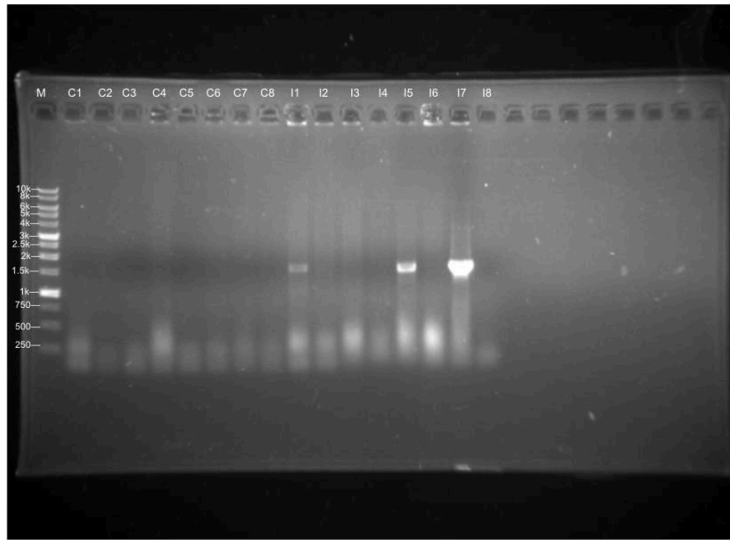

M: 1kb DNA ladder
C1-C8:pSTK-1B-PlepA-Ccasm3-B0015
I1-I8:pSTK-1D-PcpcG2-veg-Csn:ACCD-B0015

20260719 colony pcr G

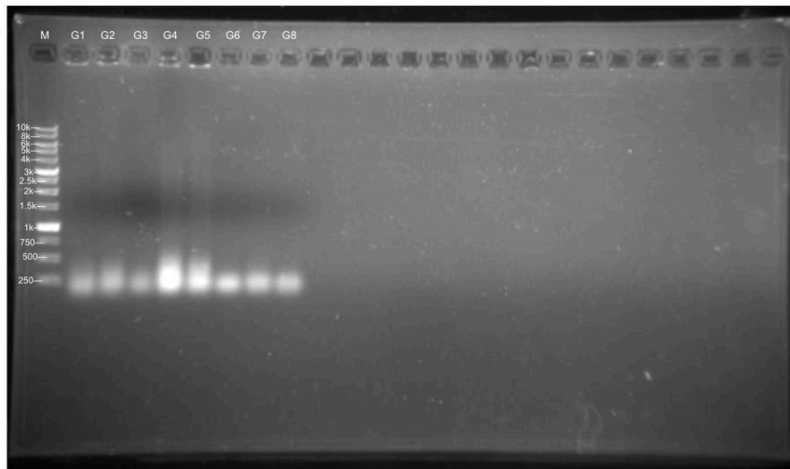

M: 1kb DNA ladder
G1-G8:pSTK-1D-PcpcG2-veg-csn:LEA-B0015

Done by Ryan Y, Olivia L

Date: 20260719

Title: Transformation

Purpose: Transform and cultivate engineered bacteria on agar plate

Procedure:

1. Turn on the water bath and set it to 42°C
2. Prepare an ice bucket
3. Thaw the competent cells on ice until 50% thawed
4. Add 5 µl of plasmid to 100 µl of competent cell(do not pipette)
5. Place the mixture on ice for 20 minutes
6. Place the tubes in a 42°C at water bath for 45 seconds-1 minute

7. Place tubes back on ice for 3 minutes to recover
8. Add 1 mL of soc medium under a bunsen burner
9. Incubate at 37°C for 1 hour
10. Centrifuge for 5 minutes at 5000 g
11. Remove 1 mL of supernatant
12. Resuspend the cell pellet in the remaining supernatant
13. Spread the mixture on LBA plate
14. Incubate at 37°C for 16 hours

Done by Abby T, Audrey H

Date: 20260719

Title: Plant seeds for Exp set 6 and exp set 7 vernalization

Purpose: plant 4 neutral plates (27 seeds/plate)

Procedure:

1. Prepare 4 12x12cm, ½ MS medium, pH 5.8, 1% square agar plates (neutral)
2. Sanitize your hands and all equipments:
 - a. ½ MS medium square agar plates
 - b. Autoclaved toothpicks
 - c. Sterile water
 - d. Planting grid line paper then place in BSC
3. Open the aluminum foil of the autoclaved toothpicks and take one out
4. Dip the tip of the toothpick in the sterile water
5. Attach seeds one at a time and tip it on neutral plates base on the template at the bottom
6. Seal with tape and label the date of vernalization

Done by: Abby T, Audrey H

Recorded by: Audrey H

Date: 20260719

Title: Plant seeds for set A and vernalization in hydroponic system

Purpose: stabilizing plant growth

Procedure:

1. Prepare 2000mL ½ MS medium using the following protocol
 - a.

Material	volume
½ MS medium powder (w/o sucrose)	4.6g
MES buffer	1.0g
Distilled water	2000mL

- b. pH: 5.72

2. Microwave the centrifuge tubes with 1X ½ MS agar from 7/12 and heat it until it melts (agar is 1% agar)
3. Wait until the agar cools down a little bit
4. Prepare 48 eppendorf lids
5. Add 250 µl agar into each eppendorf lids
6. Wait until the agar solidifies
7. Add 2.2mL of medium into 48 tube rack
8. Put the lids with agar onto the tube rack with medium (hole facing up)
9. Use an autoclaved toothpick to pick up one seed and add the seed into the hole of the lid
10. Repeat step 9 until all lids have seeds
11. Wrap the tube rack with aluminium foil
12. Place the tube rack with aluminium foil into the refrigerator in 4°C
13. Wait 4 days for seeds to vernalize

Done by: Alex

Recorded by: Audrey

2026/7/20 

Date: 20260720

Title: Golden Gate Assembly of plasmid into Lv.0 vector

Purpose: to assemble the following parts to STK level 0 vector

Procedure:

- A. pTwist-SamyQ-LEA-v2t1
- B. pTwist-SamyQ-LEA-v2t2
- C. pTwist-SamyQ-Accd-v2t1
- D. pTwist-SamQ-Accd-v3t1
- E. pTwist-SamyQ-AmilCP-v2t1
- F. pTwist-SamyQ-AmilCP-v3t1

1. Prepare the reaction as follow:

A. SamyQ-LEA-v2t1

Component	Amount	Dilution
Vector STK001_pSTK-0-sfGFP - concentration(ng/μL): 346 - bp: 3090	0.05804913295 μ L	1.16 μ L (dilute 1/20) 0.5 μ L vector 9.5 μ L ddH ₂ O
Insert SamyQ-LEA-v2t1 - concentration(ng/μL): 141 - bp: 2840	0.2618439716 μ L	1.047375886 μ L (dilute 1/4) 1 μ L vector 3 μ L ddH ₂ O
T4 DNA ligase buffer	2 μ L	
Sterile water	13.79 μ L	
T4 DNA ligase	1 μ L	
BsaI-HFv2	1 μ L	
Total	20 μ L	

B. SamyQ-LEA-v2t2

Component	Amount	Dilution
Vector STK001_pSTK-0-sfGFP - concentration(ng/μL): 346 - bp: 3090	0.05804913295 μL	1.16 μL (dilute 1/20) 0.5μL vector 9.5μL ddH ₂ O
Insert SamyQ-LEA-v2t2 - concentration(ng/μL): 111 - bp: 2840	0.333	1.33 μL (dilute 1/4) 1μL insert 3μL ddH ₂ O
T4 DNA ligase buffer	2 μL	
Sterile water	13.51 μL	
T4 DNA ligase	1 μL	
Bsal-HFv2	1 μL	
Total	20 μL	

C. pTwist-SamyQ-Accd-v2t1

Component	Amount	Dilution
Vector STK001_pSTK-0-sfGFP - concentration(ng/μL): 346 - bp: 3090	0.05804913295 μL	1.16 μL (dilute 1/20) 0.5μL vector 9.5μL ddH ₂ O
Insert SamyQ-Accd-v2t1 - concentration(ng/μL): 125 - bp: 3401	0.353704 μL	1.061112 μL (dilute 1/3) 1μL insert 2 μL ddH ₂ O
T4 DNA ligase buffer	2 μL	
Sterile water	13.78μL	

T4 DNA ligase	1 μ L	
Bsal-HFv2	1 μ L	
Total	20 μ L	

D. pTwist-SamQ-Accd-v3t1

Component	Amount	Dilution
Vector STK001_pSTK-0-sfGFP - concentration(ng/μL): 346 - bp: 3090	0.05804913295 μ L	1.16 μ L (dilute 1/20) 0.5 μ L vector 9.5 μ L ddH ₂ O
Insert SamQ-Accd-v3t1 - concentration(ng/μL): 142 - bp: 3401	0.311 μ L	1.24 μ L (dilute 1/4) 1 μ L insert 3 μ L ddH ₂ O
T4 DNA ligase buffer	2 μ L	
Sterile water	13.6 μ L	
T4 DNA ligase	1 μ L	
Bsal-HFv2	1 μ L	
Total	20 μ L	

E. SamyQ-AmilCP-v2t1

Component	Amount	Dilution
Vector STK001_pSTK-0-sfGFP concentration(ng/μL): 346	0.05804913295 μ L	1.16 μ L (dilute 1/20)

bp: 3090		0.5µL vector 9.5µL ddH ₂ O
Insert SamyQ-AmilCP-v2t1 - concentration(ng/µL): 138 - bp: 3029	0.2853405797	1.141362319µL (dilute 1/4) 1µL insert 3µL ddH ₂ O
T4 DNA ligase buffer	2 µL	
Sterile water	13.7µL	
T4 DNA ligase	1 µL	
BsaI-HFv2	1 µL	
Total	20 µL	

F. SamyQ-AmilCP-v3t1

Component	Amount	Dilution
Vector STK001_pSTK-0-sfGFP - concentration(ng/µL): 346 - bp: 3090	0.05804913295 µL	1.16 µL (dilute 1/20) 0.5µL vector 9.5µL ddH ₂ O
Insert SamyQ-AmilCP-v3t1 - concentration(ng/µL): 150 - bp: 3029	0.263 µL	1.052µL (dilute 1/4) 1µL insert 3µL ddH ₂ O
T4 DNA ligase buffer	2 µL	
Sterile water	13.79µL	
T4 DNA ligase	1 µL	

Bsal-HFv2	1 μ L	
Total	20 μ L	

Done by Naomi, Sophie H.

Date: 20260720

Title: Transfer plants from tube rack to Box 1,2,3

Purpose: to assemble the following parts to STK level 0 vector

Procedure:

1. Prepare 2000ml of 10% bleach and add into 2000ml beaker
2. Place the 3D boards into the beaker and use plastic wrap to cover it to prevent the 3D board from floating too high
3. Wait for 10 minutes
4. Use ddH₂O to rinse the 3D boards until there are no smell of the bleach
5. Put one tweezer, 3 tip boxes that are autoclaved, 3 3D printer boards, and the tube rack with 41 plants into BSC
6. Put 1000ml of ½ MS medium from 7/19 into the BSC and add about 300ml medium into each tip box
7. Put the 3D printer boards onto the medium in the tip boxes
8. Use the tweezer to pick up one lid and transfer it to the 3D printer boards in the tip boxes. Remember to choose the lids with plants that are growing well.
9. Take everything out of BSC
10. Label the 3 tip boxes as Box 1, Box 2, Box 3
11. Put Box 1, 2, and 3 into the growth chamber at 22.2%

Done by Alex, Ryan

Recorded by: Audrey H

2026/7/21 

Date: 20260721

Title: Prepare antibiotic - erythromycin

Purpose: to prepare erythromycin for growing colonies after electroporation

Procedure:

1. Add 0.01g of erythromycin with 10mL 100% alcohol
2. Mix until erythromycin completely dissolve in the alcohol
3. Distribute 1000µL to each tubes

Done by Sarah, Naomi, Sophie H.

Date: 20260721

Title: Miniprep of colour highlighted CDS in STK library

Purpose: Collect concentrated bacterial liquid

Procedure:

1.

		Concentration (ng/ML)	A260/230	A260/280
pSTK 052	pSTK-OC-eCFP codon	108	1.200	1.964
pSTK 053	pSTK-OC-GFPmut3b	82	1.155	1.952
pSTK 054	pSTK-OC-scarlett	95	1.532	2.065
pSTK 055	pSTK-OC-ECFP codon	100	1.220	1.961
pSTK 125	PSTK-OC-mcherry	33	2.357	2.063
pSTK 127	pSTK-OC-morange	128	1.422	2.000
pSTK 158	pSTK-OC-mRFP1	108	1.113	1.929

2. Transfer 1/8mL of liquid culture to 2mL microtube
3. Centrifuge at 16000 g for 60 seconds. Discard the supernatant.
4. Repeat steps 1-2 until collected all medium
5. Add 200ML MX1 buffer
6. Resuspend by vortexing or pipetting
7. Add 250 ML MX2 buffer
8. Gently mix by inverting 8-10 times. Incubate at RT for 2 minutes
9. Add 350µL MX3 buffer
10. Invert slowly 4-5 times
11. Centrifuge for 5 minutes at 1600 g
12. Transfer ~800µL of the supernatant into assembled column (till all transferred)
13. Centrifuge at 7000 g for 30 seconds and discard flow through
14. Wash the column with 500 µL of WN buffer. Centrifuge at 7000 g for 30 seconds. discard flow through

15. Wash the column with 700µL of WS buffer. Centrifuge at 7000 g for 30 seconds. Discard flow through.
16. Centrifuge for 8 mins at 10000 g
17. Place column into new 1.5mL microtube
18. Add 50µL EB to the center of the membrane
19. Incubate for 3 mins at RT
20. Centrifuge at 10000 g for 3 minutes.

Done by: Abby K. Sophie C

Recorded by: Audrey H

Date: 20260721

Title: Preparation of LBE plates with glycerol

Purpose: Culture electroporation results

Procedure:

1. Measure 3.75 of agar powder and 6.25g of LB powder
2. Pour them into a flask with 250ML of ddH₂O and 2mL of glycerol
3. Use a string bar to mix until even
4. Wrap with aluminium foil and stick indicator tape on top
5. Autoclave for 20 minutes at 127°C.
6. Add 250mL 1 mg/mL erythromycin
7. Use a pipette gun to fill 10cm circular agar plates, 15mL per plate
8. Set at room temperature under the BSC

Done by: Sophie C, Abby K

Recorded by: AudreyH

Date: 20260721

Title: Transformation of level 0 SamyQ and 0720 GGA of G, H, J, K

Purpose: Level 0 assembly

Procedure:

1.

C	PlepA-RBS0-ccasm3
G	pepcG2veg-csn:LEA
J	pveg-csn:ACCD
H	pveg-csn:LEA

2. Level 0 assembly products:
3. Turn on water bath to 42°C, prepare ice bucket
4. Get LBA with glycerol plates
5. Inaw competent cells until 50% thawed and add 10 µL plasmid

6. Incubate on ice for 20 mins
7. Heat shock at 42°C for 45 seconds
8. Incubate on ice for 3 mins
9. Add 1mL SOC medium and incubate at 37°C for 1 hour
10. Centrifuge 5000 g for 5 mins and take out 1200 µL
11. Resuspend cell, spread on LBA plate and incubate at 57°C overnight

Done by: Chloe and Sophie C

Recorded by: Audrey H

Date: 20260721

Title: Electrophoration

Purpose: insert the plasmid

pSTK-EXP1-Pveg-EEBD-RBS0-SamyQ-TAT-HA-AmilCP-B0015 into B.sub

Procedure:

1. Prewarm SOC recovery medium in 37°C
2. Precold 0.2cm electrophoration cuvette
3. Mix 4µl plasmid (con 2050 ng/µl) with 96µl electrocompetent cell in 1-5mL microfuge
4. Set up electroporator (gene pulsatory cell)
 - a. Go to home page
 - b. When using the 0.2cm cuvette, press 2 or 3 and enter
 - c. Set voltage resistance, voltage, and capacity
 - i. Voltage: 2.3kV
 - ii. Resistance: 300Ω
 - iii. capacitance: 25F
5. Transfer mixture of cell and DNA to cold electrophoration cuvette and tap the suspension to the bottom
6. Place the cuvette into the shock pad and close the lid
7. Wait until beep sound
8. Remove the cuvette quickly from the shock pad and immediately add 1mL SOC medium to the cuvette
 - a. Purification of process should be in 10 seconds
9. Transfer the cell suspension to a polypropylene tube and incubate at 37°C for three hours, shaking at 200rpm
10. Check and record the parameter
11. Centrifuge the tubes at 5000 g for 5 mins
12. Spread the LB plate with erythromycin

Done by: Sarah chou, phoebe

Supervised by: Gabriel, abby K, Sophie C

Recorded by: Sophie H

Date: 20260721

Title: Redone Golden Gate Assembly of L

Purpose: Reassemble GanA-Gala

Procedure:

1. Assemble the following parts:

L: pSTK-EXP-01-Pveg-EEBD-RBS0-GanA-GalA-B0015

	Concentration	Bp length	Original (ML)	Diluted (ML)
Vector (EXP1)	162ng/ML	6929	0.278015	1.39 (1/5)
Promoter(Pveg)	162ng/ML	2361	0.109618	1.1(1/10)
RBS(EEBD-RBS0)	150ng/ML	2113	0.207306	2.1(1/10)
CDS(GanA-GalA)	177ng/ML	2080	0.152768	1.53(1/10)
Terminator (B0015)	107ng/ML	2213	0.268869	1.34(1/5)
Spacer(pSTK-0CX-spacer)	127ng/ML	2095	0.10722	1/07(1/10)
Histag(pSTK-0CS-HisGX)	135ng/ML	2402	0.11561	1.10(1/10)
T4 DNA ligase buffer	X	X	2	2
Sterile water	X	X		8.54
T4 DNA ligase	X	X	1	1
Bsal-HFv2	X	X	1	1
Total	X	X	20	20

Done by: Olivia L

Supervised by: Gabriel

Recorded by: Audrey H

Date: 20260720

Title: Transfer Arabidopsis to salinity plates

Purpose: Start experiment set 5 using seeds from 0719

Procedure:

1. Take the stack of agar plates wrapped in aluminum foil out of the refrigerator
2. Open the BSC
3. Spray alcohol on the agar plates before placing them inside the BSC
4. Use a marker to draw a straight line along a ruler 3cm from the top of the plate using
5. Draw 8 dots at the ends of the roots

6. Place the square plates vertically at a slight angle
7. Choose 8 germinated seeds after vernalization that have similar root lengths or germinated status and align the endpoints of their roots to the line drawn on the plate
8. Seal the plates with breathable tape and label the incubation date with a marker before placing the vertical plates on the stacks in the growth chamber, and slightly tilted towards the back

Done by Audrey H, Abby K

Recorded by: Audrey H

2026/7/22 

Date: 20260722

Title: Colony PCR of C, G, and K

Purpose: To amplify the desired DNA fragments

Procedure:

1. Prepare the master mix using the following reaction:

		[K]	[G]	[C]
	1 rxn	4rxn	2rxn	11rxn
FW primer	1μL	4μL	2μL	11μL
RV primer	1μL	4μL	2μL	11μL
Q-amp master mix	10μL	40μL	20μL	110μL
ddH ₂ O	8μL	32μL	16μL	88μL

K: pExp1-Pveg-SamyQ-TatHA-AmilCP-B.sub

G: (STK) Module4: pcpcG2-veg-EEBD-RBS0-Csn:LEA-B0015

C: (STK) Module 2: PlepA-RBS0-Ccasm3-B0015 (bp:1667)

[K]:FW&RV primer:Exp1

[G]:FW primer: biobrick prefix/ RV primer:LEA

[C]: Fw primer: biobrick prefix/ RV primer:Ccasm3 rv

2. Distribute 20µL per microtube
3. Pick the circled colonies with the tip and streak the gridded plates
4. Dispense the tips into one of the microtubes
5. Repeat steps 3 & 4 until all microtubes are filled
6. Set up a 5min timer and begin after last tip
7. Remove tips and close lids after 5min and centrifuge
8. Put the PCR tubes into PCR machine and set up thermal cycle as below:

95°C	5 min	1 cycle
95°C	20 secs	35 cycles
65°C	30 secs	//
72°C	2 min (bp 1800)	//
72°C	5 min	1 cycle

Done by Sophie C.

ELN by Sophie H.

Date: 20260722

Title: 0.8% gel preparation

Purpose: Make 0.8% gel for electrophoresis

Procedure:

1. Assemble mold, tray, and comb
2. Weight 0.64g of agarose powder and prepare 80mL of 1X TAE buffer
3. Add 1X TAE buffer and Agarose powder into conical flask
4. Cover flask with plastic wrap and use tip to puncture holes
5. Microwave for 30sec than agitate, repeat this step until all particles were dissolved
6. When gel stops smoking but is still hot, add 2 μ L of safety dye and agitate it
7. Pour the melted gel into two trays
8. Leave to solidify in room temperature

Done by Sophie C, Sophie Liu

ELN by Sophie H

Date: 20260722

Title: Electrophoresis

Purpose: Check colony PCR product of C,G,K,A,E,

Procedure:

1. Load the samples in the following order:

Gel1: volume(μ L)

<u>M</u>	<u>C1</u>	<u>C2</u>	<u>C3</u>	<u>C4</u>	<u>C5</u>	<u>C6</u>	<u>C7</u>	<u>C8</u>	<u>C9</u>	<u>C10</u>	<u>G1</u>	<u>K1</u>	<u>K2</u>	<u>K3</u>
4	5	5	5	5	5	5	5	5	5	5	5	5	5	5

Gel2: volume(μ L)

<u>M</u>	<u>A1</u>	<u>A2</u>	<u>A3</u>	<u>E1</u>	<u>E2</u>	<u>E3</u>	<u>B1</u>	<u>B2</u>	<u>B3</u>	<u>F1</u>	<u>F2</u>	<u>F3</u>	<u>D1</u>	<u>D2</u>	<u>D3</u>	<u>H1</u>
5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5
<u>H2</u>	<u>H3</u>															
5	5															

M: 1kb DNA ladder

A: pSTK-SamyQ-LEA-v2t1

E: pSTK-SamyQ-LEA-v2t2

B: pSTK-SamyQ-AccD-v2t1

F: pSTK-SamyQ-AccD-v3t1

D: pSTK-SamyQ-AmilCP-v2t1

H: pSTK-SamyQ-AmilCP-v3t1

C: PlepA-RBS0-ccdsm3-B0015


G: pcpcG2veg-EEBD-RBS0-Csn::LEA-B0015

K: Pveg-EEBD-RBS0-SamyQ-TAT-HA-AmilCP-B0015

2. Resolve the DNA at 10v for 30min


Result:

20260722 colony pcr PlepA-EEBD-RBS0-CcaSm3-B0015

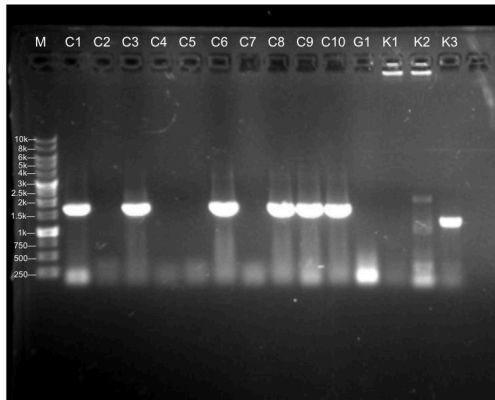

M: 1kb DNA ladder
C: PlepA-EEBD-RBS0-CcaSm3-B0015
G: PcpG2-veg-EEBD-RBS0-csn::LEA-B0015
K: Pveg-EEBD-RBS0-SamyQ-TAT-HA-AmilCP-B0015

20260722 Bsa1 Enzyme digestion

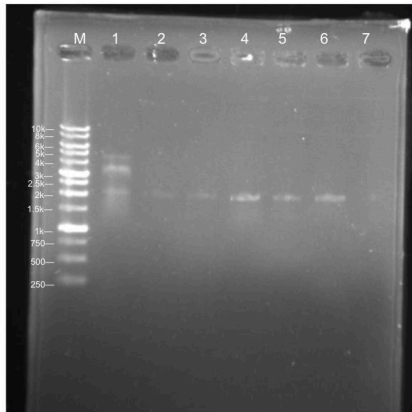

M: 1kb DNA ladder
1: Ho1:pcyA
2: pveg
3: spacer (old enzyme)
4: GanA-Gala
5: HIStag
6: B0015
7: spacer (new enzyme)

Done by Sophie C, Abby K
ELN by Sophie H

Date: 20260722

Title: Transformation

Purpose: To transform the desired bacteria in medium to agar plate

Procedure:

pSTK-0A-Pveg

pSTK-0D-B0015

pSTK-0B-EEBD

pGEM-Teasy-0C-Ho1:pcyA

1. Turn on water bath at 52°C and prepare ice
2. Get LBC with glycerol plates

3. Thaw competent cell until 100% thawed and split them into 3 tubes of 30 and 1 tube of 10
4. Add 2µL plasmid into the tubes of 30 and add 1µL plasmid into the tube of 10
5. Incubate on ice for 20min
6. Heat shock at 42°C for 45sec and then recover on ice for 3min
7. Spread on LBC plates and incubate at 37°C overnight

Done by Chole, Sophie C

ELN by Sophie H

Date: 20260722

Title: Gel Electrophoresis for Enzyme Digestion

Purpose: Check the approximate length of enzyme digestion product

Procedure:

M: 1kb DNA ladder

1: Ho1:pcyA

2: Pveg

3: spacer(old enzyme)

4: GanA-Gala

5: HIStag

6: B0015


7: spacer (new enzyme)

1. Run the gel in the following order:

<u>M</u>	<u>1</u>	<u>2</u>	<u>3</u>	<u>4</u>	<u>5</u>	<u>6</u>	<u>7</u>
5	5	5	5	5	5	5	5

Result:

20260722 Bsa1 Enzyme digestion

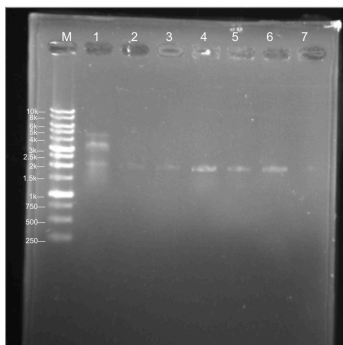

M: 1kb DNA ladder
 1: Ho1:pcyA
 2: pveg
 3: spacer (old enzyme)
 4: GanA-Gala
 5: HIStag
 6: B0015
 7: spacer (new enzyme)

Done by Sophie L

Date: 20260722

Title: Issue of Box B1

Purpose: Summarize and analyze the problem of hydroponic box and find these problem before next round of testing

Problems:

Problem 1: The 3D printer board in Box B1 started to sink into the medium, This is because the medium leads into the 3D printer board. This caused the lids to float on medium, and not on the board.

Problem 2: Two of the lids seemed to be contaminated by red bacteria or fungi

Problem 3: The medium is contaminated and it looks turbid

Problem 4: 8 out of 13 plants are not growing, including the 2 that had been contaminated by red bacteria

Done by plant model team

2026/7/23



Date: 20260723

Title: Liquid culture

Purpose: To transfer colony inside TB medium

Procedure:

1. Label 10 liquid culture tubes
C1, C3, C6, C8, C9, C10: PlepA-RBS0-CcaSm3-B0015
E1: pSTK-SamyQ-LEA-v2t2
D1, D2, D3: pSTK-SamyQ-AmilCP-v3t1
2. Turn on bunsen burner
3. Add 3 μ L of ampicillin and 3mL of TB medium into each tubes
4. Dip a pipette tip into colonies
5. Dispense tips in their respective tubes
6. Place the tubes in incubator at a slanted angle

Done by Sophie C

Supervised by Gabriel

Date: 20260723

Title: Colony PCR of Ho1:pcyA

Purpose: Screen for successful bacterial inserts from the transformation of Ho1:pcyA

Procedure:

1. Prepare the reaction as follows:
A: pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015

	1 rxn	58rxn
Forward Primer (pSTK-seqFW)	1 μ L	58 μ L
Reverse Primer (pcyA R)	1 μ L	58 μ L
Q-amp master mix	10 μ L	580 μ L
ddH ₂ O	8 μ L	464 μ L
total	20 μ L	1160 μ L

2. Circle 24 colonies each from plate A1 and A2, and 8 colonies from plate B

A1 & A2: Pveg-RBS0-Ho1-pcyA-B0015

B: pGEM-Teasy-Ho1-pcyA

3. Pick up 24 colonies from each agar plate with tip and streak the labelled grid 3 times
4. Dispense the tips with colonies into the PCR strips and wait 5 minutes
5. Take out the tips and close the lid
6. Place into PCR machine for thermal cycle
95°C – 2 min (1 cycle)

95°C – 30 sec //
 55°C – 20 sec } 30 cycles
 72°C – 2 min //
 72°C – 5 min (1 cycle)
Done by Abby T, Abby K, Sophie C
Supervised by Chars

Date: 20260723

Title: Golden Gate Assembly of mCherry, mScarlet

Purpose: Assemble parts

Procedure:

1. Assemble the following parts:	original(μL)	diluted(μL)
R: Pveg-RBS0-mCherry-B0015		
Vector: EXP1 (162ng/μL, 6929 bp)	0.278	1.39($\frac{1}{5}$)
I promoter: Pveg (168ng/μL, 2361 bp)	0.182	1.82(1/10)
I RBS: RBS0	0.146	1.46(1/10)
I CDS: mCherry (83ng/μL, 2795bp)	0.437	0.44
I terminator: B0015 (148ng/μL, 2213bp)	0.194	1.94(1/10)
T4 DNA ligase buffer	2	2
Sterile water	last added, add to 20μL	
T4 DNA ligase	1	1
Bsal-HFv2	1	1
total	20	20

Z: Pveg-RBS0-mScarlet-B0015

Vector: EXP1 (162ng/μL, 6929 bp)	0.278	1.39($\frac{1}{5}$)
I promoter: Pveg (168ng/μL, 2361 bp)	0.182	1.82(1/10)
I RBS: RBS0	0.146	1.46(1/10)

I CDS: mScarlet (95ng/μL, 2780bp)	0.38	0.4
I terminator: B0015 (148ng/μL, 2213bp)	0.194	1.94(1/10)
T4 DNA ligase buffer	2	2
Sterile water	last added, add to 20μL	
T4 DNA ligase	1	1
Bsal-HFv2	1	1
total	20	20

J: (STK) Module4-Pveg-EEBD-RBS0-CSN::ACCD

Vector: EXP1 (162ng/μL, 6929 bp)	0.278	1.39(1/5)
I promoter: Pveg (168ng/μL, 2361 bp)	0.182	1.82(1/10)
I RBS: EEBD-RBS0	0.205	1.02(1/5)
I CDS: Csn::ACCD (142ng/μL, 4212bp)	0.385	1.9(1/10)
I terminator: B0015 (148ng/μL, 2213bp)	0.194	1.94(1/10)
T4 DNA ligase buffer	2	2
Sterile water	last added, add to 20μL	
T4 DNA ligase	1	1
Bsal-HFv2	1	1
total	20	20

I:(STK) Module 4-PcpcG2-veg-EEBD-RBS0-CSN::ACCD-B0015

Vector: 1D (132ng/μL, 3120 bp)	0.15	1.5(1/10)
I promoter: Pveg (65ng/μL, 3252bp)	0.6504	0.65
I RBS: RBS0 (151ng/μL, 2392bp)	0.205	1.02(1/5)

I CDS: CSN::ACCD (142ng/μL, 4212bp)	0.385	1.9(%)
I terminator: B0015 (148ng/μL, 2213bp)	0.194	1.94(1/10)
T4 DNA ligase buffer	2	2
Sterile water	last added, add to 20μL	
T4 DNA ligase	1	1
Bsal-HFv2	1	1
total	20	20

Date: 20260723

Title: Prepare 2.0% Agarose Gel

Purpose: prepare gel for electrophoresis

Procedure:

1. Assemble mad, tray, and comb
2. Weight 0.40g of Agarose powder and prepare 20mL of 1X TAE buffer
3. Add 1X TAE and Agarose powder into conical flask
4. Cover the flask with plastic wrap and puncture holes with tip(middle one would be the best)
5. Microwave for 30sec then agitate. Repeat this step until all particles were dissolved
6. When gel stops smoking, add 1μL of safety dye and agitate it
7. Pour the melted gel into tray
8. Leave to solidify in RT

Done by Sophie L

Date: 20260723

Title: Make 1% agarose gel*2

Purpose: Prepare agarose gel for electrophoresis

Procedure:

1. Assemble mad, tray, and comb
2. Weight 0.8g of Agarose powder and prepare 80mL of 1X TAE buffer
3. Add 1X TAE and Agarose powder into conical flask
4. Cover the flask with plastic wrap and puncture holes with tip(middle one would be the best)
5. Microwave for 30sec then agitate. Repeat this step until all particles were dissolved
6. When gel stops smoking, add 1μL of safety dye and agitate it
7. Pour the melted gel into tray
8. Leave to solidify in RT

Done by Abby K, Sophie C

Date: 20260723

Title: Gel electrophoresis

Purpose: Verify the approximate length of basepair of Pveg-RBS0-Ho1:pcyA-B0015

Procedure:

1. Load the samples in the following order: (A1/A2:Pveg-RBS0-Ho1:pcyA-B0015)

Gel1:


M	A 1 1	A 1 2	A 1 3	A 1 4	A 1 5	A 1 6	A 1 7	A 1 8	A 1 9	A 1 1 0	A 1 1 1	A 1 1 2	A 1 1 3	A 1 1 4	A 1 1 5	A 1 1 6	A 1 1 7	A 1 1 8	A 1 1 9
5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5
A 1 2 0	A 1 2 1	A 1 2 2	A 1 2 3	A 1 2 4															
5	5	5	5	5															

Gel2:

M	A 2 1	A 2 2	A 2 3	A 2 4	A 2 5	A 2 6	A 2 7	A 2 8	A 2 9	A 2 1 0	A 2 1 1	A 2 1 2	A 2 1 3	A 2 1 4	A 2 1 5	A 2 1 6	A 2 1 7	A 2 1 8	A 2 1 9
5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5
A 2 2 0	A 2 2 1	A 2 2 2	A 2 2 3	A 2 2 4															
5	5	5	5	5															

2. Resolve the DNA at 100V for 30mins


20260723 Colony PCR Module 1 Gel 2

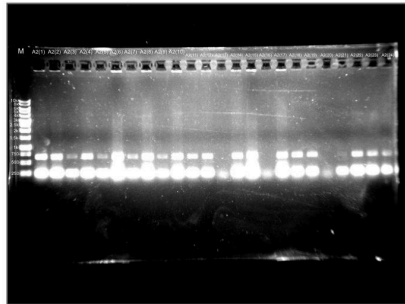

M: 1kb DNA ladder
A2(1)-A2(24):
Pveg-RBS0-Ho1:pcyA-B0015
• forward primer(pSTK-seqFW)
• reverse primer(pcyA R)
• bp: 1871

Result:

20260723 Colony PCR Module 1 Gel 1

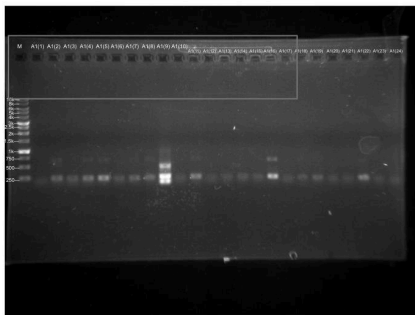

M: 1kb DNA ladder
A1(1)-A1(24):
Pveg-RBS0-Ho1:pcyA-B0015
• forward primer(pSTK-seqFW)
• reverse primer(pcyA R)
• bp: 1871

Primers: pSTK-seqFW and pcyA rv

Done by Sophie C.

2026/7/24 

Date:20260724

Title:PCR for primer quality check

Purpose:Verify primer pairs for colony PCR is work

Procedure:

1. Prepare reaction master mix:

template	1.2	1.2
	pSTK-0-spGFP	pSTK-Exp2-sPEFP
ddH ₂ O	46.8	46.8
Q-AMP mix	60	60

2. Dispense reaction by the following layout:

		1	2
pSTK seq	A	fw fw	fw fw
	B	rv rv	rv rv
	C	fw rv	fw rv
pSTK EXP seq	D	fw fw	fw fw
	E	rv rv	rv rv
	F	fw rv	fw rv
pSTK- seq	G	fw rv	fw rv
pSTK- EXP- seq	H	fw rv	fw rv
		old primer	new primer

3. Thermal cycle:

95°C	5 min	
95°C	20 sec	30 cycles
60°C	20 sec	
72°C	1.5 min	
72°C	5 min	
RT		

Done by: Chars

Recorded by: Audrey H

Date:20260724

Title:Make 1% agarose gel

Purpose: Make gel for electrophoresis

Procedure:

1. Assemble mold, tray, and comb
2. Weight 0.4g agarose powder and prepare 40mL 1XTAE buffer
3. Add into a flask
4. Cover with plastic wrap and poke holes
5. Place into microwave for 30 secs
6. Agitate
7. Put into microwave for 5 sec then agitate
8. Repeat step 7 until no sediments
9. When it stops smoking, add 2 μ L safety dye, agitate
10. Pour into the mold
11. Let solidify at RT

Done by: Abby K, Sophie C

Supervised by: Sophie C

Recorded by: Audrey H

Date:20260724

Title:Redoing 20260723 Colony PCR

Purpose: Amplify plasmid number

Procedure:

1. Prepare the following reaction (master mix) x26 reactions and dispense solution into PCR tubes

	1 reaction	26 reactions
Fw primer (biobrick prefix)	1 μ L	26 μ L

Rv primer (stk reverse)	1µL	26µL
Q-amp master mix	10µL	260µL
ddH ₂ O	8µL	208µL

A1 & A2: Pveg-RBSO-Ho1: pcyA-B0015

2. Distribute 20 µL per microtube
3. Pick the circled colonies with the top and streak the gridded plates
4. Dispense the tops into one of the microtubes
5. Repeat steps 3 and 4 until all microtubes are filled with 1 tip
6. Set up a 5 mins timer and begin after last top
7. Remove tips and close lid after 5 min and centrifuge
8. PCR

95°C	5 mins	1
95°C	20 sec	35 cycles
65°C	20sec	
72°C	2 mins 30 sec	
72°C	5 min	1
25°C	0 min	1

Done by: AbbyK, Sophie C

Recorded by: Audrey H

Date:20260724

Title: Transformation

Purpose: insert plasmids into competent cell

Procedure:

G:(STK) module 4- PcpcG2-veg-EEBD-RBS0-CSN::LEA-B0015

H:(STK) module 4- Pveg-EEBD-RBS0-CSN::LEA-B0015

I: (STK) module 4-PcpcG2-veg-EEBD-RBS0-CSN::ACCD-B0015

Y: Pveg- RBS-mcherry -B0015

J: (STK) Module 4 - Pveg - EEBD - RBSO - CSN: ACCD- B0015

Z: Pveg-RB50-mScarlet -B0015

1. set water bath to 42°C and prepare an ice bucket
2. Prepare LBA plates (6) and place them (with soc medium) in the incubation
3. Thaw the competent cell until 50% melted, then add 10µL of plasmid.
4. Incubate on ice for 20 minutes
5. Heat shock at 42°C for 45 seconds.
6. Recover on ice for 8 minutes.

7. Add 1mL soc medium
8. Incubate for 1 hour at 31°C.
9. centrifuge at 5000 g for 5 minutes
10. Remove 1mL of supernatant then vortex
11. Spread on LBA plates
12. Incubate at 31°C for 16 hours.

Done by: Abby K, Sophie C

Recorded by: Audrey H

Date:20260724

Title:Gel electrophoresis

Purpose: Check base pair of colony PCR product

Procedure:

1. Load the sample in the following sample:

<u>M</u>	<u>A1 1</u>	<u>A1 2</u>	<u>A1 3</u>	<u>A1 4</u>	<u>A1 5</u>	<u>A1 6</u>	<u>A1 7</u>	<u>A1 8</u>	<u>A1 9</u>	<u>A1 10</u>	<u>A1 11</u>	<u>A1 12</u>
5	5	5	5	5	5	5	5	5	5	5	5	5
<u>A2 1</u>	<u>A2 2</u>	<u>A2 3</u>	<u>A2 4</u>	<u>A2 5</u>	<u>A2 6</u>	<u>A2 7</u>	<u>A2 8</u>	<u>A2 9</u>	<u>A2 10</u>	<u>A2 11</u>	<u>A2 12</u>	
5	5	5	5	5	5	5	5	5	5	5	5	5

2. Resolve at 100V for 30 minutes
3. A1 & A2: Pveg-RBS0-Ho1:pcyA-B0015

Done by Olivia Lin, Ryan Y

Supervised by:Chars

Recorded by: Audrey H

Date:20260724

Title:GGA of A, L

Purpose: Assemble parts

Procedure:

A: (STK) Module 1(1A): Pveg-RBS0-Ho1:pcyA-B0015

1. prepare the reaction as follow:

Part	Original Amount	Dilution
Vector (1A)	0.176348 µL	1.76 µL (dilute 1/10)
Insert promoter(Pveg)	0.182 µL	1.82 µL (dilute 1/10)
Insert RBS(RBS0)	0.1461117 µL	1.46 µL (dilute 1/10)

Insert CDS(Ho1::PcyA)	0.0849916 μL	1.7 μL (dilute 1/20)
Insert Terminator (B0015)	0.194 μL	1.94 μL (dilute 1/10)
T4 DNA ligase buffer	2 μL	2 μL
Sterile water		7.32 μL
T4 DNA ligase	1 μL	1 μL
Bsal-HFv2	1 μL	1 μL
Total	20 μL	20 μL

L: Pveg-EEBD-RBS0-GanA-Gala-B0015

Part	Original	Dilution
Vector (shuttle vector): pSTK-EXP-1-sfGFP	0.278015 μL	1.39 μL (dilute $\frac{1}{6}$)
Insert promoter(Pveg)	0.109618 μL	1.1 μL (dilute 1/10)
Insert RBS(EEBD-RBS0)	0.205 μL	2 μL (Dilute 1/10)
Insert CDS(GanA-Gala)	0.30465 μL	1.53 (dilute 1/5) μL
Insert Terminator (B0015)	0.194 μL	1.94 μL (dilute 1/10)
Spacer (pSTK-OCX-spacer)	0.10722440 μL	1.07 μL (dilute 1/10)
Histag (pSTK-OCZ-HisGX)	0.1156185 μL	1.16 μL (dilute 1/10)
T4 DNA ligase buffer	2 μL	2 μL
Sterile water	μL	5.81 μL
T4 DNA ligase	1 μL	1 μL
Bsal-HFv2	1 μL	1 μL
Total	20 μL	20 μL

Done by: Chloe

Supervised by: Chars

Recorded by Audrey H

Date:20260724

Title:Liquid culture

Purpose: Transfer colony to medium

Procedure:

Pveg

RBS0

EEBD

B0015

PGEM-Teasy-Ho1:pcyA

1. Label the five liquid culture tubes
2. Turn on the bunsen burner and add 3mL iB medium. Add 3 μ L ampicillin to Pveg, RBS0, EEBD, B0015 and add 3 μ L chlorophenicol to Ho1:pcyA
3. Dip pipette into colonies and dispense tips into respective tubes
4. Incubate in 25°C overnight

Done by: Chars

Supervised by: Chloe

Recorded by: Audrey H

Date:20260724

Title: Prepare 2% agarose gel

Purpose: Make gel for electrophoresis

Procedure:

1. Assemble mold, tray, and comb
2. Weight 1.2 g of agarose powder and measure 60mL of 1xTAE buffer
3. Add the agarose powder and the 1x TAE buffer into conical flask
4. Microwave for 30 sec then agitate it
5. Repeat this step until no sediments
6. Wait until no smoke, add 3 μ L of safety dye and agitate it as well
7. Pour the melted gel into the tray and use tip to squeeze the air and pop the bubble
8. Wait until the gel solidifies

Done by Abby K

2026/7/25 

Date: 20260725

Title: Planting Arabidopsis in soil

Purpose: To plant the Arabidopsis seed in the soil

Procedure:

1. Soak pots and water trays in 1:100 diluted bleach for 10mins, then rinse it thoroughly with ddH₂O (make sure every corner was rinsed thoroughly until no bleach smell left)
2. Mix peat soil into a 2L flask and autoclave it for 20 minutes at 127°C.
3. it. Fill 3-inch pots until 80% full, lightly press down and thoroughly water from the surface to fully moisturize the soil before transplanting the Arabidopsis seeds.
4. Make a small hole (around one knuckle deep) in the soil using finger (the holes should be spaced evenly to prevent plants interrupting each other from growing)
5. Use tweezers to gently lift the healthy seeds (put for 7 days)
6. Place the root against the sidewall of the hole
7. Use a wash bottle to gently rinse soil over roots so they are fully covered by the soil (must avoid direct contact of the root with fingers or tweezers, support seedlings only with tweezers only from beneath the cotyledons)

Done by Alex

Record by Sophie H.

Date: 20260725

Title: Colony PCR

Purpose: To verify bacterial colonies after a transformation step in molecular cloning

Procedure:

1. Prepare the following master mix

Z: Pveg-RBS0-mScarlet-B0015

Y: Pveg-RBS0-mCherry-B0015

	1 reaction	26 reactions
Forward(Pveg-fw)	1 µl	26 µl
Reverse(EXP-rv)	1 µl	26 µl
Q-amp master mix	10 µl	260 µl
ddH ₂ O	8 µl	208 µl

J: Pveg-EEBD-RBS0-csn::ACCD-B0015

H: Pveg-EEBD-RBS0-csn::LEA-B0015

	1 reaction	26 reactions
Forward primer (EXP fw)	1µL	26 µL

Reverse primer (EXP seq 4 rv)	1µL	26 µL
Q-amp master mix	10µL	260 µL
ddH ₂ O	8µL	208 µL

G: PcpG2-veg-EEBD-RBS0-CSN::LEA-B0015
I: PcpG2-veg-EEBD-RBS0-CSN::ACCD-B0015(1836)

	1 reaction	26 reactions
Forward primer (BioBrick fw)	1µL	26 µL
Reverse primer (STK rv)	1µL	26 µL
Q-amp master mix	10µL	260 µL
ddH ₂ O	8µL	208 µL

1. distribute 20 µL per micro tube
2. pick the circled colonies with the tip and streak the gridded plates
3. Dispense the tips into one of the microtubes
4. Repeat steps stand for until all micro tubes are filled with one
5. Set up a five minute timer and begin after last tip
6. Remove tips and close lids after five minutes and spin down
7. Place into PCR machine for thermal cycle

95°C	2 min	1 cycle
95°C	30 sec	30 cycles
65°C	40 sec	
72°C	2 min	
72°C	1 min	1 cycle

Done by Sophie C, Abby K, Ryan Y
Supervised by Chars

Date: 20260725
Title: Miniprep
Purpose: Collect concentrated bacterial liquid
Procedure:

1. Transfer 1.8mL of liquid culture to a 2 mL tube
2. Centrifuge at 16,000 g for 60 seconds
3. Repeat steps 1-2 until collected all medium
4. Add 200 µl of Mx1 buffer
5. Resuspend by vortexing or pipetting
6. Add 250 µl of Mx2 buffer
7. Invert 8-10 times. Incubate at RT for 2 minutes
8. Add 350 µl of Mx3 buffer
9. Invert 3-5 times
10. Centrifuge for 5 minutes at 23,444 g
11. Add 800 µl of supernatant into assembled column
12. Centrifuge at 7000 g for 30 seconds and discard the flow through
13. Add 500 µl of WN buffer and centrifuge at 7000 g for 30 seconds and discard the flow through
14. Add 700 µl of WS buffer and centrifuge at 7000 g for 30 seconds and discard the flow through
15. Centrifuge the empty column at 10,000 g for 8 minutes
16. Assemble 1.5mL microtube w/ column
17. Add 50 µl of EB to the center of the membrane.
18. Incubate at RT for 3 minutes
19. Centrifuge at 10,000 g for 3 minutes

Biodrop results:

	concentration	A260/280
Pveg	129 ng/ µl	2.016
EEBD	123 ng/ µl	1.952
RBS0	216 ng/ µl	2.057
B0015	139 ng/ µl	1.986

Done by Olivia L and Chloe

Date: 20260725

Title: Making agarose gel(1%) (5 gels)

Purpose: Make gel for electrophoresis

Procedure:

1. Add 2g of agarose powder and 200ml of 1xTAE buffer into a flask
2. Prepare the mold, tray and comb
3. Cover with plastic wrap and poke 4-5 holes
4. Microwave for 30 seconds then agitate
5. Microwave for 5 seconds then agitate

6. Repeat step 5 until no sediments left
7. When the mixture stops smoking, add 10 µl of safety dye
8. Pipette then agitate
9. Pour into the molds
10. Let the gel solidify at room temperature

Done by Abby K and Sophie C

Date: 20260725

Title: Transformation

Purpose: Transform bacteria to agar plates

Procedure:

1. Set the water bath to 42°C and prepare ice bucket
2. Get LBA plates with glycerol
3. Thaw competent cells until they're 50% thawed and add 10µl plasmid
4. Incubate on ice for 20 minutes
5. Heat shock at 42°C for 45 seconds
6. Incubate on ice for 3 minutes
7. Add 1 ml SOC medium
8. Incubate for 1 hour at 37°C
9. Centrifuge at 5000 g for 5 minutes
10. Remove 1 ml of supernatant then vortex
11. Spread on LBA plates
12. Incubate at 37°C for 16 hours

Done by Chloe, Sophie C, Abby T

Date: 20260725

Title: Gel electrophoresis of 0725 colony PCR

Purpose: Verify plasmid sizes of 0725 colony PCR

Procedure:

1. Load 5 µl into each of the wells into each of the gels in the following layout


Gel 1:

<u>M</u>	<u>Y1</u>	<u>Y2</u>	<u>Y3</u>	<u>Y4</u>	<u>Y5</u>	<u>Y6</u>	<u>Y7</u>	<u>Y8</u>	<u>Y9</u>	<u>Y10</u>	<u>Y11</u>	<u>Y12</u>	<u>Z1</u>	<u>Z2</u>	<u>Z3</u>	<u>Z4</u>	<u>Z5</u>
<u>Z6</u>	<u>Z7</u>	<u>Z8</u>	<u>Z9</u>	<u>Z10</u>	<u>Z11</u>	<u>Z12</u>	(volume: µl)										

M: 1kb DNA ladder

Y1-Y12: Pveg-RBS0-mCherry-B0015 (bp:1164)

Z1-Z12: Pveg-RBS0-mScarlet-B0015 (bp: 1149)

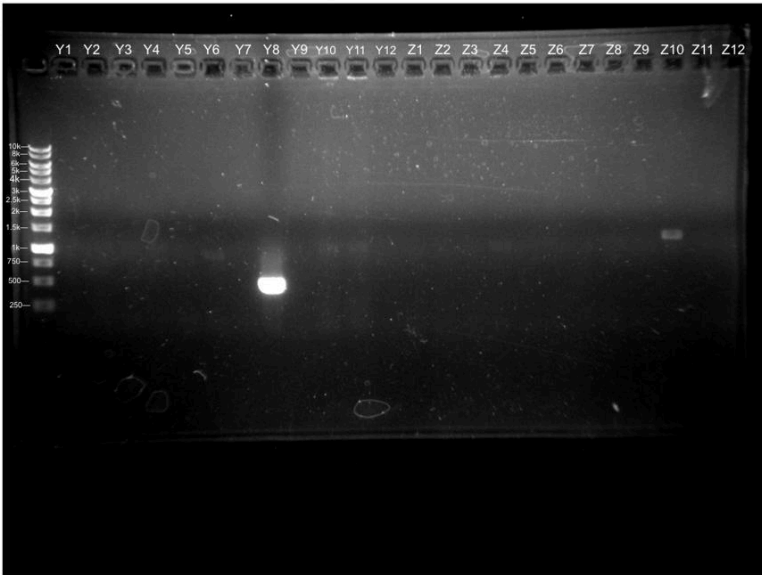


Gel 2:

M G1 G2 G3 G4 G5 G6 G7 G8 G9 G10 G11 G12 H1 H2 H3 H4
H5 H6 H7 H8 H9 H10 (volume: μ l)

M: 1kb DNA ladder

G1-G12: PcpcG2-veg-EEBD-RBS0-CSN::LEA-B0015(base pair: 1375)

H1-H10: Pveg-EEBD-RBS0-CSN::LEA-B0015(base pair:1372)

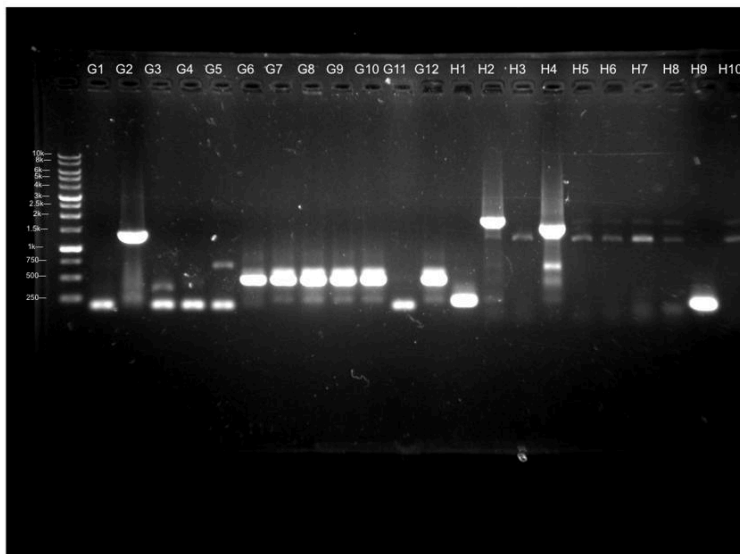


Gel 3:

M I1 I2 I3 I4 I5 I6 I7 I8 I9 I10 I11 I12 J1 J2 J3 J4 J5 J6 J7
J8 J9 J10 (volume: μ l)

M: 1kb DNA ladder

I1-I12: PcpcG2-veg-EEBD-RBS0-CSN::ACCD-B0015(base pair: 1836)

J1-J10: Pveg-EEBD-RBS0-CSN::ACCD-B0015(base pair: 2109)

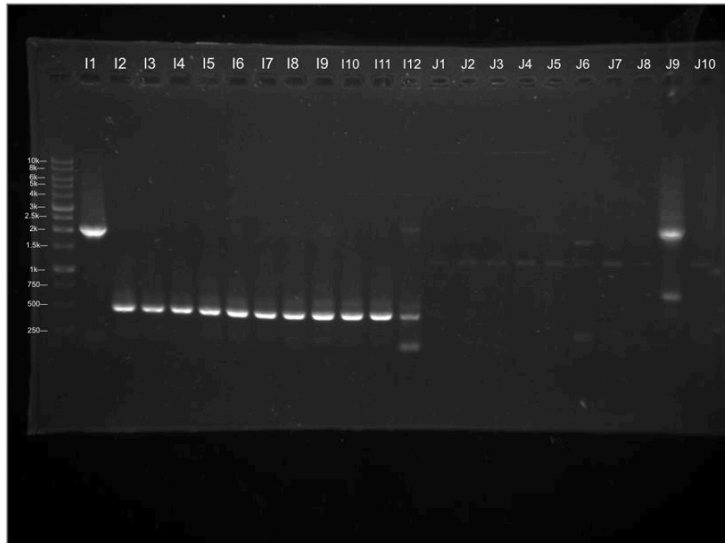

2. Resolve DNA at 100V for 30 minutes

Done by Sophie C and Abby K

Date: 20260725

Title: Colony PCR

Purpose: Verify level 0 parts(mutants)

Procedure:

1. Prepare the following reaction and master mix

B: pSTK-SamyQ-ACCD-v2t1

F: pSTK-SamyQ-ACCD-v3t1

	1 reaction	20 reactions
Forward primer (BioBrick-fv)	1 μ L	20 μ L
Reverse primer (ACCD R)	1 μ L	20 μ L
Q-amp master mix	10 μ L	200 μ L
ddH ₂ O	8 μ L	16/0 μ L

A: pSTK-SamyQ-LEA-v2t1

	1 reaction	10 reactions
Forward primer (BioBrick-fv)	1 μ L	10 μ L

Reverse primer X(LEA R)	1µL	10 µL
Q-amp master mix	10µL	100 µL
ddH ₂ O	8µL	80 µL

H: pSTK-SamyQ-AmilCP-v3t1

	1 reaction	10 reactions
Forward primer (BioBrick-fv)	1µL	10 µL
Reverse primer (AmilCP R)	1µL	10 µL
Q-amp master mix	10µL	100 µL
ddH ₂ O	8µL	80 µL

2. distribute 20 µL per micro tube
3. pick the circled colonies with the tip and streak the gridded plates
4. Dispense the tips into one of the microtubes
5. Repeat steps stand for until all micro tubes are filled with one
6. Set up a five minute timer and begin after last tip
7. Remove tips and close lids after five minutes and spin down
8. Place into PCR machine for thermal cycle

95°C	2 min	1 cycle
95°C	30 sec	30 cycles
65°C	40 sec	
72°C	1 min 30/ sec	
72°C	1 min	1 cycle

Done by Ryan Y, Olivia L, and Chloe
Supervised by Chars and Gabriel

Date: 20260725

Title: ½ MS medium preparation (2000mL)

Purpose: Prepare ½ MS medium for further Arabidopsis planting(hydroponics)

Procedure:

1. Measure and add to the flask
 - a. ½ MS medium(without sucrose): 4.6g(2.3*2)
 - b. MES buffer: 1g(0.5*2)
 - c. ddH₂O: 2L(1*2)
2. Use the hot plate to stir(with a stirring bar) until no sediments are left in the solution

3. Determine pH of the solution and adjust the pH of the ½ MS medium using a pH probe to around 5.7-5.8pH
 - a. Add 1M KOH to increase the pH
 - i. Original pH: 4.1
 - ii. Final pH: 5.7
4. Split the 2000mL medium into 2 500mL flask and the original beaker
5. Wrap the neck with aluminum foil and place the indicator tape on top
6. Autoclave at 127°C for 20 minutes
7. Wait for the pressure to reach 0 before opening the autoclave
8. Pour 300mL of medium into each of the hydroponic boxes

Done by Sara and Sophie H

Supervised by Alex

Date: 20260725

Title: Gel electrophoresis

Purpose: Check base pair of colony PCR product

Procedure:

1. Load samples into the gels in the following layout

Gel 1:


<u>M</u>	<u>A1</u>	<u>A2</u>	<u>A3</u>	<u>A4</u>	<u>A5</u>	<u>A6</u>	<u>A7</u>	<u>A8</u>	<u>B1</u>	<u>B2</u>	<u>B3</u>	<u>B4</u>	<u>B5</u>	<u>B6</u>	<u>B7</u>	<u>B8</u>
5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5	5

(volume: µl)

M: 1kb DNA ladder

A: pSTK-SamyQ-Lea-v2t1

B: pSTK-SamyQ-ACCD-v2t1

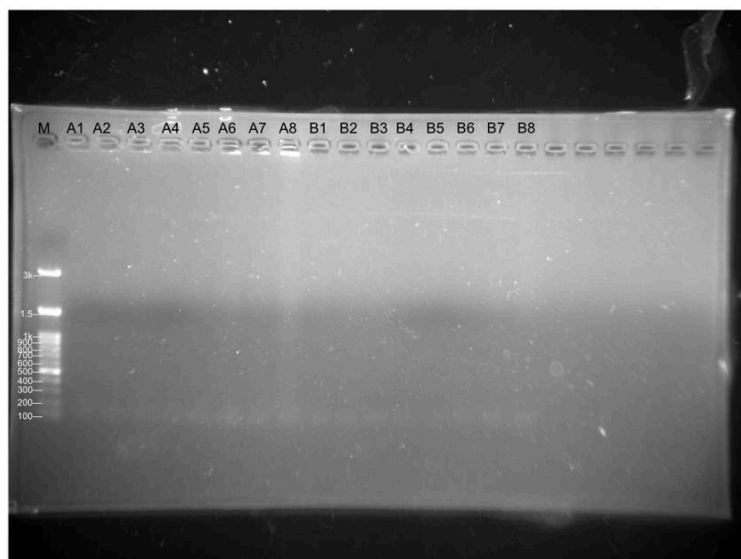

Gel 2:

<u>M</u>	<u>H1</u>	<u>H2</u>	<u>H3</u>	<u>H4</u>	<u>H5</u>	<u>H6</u>	<u>H7</u>	<u>H8</u>	<u>F1</u>	<u>F2</u>	<u>F3</u>	<u>F4</u>	<u>F5</u>	<u>F6</u>	<u>F7</u>	<u>F8</u>
----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------

5 5 5 5 5 5 5 5 5 5 5 5 5 5 5 5
(volume: μ l)

M: 1kb DNA ladder

H: pSTK-SamyQ-AmilCP-v3t1

F: pSTK-SamyQ-ACCD-v3t1



2. Resolve at 100V for 30 minutes

Done by Sophie C and Abby T

Date: 20260725

Title: Transformation

Purpose: Transform bacterial liquid on agar plates

Procedure:

A: Pveg-RBS0-Ho1:pcyA-B0015

L: Pveg-EEBD-RBS0-GanA-Gala-B0015

1. Set the water bath to 42°C and prepare ice
2. Get LBA plates with glycerol
3. Thaw competent cells until they're 50% thawed and add 10 μ l of plasmid
4. Incubate on ice for 20 min
5. Heatshock at 42°C for 45 seconds
6. Incubate on ice for 3 minutes
7. Add 1 mL SOC medium
8. Incubate for 1 hour at 37°C
9. Centrifuge at 5000 g for 5 minutes
10. Remove 1 ml of supernatant then vortex
11. Spread on LBA plates
12. Incubate at 37°C for 16 hours

Done by Chloe and Sophie C

Date: 20260725

Title: Liquid Culture

Purpose: Transfer colony to medium

Procedure:

A1-3, A1-4, A1-8, A1-9, A2-6 : Pveg-RBS0-Ho1:PcyA-B0015

H3, H7 : Pveg-EEBD-RBS0-csn:LEA-B0015

G2:PcpcG2-veg-EEBD-RBS0-csn:LEA

Z10:Pveg-RBS0-mScarlet-B0015

J9:Pveg-EEBD-RBS0-csn:ACCD- B0015

I1:PcpcG2-veg-EEBD-RBS0-csn:ACCD-B0015

1. Label 11 liquid culture tubes
2. Turn on Bunsen Burner, add 3ml TB medium and 3 μ l of ampicilin to all 11 tubes
3. Dip pipette into colonies and dispense tips into respective tubes
4. Incubate at 37°C overnight

Done by Abby T and Sophie C

Supervised by Chars


Date: 20260725

Title: CH biotech protectant testing on box 1, 2, 3, 4

Purpose: test on when we should add protectants

Procedure:

1. Autoclave one tip box for 20 minutes at 127°C and cleanse the 3D printed board with 10% bleach
2. Add 300 mL of $\frac{1}{2}$ MS medium into the tip box, and label this box, and label this box as box 4
3. Add 600 μ L of CH biotech protectant into the medium of box 4
4. Place the 3d printed board onto the medium in box 4
5. Transfer 3 plants each from box 1 and box 3 to box 4
6. Transfer 4 plants from box 2 to box 4

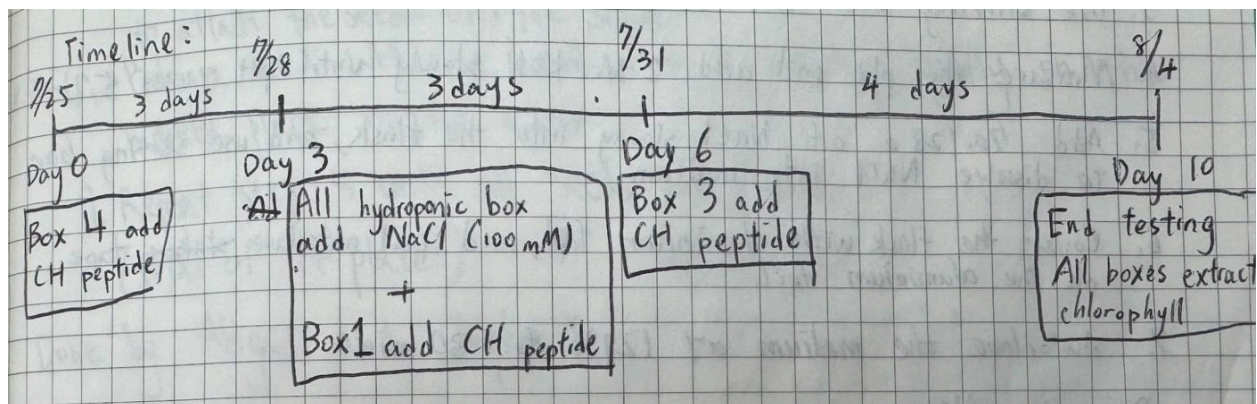

Done by Alex
Supervised by Gina

Date: 20260725

Title: Prepare 3M NaCl with ½ MS medium

Purpose: hydroponics testing with CH biotech protectant

Procedure:

1. Prepare 400ml ½ MS medium with the following protocol

½ MS powder (w/o sucrose)	0.92g
MES buffer	0.2g
Distilled water	400ml

2. Add the components into 50ml flask
3. Use stirring bar to dissolve the powder into water
4. Measure the pH and add 1M KOH slowly until pH reaches 5.7-5.8
5. Add 70.128g of NaCl slowly into the flask, and use stirring bar to dissolve NaCl into medium
6. Cover the flask with aluminum foil, and put autoclave tape on the aluminum foil
7. Autoclave for 20 minutes at 127°C

Done by Alex
Supervised by Gina

2026/7/26 

Date: 20260726

Title: Miniprep

Purpose: extract plasmid DNA from liquid culture of colony PCR products(0725)

Procedure:

A1-3, A1-4, A1-8, A1-9, A2-6 : Pveg-RBS0-Ho1:PcyA-B0015

H3, H7 : Pveg-EEBD-RBS0-csn:LEA-B0015

G2:PcpcG2-veg-EEBD-RBS0-csn:LEA

Z10:Pveg-RBS0-mScarlet-B0015

J9:Pveg-EEBD-RBS0-csn:ACCD- B0015

I1:PcpcG2-veg-EEBD-RBS0-csn:ACCD-B0015

1. Pellet bacterial overnight culture by centrifugation at 23444 g for 30 seconds at RT
2. Resuspend pelleted bacterial cells in 250 μ l Buffer P1 and transfer to a microcentrifuge tube
3. Add 250 μ l Buffer P2 and mix thoroughly by inverting the tube 4-6 times until the solution becomes clear(turn blue)
4. Add 350 μ l Buffer N3 and mix immediately and thoroughly by inverting the tube 4-6 times(turn colorless)
5. Centrifuge for 10 min at 17900 g in a table top microcentrifuge
6. Apply supernatant from step 5 to QIA prep 2.0 spin column by pipetting. Centrifuge for 30 seconds and discard the flow through
7. Wash by adding 0.5ml Buffer PB. Centrifuge for 30 seconds(17900 g) and discard the flow through
8. Wash by adding 0.75ml Buffer PE. Centrifuge for 30 seconds(17900 g) and discard the flow through
9. Centrifuge for 1 min to remove residual wash buffer
10. Place QIA prep 2.0 column in a clean 1.5 ml centrifuge tube. Add 50 μ l sterilized water to the center of membrane
11. Stand for 1 min
12. Centrifuge at 17900 g for 1 min

Biodrop results

	concentration(ng/ μ l)	A260/A280
A1-3	312.0	1.914
A1-4	420.0	1.918
A1-8	464.0	1.909
A1-9	467.0	1.914
A2-6	377.0	1.914

G2	133.0	1.928
H3	143.0	1.882
H7	108.0	1.929
I1	162.0	1.884
J9	410.0	1.898
Z10	274.0	1.903

Done by Abigail and Ryan W
Supervised by Chars

Date: 20260726

Title: Colony PCR of A, L

Purpose: screen for successful bacterial inserts from 7/25 transformation

Procedure:

A: Pveg-RBS0-Ho1:pcyA-B0015

L: Pveg-EEBD-RBS0-GanA-Gala-B0015

1. Prepare the reaction as follows

A. pSTK-1A-Pveg-RBS0-Ho1:pcyA-B0015 (pcr product bp: 2111)

	1 reaction	12 reactions
Forward primer (Biobrick Prefix)	1µL	12 µL
Reverse primer (pSTK_Seq.Rv)	1µL	12 µL
Q-amp master mix	10µL	120 µL
ddH ₂ O	8µL	96 µL

L: pSTK-EXP1-Pveg-EEBD-RBS0-Spacer-GanA-Gala-B0015 (pcr product: 2993)

	1 reaction	6 reactions
Forward primer (Exp-seq-fw)	1µL	6 µL
Reverse primer (Exp-seq-rv)	1µL	6 µL
Q-amp master mix	10µL	60 µL
ddH ₂ O	8µL	48 µL

1. Circle 11 colonies from A, 5 colonies from L
2. Turn on Bunsen Burner and Pick up 11 colonies from A, 5 colonies from L with tip and streak the labelled grid
3. Dispense the tips with colonies into the PCR strips and wait 5 minutes
4. Take out the tips and close the lid
5. Place into PCR machine for thermal cycle

95°C	2 min	1 cycle
95°C	30 sec	30 cycles
55°C	20 sec	
72°C	3 min (largest bp:2993)	
72°C	5 min	1 cycle

Done by Sophie L and Abby T
Supervised by Chars

Date: 20260726

Title: Prepare 0.8% Agarose Gel

Purpose: make gels for gel electrophoresis

Procedure:

1. Assemble mold, tray and comb
2. Measure 0.32g of agarose powder and 40 ml of 1xTAE buffer
3. Pour agarose powder and 1XTAE buffer into conical flask
4. Microwave for 30 seconds then agitate it. Repeat this step until no sediments remain
5. Wait until no longer smoking then add 2 µl of safety dye then agitate
6. Pour melted gel into mold
7. Store in RT until solidified

Done by Sophie L

Date: 20260726

Title: Enzyme Digestion

Purpose: Check plasmid size before DNA sequencing

Procedure:

A1-3, A1-4, A1-8, A1-9, A2-6 : Pveg-RBS0-Ho1:PcyA-B0015

H3, H7 : Pveg-EEBD-RBS0-csn:LEA-B0015

G2:PcpcG2-veg-EEBD-RBS0-csn:LEA

Z10:Pveg-RBS0-mScarlet-B0015

J9:Pveg-EEBD-RBS0-csn:ACCD- B0015

I1:PcpcG2-veg-EEBD-RBS0-csn:ACCD-B0015

1. Prepare the following reaction:

A1-3(full length base pair: 4032):

10x CutSmart	2 μ L
DNA	1.603 μ L
ddH ₂ O	15.397 μ L
EcoRI(1 site)	0.5 μ L

A1-4(full length base pair: 4032):

10x CutSmart	2 μ L
DNA	1.190 μ L
ddH ₂ O	15.81 μ L
EcoRI(1 site)	0.5 μ L

A1-8(full length base pair: 4032):

10x CutSmart	2 μ L
DNA	1.078 μ L
ddH ₂ O	15.922 μ L
EcoRI(1 site)	0.5 μ L

A1-9(full length base pair: 4032):

10x CutSmart	2 μ L
DNA	1.071 μ L
ddH ₂ O	15.929
EcoRI(1 site)	0.5 μ L

A2-6(full length base pair: 4032):

10x CutSmart	2 μ L
DNA	1.326 μ L
ddH ₂ O	15.674 μ L
EcoRI(1 site)	0.5 μ L

H3(full length base pair: 7030):

NE Buffer	2μL
DNA	3.497μL
ddH ₂ O	13.503μL
Bsmbl(1 site)	0.5 μL

H7(full length base pair: 7030):

NE Buffer	2μL
DNA	4.630μL
ddH ₂ O	12.37μL
Bsmbl(1 site)	0.5 μL

G2(full length base pair: 3296):

NE Buffer	2μL
DNA	3.759μL
ddH ₂ O	13.241μL
Bsmbl(2 sites)	0.5 μL

Z10(full length base pair: 2780):

10x CutSmart	2μL
DNA	1.825μL
ddH ₂ O	15.175μL
EcoRI(1 site)	0.5 μL

J9(full length base pair: 7767):

10x CutSmart	2μL
DNA	1.220μL
ddH ₂ O	15.78μL
BbsI(1 site)	0.5 μL

l1(full length base pair 3757):

10x CutSmart	2 μ L
DNA	3.086 μ L
ddH ₂ O	13.914 μ L
EcoRI(1 site)	0.5 μ L

2. Prepared the reaction according to the formula: add water, DNA, buffer first before enzymes
3. Incubated the formula in PCR machine

55°C	15 hours
80°C	20 min

Done by Abby T, Ryan W, and Abigail
Supervised by Chars

Date: 20260726

Title: PCR

Purpose: Amplify plasmid number

Procedure:

1. Prepare the following reaction for:

L3SIPS2-L3S2P56

MF001

MF002

Reagent	1 rxn (50 μ L)	3 rxn (150 μ L) *3 (master mix)
10x ExTaq Buffer	5 μ L	15 μ L
2.5 μ M dNTP mix	4 μ L	12 μ L
10 μ M forward primer	1 μ L	3 μ L
10 μ M reverse primer	1 μ L	3 μ L

DNA temp	1 μ L	3 μ L
PCR grade water	37.75 μ L	113.25 μ L
Ex Taw polymerase	0.25 μ L	0.75 μ L

2. Distribute 50 μ L of the 3 rxn mixture into each tube for the corresponding part
3. Place into PCR machine, applying the thermalcycle below:

98 °C	30 seconds	1x
98 °C	10 seconds	25x
65 °C	30 seconds	
72 °C	30 seconds	
72 °C	2 minutes	1x
25 °C	0 minutes	

Done by Chloe

Date: 20260726

Title: SOC Medium Preparation

Purpose: prepare soc medium for transformation

Procedure:

1. Prepare 3.17 g of SOC broth powder.
2. Prepare 100 mL ddH₂O
3. Mix together by using a stirring rod.

Done by Abby T

Supervised by Chars

Date: 20260726

Title: Transformation

Purpose: redo transformation of 0725

Procedure:

A: Pveg-RBS0-Ho1:pcyA-B0015

L: Pveg-EEBD-RBS0-GanA-Gala-B0015

1. Set the water bath to 42°C and prepare ice
2. Get LBA plates with glycerol
3. Thaw competent cells until they're 50% thawed and add 10 μ L of plasmid
4. Incubate on ice for 20 min

5. Heatshock at 42°C for 45 seconds
6. Incubate on ice for 3 minutes
7. Add 1 mL SOC medium
8. Incubate for 1 hour at 37°C
9. Centrifuge at 5000 g for 5 minutes
10. Remove 1 ml of supernatant then vortex
11. Spread on LBA plates
12. Incubate at 37°C for 16 hours

Done by Chloe and Abigail

Supervised by Chars

2026/7/27 

Date: 20260727

Title: Miniprep

Purpose: Purify desired plasmid from liquid culture

Procedure:

D10: pSTK-0AB-PtagC+RBS

H2: pSTK-0D-bkdB Cold double terminator

G12: pSTK-0D-L3S2P21

E8: pSTK0-B-RBS2

D11: pSTK-0B-3P promoter

D7: pSTK-0A-Phyperspank

1. Transfer 1 mL of liquid culture to a 2mL tube
2. Centrifuge at 16000 g for 60 seconds
3. Discard the supernatant
4. Repeat step 1-3 until all medium are collected
5. Add 200 µl Mx1 buffer
6. Resuspend by vortexing or pipetting
7. Add 250 µl Mx2 buffer
8. Invert slowly 8-10 times, incubate at room temperature for 2 min
9. Add 350 µl Mx3 buffer
10. Invert slowly 4-5 times
11. Centrifuge at 16000 g for 5 min
12. Transfer around 800 µl supernatant into the column
13. Centrifuge at 7000 g for 30 min
14. Discard the flow through
15. Add 500 µl of WN buffer, centrifuge 7000 g for 30 seconds, discard the flow through
16. Add 500 µl of WS buffer, centrifuge 7000 g for 30 seconds, discard the flow through
17. Centrifuge 1000 g for 8min
18. Assemble 1.5mL microtube
19. Add 50 µl of EB buffer on the center of the membrane
20. Incubate in RT for 3 min
21. Centrifuge 10000 g for 3min

Biodrop results

	concentration(ng/ µl)	A260/A280	A260/A230
D10	119	1.725	0.967
H2	39	1.857	1.114
G12	119	1.653	0.902
E8	101	1.683	1.01

D11	123	1.662	0.932
D7	110	1.774	0.887

Done by Sarah, Naomi, Sophie L

Date: 20260727

Title: Make 0.8% Agarose gel(small+big)

Purpose: electrophoresis for 0726 enzyme digestion

Procedure:

1. Assemble mold, tray and comb
2. Measure 0.48g of agarose powder and 60mL of 1x TAE buffer
3. Pour agarose powder and 1xTAE buffer into conical flask
4. Microwave for 30 seconds then agitate it
5. Repeat this step until no sediments
6. Wait until no smoke and add 3 µl of safety dye
7. Agitate it
8. Pour the melted gel into the mold
9. Wait until the gel solidifies in RT

Done by Sophie C

Date: 20260727

Title: Colony PCR of I, L

Purpose: Amplify plasmid number

Procedure:

I: pSTK-PcpcG2-veg-EEBD-RBS0-csn:ACCD-B0015

L: pSTK-Exp1-Pveg-GanA:Galactosidase DH5a Assembly

1. Prepare the reaction as follows

I: pSTK-PcpcG2-veg-EEBD-RBS0-csn:ACCD-B0015

	1 reaction	12 reactions
Forward primer (PcpcG2-veg fw)	1µL	12 µL
Reverse primer (ACCD rv)	1µL	12 µL
Q-amp master mix	10µL	120 µL
ddH ₂ O	8µL	96 µL

L: pSTK-Exp1-Pveg-GanA:Galactosidase DH5a Assembly(x26)

	1 reaction	26 reactions
Forward primer (Pveg Fw)	1µL	26 µL

Reverse primer (Exp seq4rv)	1µL	26 µL
Q-amp master mix	10µL	260 µL
ddH ₂ O	8µL	208 µL

1. Circle 10 colonies from I, 24 colonies from L
2. Turn on Bunsen Burner and pick up 11 colonies from A, 5 colonies from L with tip and streak the labelled grid
3. Dispense the tips with colonies into the PCR strips and wait 5 minutes
4. Take out the tips and close the lid
5. Place into PCR machine for thermal cycle

95°C	2 min	1 cycle
95°C	30 sec	30 cycles
55°C	20 sec	
72°C	3 min (largest bp:2993)	
72°C	5 min	1 cycle

Done by Sophie C

Supervised by Gabriel

Date: 20260727

Title: Electrophoresis and Gel extraction

Purpose: Isolate desired DNA fragment from enzyme digestion

Procedure:

1. Mix 50 µL DNA with 10 µL 6X loading dye.
2. Loading sample by the following layout.

Gel 1

M	A	
5	60	volume (µL)

M= 1kb DNA ladder (blue)

A= MF001

Gel 2

M	B	
5	60	volume (μL)

M= 1kb DNA ladder (blue)

B= MF002

Gel 3

M	C	
5	60	volume (μL)

M= 1kb DNA ladder (blue)

C= L3SIPS

3. Resolve the DNA by 100V 30 minutes.

Done by Abby T, Chloe

Supervised by Gabriel

Date: 20260727

Title: Gel Extraction

Purpose: DNA fragment isolation and purification.

Procedure:

1. Use a sanitized sharp knife to cut the gel band containing the desired DNA fragment under the UV light box.
2. Place each 3 gel bands into separated and labeled 2ml tubes (Mf002, Mf001, and L3SIPS).
3. Add 700 μL of GP Buffer into each of the tubes.
4. Place the tubes in a dry bath on 60°C to dissolve the gel slices, check every 2 minutes until entirely dissolved.
5. Load no more than 700 microliter of melted gel into the collection tube
6. Centrifuge for 30seconds at 7000 g. Discard the flow-through. Repeat step 5-6 step for the rest of the mixture.
7. Wash the column once with 0.5ml of WN Buffer by centrifuging for 30 seconds. Discard the flow through.
8. Wash the column once with 0.5ml of WS Buffer by centrifuging for 30 seconds. Discard the flow through.
9. Centrifuge the column at full speed 23444g for 3 minutes to remove residual ethanol.

10. Place the column onto a new centrifuge tube. Add 15 μ L of EB buffer onto the center of the membrane.
11. Stand the column for 2 minutes and centrifuge at full speed 23444g for 2 minutes to elute DNA.
12. Turn on BioDrop, adjust the setting to pathlength unit (concentration).
13. Clean the sample port with 1 μ L of Nucleus free water.
14. Calibrate.
15. Load 1 μ L of DNA from the 2 tubes each time and record down the results.
16. Clean the port with Nucleus Free water and close BioDrop.

Done by Abby T, Chloe

Supervised by Gabriel

Date: 20260727

Title: Plant seeds in circular $\frac{1}{2}$ MS medium (neutral) plates

Purpose: plant seeds to vernalize for 4 days before testing on 正瀚 protectants

Procedure:

1. Distribute (columbia-0 kyle seeds into individual tubes
2. Open bsc and put the circular 0.8% $\frac{1}{2}$ ms medium plates from 0726, see template (horizontal 5 seeds x vertical 5 seeds = 25 seeds / plate), seeds in tubesm sterilized water and bamboo sticks in the bsc
3. Place the plate on the seed template
4. Stick the bamboo stick in sterilized water so it is easier to attach the seeds on the stick
5. Pick a seed with the bamboo stick and tap lightly on the marked place on the agar
6. Repeat step 5 until all seeds are planted (total of 14 plates)

Done by Abigail and Sophia Y.

Date: 20260727

Title: Make experiment set 6 salinity plates

Purpose: to prepare salinity plates for experiment set 6

Procedure:

1. Melt the solidified $\frac{1}{2}$ MS agar (in tubes 30mL) by putting the 2 tubes into microwave

For 300 mL of agar:

6 agar plates — only made 2 agar plates (control X control and 150 mM NaCl x control with the stored $\frac{1}{2}$ MS solidified agar, 3M NaCl liquid and 100 mM trehalose liquid)

- Control (no trehalose) $\times$ 75 mM NaCl
- Control (no trehalose) $\times$ 100 mM NaCl
- 5 mM tre $\times$ control (no NaCl)
- 5 mM tre $\times$ 75 mM NaCl
- 5 mM tre $\times$ 100 mM NaCl
- 5 mM tre $\times$ 150 mM NaCl

1. Prepare the following
 - a. $\frac{1}{2}$ MS medium powder (w/o sucrose): 0.6 g
 - b. MES buffer: 0.15 g
 - c. Distilled water: 300 mL
 - d. Sucrose (1%): 3 g
 - e. Agar powder (1%): 3 g
2. Add the contents into a 1000 mL flask
3. Use a hot plate to stir (with a magnetic stirring rod) until no sediments are left in the solution
4. Use the pH probe and measure the pH of the solution, make sure the pH is 5.7–5.8
 - a. If pH too low, add 0.1 M KOH slowly to increase it
5. Pour 300 mL of mixture out into another flask to separate into six flasks, label with:
 - a. [0 mM] NaCl, 5 mM trehalose
 - b. [75 mM] NaCl, 5 mM trehalose
 - c. [100 mM] NaCl, 5 mM trehalose
 - d. [150 mM] NaCl, 5 mM trehalose
 - e. [75 mM] NaCl, 0 mM trehalose
 - f. [100 mM] NaCl, 0 mM trehalose
6. Add 219.2 mg of NaCl into the flask labelled [75 mM NaCl]
7. Add 292.2 mg of NaCl into the flask labelled [100 mM NaCl]
8. Add 438.3 mg of NaCl into the flask labelled [150 mM NaCl]
9. To their respective trehalose concentrations, add:

Trehalose concentration	Amount of trehalose
0mM	0mg
5mM	85.575mg

10. Wrap the necks of the flasks with aluminum foil and stick a piece of indicator tape on top
11. Autoclave for 20 min at 127°C
12. Wait for pressure to reach 0 before retrieving flasks

Done by: Abby, Abigail

Supervised by: Gabriel

ELN by Naomi

2026/7/28 

Date: 20260728

Title: Prepare ½ MS medium (Hydroponic medium)

Purpose: prepare medium for hydroponic testing (正瀚protectant + NaCl)

Procedure:

1. The following protocol is for 2000mL of ½ MS medium
 - a. ½ MS medium powder without sucrose: 4.6g
 - b. MES buffer: 1g
 - c. Distilled water: 2000mL
2. Add the contents in the beaker
3. Use a hot plate to stir (with a magnetic stir bar) until no sediments are left in the solution
4. Use the pH probe and measure the pH of the solution, ensure the pH is 5.7~5.8
 - a. Actual pH: 4.1, pH after adjustment: 5.7
5. If the pH is too low, add 0.1M KOH slowly to increase it
6. Wrap the beaker with aluminum foil and stick a piece of autoclave indicator tape on the top
7. Autoclave for 20 minutes at 127°C.

Done by Abigail, AudreyH, Sophie H

Date: 20260728

Title: PCR

Purpose: To amplify the number of DNA circuit

Procedure:

1. Prepare the following reaction for L3S1PS2-L3S2P56, MF001, MF002

Reagent	1 rxn	3 rxn(150µL)
10x EX Taq buffer	5µL	15µL
2.5um dNTP mix	4µL	12µL
10 µL fw	1µL	3µL
10µL rv	1µL	3µL
DNA temp	1µL	3µL
PCR-grace H ₂ O	37.75µL	113.25µL
Ex taq poly	0.25µL	0.75µL

2. Distribute 50µL of the 3rxn mixtures into each tube for the corresponding parts
3. Place in PCR machine and assemble the thermalcycle as below:

98°C	30sec	x1
98°C	10sec	
65°C	30sec	}x25
72°C	30sec	
72°C	2min	
25°C	0sec	x1

Done by Sophie C, Chloe

Date: 20260728

Title: Making 2% agarose gel x3

Purpose: Make gel for electrophoresis

Procedure:

1. Assemble mold, tray, and comb
2. Weight 1.2 g of agarose powder and measure 60mL of 1xTAE buffer
3. Add the agarose powder and the 1x TAE buffer into conical flask
4. Microwave for 30 sec then agitate it
5. Repeat this step until no sediments
6. Wait until no smoke, add 3µL of safety dye and agitate it as well
7. Pour the melted gel into the tray and use tip to squeeze the air and pop the bubble
8. Wait until the gel solidifies

Done by Sophie C, Chloe

Date: 20260728

Title: Gel extraction

Purpose: recollect PCR product from gel after electrophoresis

Procedure:

1. Use a sanitized sharp knife to cut the gel band containing the desired DNA fragment under the UV light box
2. Place each 3 gel bands into separated 2ml tubes (MF001, MF002, L3S1PS)
3. Add 300µL GP buffer into each tube
4. Place tubes in a 60°C dry bath to dissolve the slices(check every 2min until fully dissolved)
5. Load no more than 700µL of melted gel in the collection tube
6. Centrifuge 7000 g for 30secs, then discard the flow through and repeat steps 5-6 for the rest of the mixture
7. Wash column with 500µL WN buffer and centrifuge 7000 g for 30secs, then discard the flow through
8. Wash column with 500µL WS buffer and centrifuge 7000 g for 30secs, then discard the flow through
9. Centrifuge empty column at 23,444 g for 3min

10. Place column on a new 1.5mL tube and add 15µL Eb buffer to the center of the membrane
11. Stand the column for 2min and centrifuge at 23,444 for 2min to elute DNA
12. Measure concentration and A260/280 of plasmid on biodrop.

Done by Chloe, Sophie C.

Date: 20260728

Title: GGA of A, L

Purpose: To assemble DNA construct

Procedure:

A: (STK) Module 1(1A): Pveg-RBS-ho1:pcyA-L3S2P21

1. prepare the reaction as follow:

Part	Original Amount	Dilution
Vector (1A) 115 ng/µL, 3120 bp	0.18 µL	1.76 µL (dilute 1/10)
PtagC+RBS 725 ng/µL, 2184 bp	0.23 µL	2.27 µL (dilute 1/10)
Ho1:pcyA 238 ng/µL, 1556 bp	0.08µL	1.7µL (dilute 1/20)
L3S2P21 170 ng/µL, 2145bp	0.16 µL	1.64 µL (dilute 1/10)
T4 DNA ligase buffer	2 µL	2 µL
Sterile water	8.63 µL	//
T4 DNA ligase	1 µL	1 µL
BsaI-HFv2	1 µL	1 µL

Total	20 μ L	20 μ L
-------	------------	------------

L: (STK) Module 5: 3P-EEBD-RBS0-spacer-GanA-Gala-spacer-B0015

Part	Original Amount	Dilution
Vector (Exp1) 162 ng/ μ L, 6929 bp	0.28 μ L	1.39 μ L (dilute 1/10)
3P 106 ng/ μ L, 2909 bp	0.36 μ L	1.8 μ L (dilute 1/10)
EEBD-RBS0 123 ng/ μ L, 2392 bp	0.25 μ L	1 μ L (dilute 1/20)
GanA-Gala 177 ng/ μ L, 4148bp	0.30 μ L	1.5 μ L (dilute 1/10)
Spacer CX	0.107 μ L	1.07 μ L
Spacer CZ	0.116 μ L	1.16 μ L
L3S2P21	0.16 μ L	1.64 μ L
T4 DNA ligase buffer	2 μ L	2 μ L
Sterile water	6.44 μ L	6.44 μ L
T4 DNA ligase	1 μ L	1 μ L
BsaI-HFv2	1 μ L	1 μ L
Total	20 μ L	20 μ L

Use the thermocycler to run the program below:

Temperature	Time	Cycles
37 °C (Restriction enzyme digestion)	5 minutes	30 cycles
16 °C (Ligation)	5 minutes	
37 °C (Restriction enzyme digestion)	10 minutes	1 cycle
80 °C (Heat inactivation of BsaI restriction enzyme)	20 minutes	1 cycle

Done by Sophie C, Chloe

Date: 20260728

Title: Liquid culture

Purpose: To transfer the colony on agar plate in medium

Procedure:

D7: pSTK-0A-Phyperspank

D10: pSTK-0AB-PtagC+RBS

D11: pSTK-0A-P3P

E8: pSTK-0B-RBS2

H2: pSTK-0D-bkaB-Co1D-double terminator

G12: pSTK-0D-L3S2P21

1. Label 6 liquid culture tubes
2. Turn on bunsen burner and add 3mL TB medium and 3µL of chloramphenicol in each tube
3. Add glycerol stock into each tube
4. Incubate in shaker overnight

Done by Sophie C, Gabriel

Date: 20260728

Title: Miniprep

Purpose: Purify desired plasmid

Procedure:

A11: Pveg-RBS0-Ho1:pcyA-B0015

1. Transfer 1.8mL of A11 liquid culture to a 2mL microtube
2. Centrifuge at 23444 g for 30 sec, then discard the supernatant
3. Repeat steps 1-2 until all medium is collected, remember to empty the medium in the tips too

4. Add 200µL of MX1 buffer, vortex or pipette to resuspend it
5. Add 250µL of MX2 buffer, vortex or pipette to resuspend it
6. Invert slowly 8-10 times, Incubate at RT for 2mins
7. Add 350µL of MX3 buffer, and invert slowly 4-5 times
8. Centrifuge for 5min at 16000 g
9. Transfer about 900µL of supernatant into assembled column
10. Centrifuge at 7000 g for 30secs and discard the flow through
11. Wash the column with 500µL WN buffer, then centrifuge at 700 g for 30sec
12. Discard the flow through
13. Wash the column with 500µL WS buffer, then centrifuge at 700 g for 30sec
14. Discard the flow through
15. Centrifuge for 8min at 10000 g
16. Place the new column into new 1.5ml microtube
17. Add 50µL Eb to the center of th membrane
18. Incubate for 30min at RT
19. Centrifuge at 10000 g for 30sec

Biodrop Results:

A11	Concentration	A260/230	A260/280
	75ng/µL	1.230	2.032

Done by Sophie C.

Date:20260728

Title:Gel electrophoresis

Purpose: To check the results of DNA circuit after GGA and PCR

Procedure:

1. **Load the samples in the following order:**


<u>M</u>	<u>L25</u>	<u>L26</u>	<u>L27</u>	<u>L28</u>	<u>L29</u>	<u>L30</u>	<u>L31</u>	<u>L32</u>	<u>L33</u>	<u>L34</u>	<u>L35</u>	<u>L36</u>	<u>L37</u>	<u>L38</u>
5	5	5	5	5	5	5	5	5	5	5	5	5	5	5

<u>L39</u>	<u>L40</u>	<u>L41</u>	<u>L42</u>	<u>L43</u>	<u>L44</u>	<u>L45</u>	<u>L46</u>	<u>L47</u>	<u>L48</u>
5	5	5	5	5	5	5	5	5	5

2. Resolve the gel at 100V for 30sec

Result:

20260728 Colony PCR (L in E. Coli)

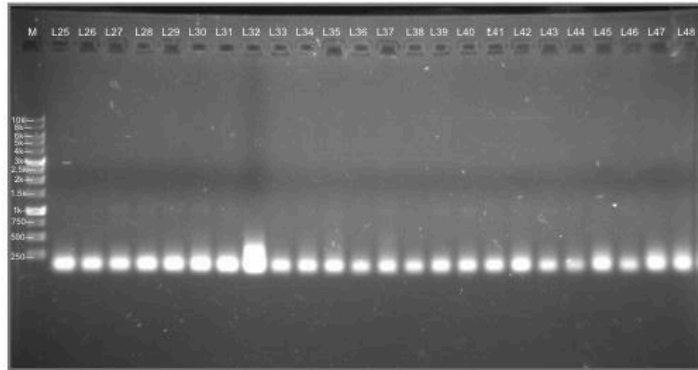

M: 1kb DNA ladder
L: Pveg-EEBD-RBS0-GanA-Gala-B0015
(bp:
• Forward primer (Pveg fw)
• Reverse primer (EXP-seq4 rv)
• bp: 2946

All fail

Done by Sophie C

Date: 20260728

Title: Transfer plants from incubator to salinity plates

Purpose: start salinity stress testing of experiment 6 (0, 75, 150, 150mM NaCl, 0, 5mM trehalose)

Procedure:

1. Take plates (incubated for 4 days) out of growth chamber
2. Open BSC
3. Place 2 vertical plate holders, tweezers, [unclear/crossed out], the salinity plates and plates that incubated for 4 days into the BSC
4. Choose plants that have similar root lengths to transfer after visually comparing the incubated plants
5. Lightly lift the chosen plant using a tweezer below the leaves
6. Place the plant on the salinity plates
7. Use a marker to draw a straight line along a ruler using the seed template
8. Align the root endpoints with the drawn line and use a marker to indicate the endpoint for further image analysis
9. Seal the plates with breathable tape and label the test group and transferred date
10. Place the vertical plates on the racks in the growth chamber at a slanted angle towards the back (front faces forward)

Done by: Sara, Sophie H, Abby T

Date: 20260728

Title: Change medium for hydroponic box 1, 2, 3, 4

Purpose: change medium to test on when we should add protectant

Procedure:

1. Take 4 new tip boxes, 3M NaCl, CH protectant, and 1000 μ L pipette, P100 pipette gun, and serological pipettes to BSC
2. Add 290 mL of medium to Box 1, 2, 3 and 4 and use the beaker with a label of 290 mL to measure the amount of medium
3. Add 600 μ L of CH protectant into Box 1
4. Add 10 mL of 3M NaCl to Box 1, 2, 3, and 4
5. Take the 3D printed rack from the corresponding boxes and put them into the corresponding boxes
6. Label the new boxes
7. Wrap your phone with plastic wrap and spray alcohol on it and use it to take test pictures of the 4 boxes

Done by: Abby T, Audrey H

2026/7/29 

Date: 20260729

Title: Gel extraction

Purpose: to isolate desired DNA fragment

Procedure:

1. Use a sanitized sharp knife to cut the gel band containing the desired DNA fragment under the UV light box.
2. Place each 3 gel bands into separated and labeled 2ml tubes
Mf002 bp: 89
Mf001 bp: 90
L3SIPS bp: 175
3. Add 700 μ L of GP Buffer into each of the tubes.
4. Place the tubes in a dry bath on 60°C to dissolve the gel slices, check every 2 minutes until entirely dissolved.
5. Load no more than 700 microliter of melted gel into the collection tube
6. Centrifuge for 30seconds at 7000 g. Discard the flow-through. Repeat step 5-6 step for the rest of the mixture.
7. Wash the column once with 0.5ml of WN Buffer by centrifuging for 30 seconds. Discard the flow through.
8. Wash the column once with 0.5ml of WS Buffer by centrifuging for 30 seconds. Discard the flow through.
9. Centrifuge the column at full speed 23444g for 3 minutes to remove residual ethanol.
10. Place the column onto a new centrifuge tube. Add 15 μ L of EB buffer onto the center of the membrane.
11. Stand the column for 2 minutes and centrifuge at full speed 23444g for 2 minutes to elute DNA.
12. Turn on BioDrop, adjust the setting to pathlength unit (concentration).
13. Clean the sample port with 1 μ L of Nucleus free water.
14. Calibrate.
15. Load 1 μ L of DNA from the 2 tubes each time and record down the results.
16. Clean the port with Nucleus Free water and close BioDrop.

Done by Abby T

Date: 20260729

Title: Making 1% Agarose gel

Purpose: Make Agarose ge; for electrophoresis

Procedure:

1. Add 0.2g of Agarose powder and 20mL of TAE buffer into a flask
2. Assemble mold, tray, and comb
3. Cover the flask with plastic wrap and poke few holes on it with tip
4. Microwave it for 30sec then agitate it, repeat this step until no sentiment
5. Wait until the flask doesn't smoke but is still hot, add 1 μ L of diluted safety dye into the flask and gently mix it
6. Pour the melted gel into the mold
7. Wait for it to solidify

Done by Ryan Yuan

Date:20260729

Title:Electroporation of K

Purpose: Transfer bio-engineered plasmid into *Bacillus subtilis* common cell

Procedure:

K: Pveg-EEBD-SamyQ-TAT-HA-B0015 354ng/μL

1. Prewarm SOC medium in 37 degree celsius
2. Precold 0.2cm electroporation cuvette
3. Mix 3μL plasmid (con 354ng/μL) with 100μL electrocompetent cell in a 1.5mL microfuge
4. Set up electroporator
 - A. Go to homepage
 - B. when using the 0.2cm corvette, press2 or 3 and enter
 - C. set voltage resistance, voltage and capacity
 - a. voltage:2.3KV
 - b. resistance:300Ω
 - c. capacitance:25F
5. Transfer mixture of cells+DNA to a cold electroporation corvette and tap the suspension to the bottom
6. Place the corvette in the shock pad and close the lid
7. Wait until beep sound
8. Remove the corvette quickly form the shock pad and immediately add 1mL SOC medium to the cuvette
9. Transfer the cell suspension to a polypropylene tube and incubate at 37°C for 3 hr, shaking at 200rpm
10. Check & record the parameter
11. Centrifuge the tubes at 5000 g for 5min
12. Spread the LB plate without erythromycin
 - TC(ms): 6.6
 - V0Hs(V):2284

Done by Gabriel

Date:20260729

Title: PCR

Purpose: Amplify DNA plasmid sequence number

Procedure:

1. Prepare the reaction as follows for MF001, MF002, L3S1P2

	1rxn	5rxn(master mix)
Ex taq buffer	5μL	25μL

2uM dNTP mix	4μL	20μL
10uM fw primer	1μL	5μL
10uM rv primer	1μL	5μL
DNA Temp	1μL	5μL
PCR grade water	37.75μL	188.75μL
Ex taq polymerase	0.25μL	1.25μL

2. Distribute 50μL of the 5rxn mixture into the corresponding tube of the part
3. Put into thermal cycle

98°C	30sec	1x cycle
98°C	10sec	35cycle
55°C	30sec	
72°C	40sec	
72°C	2min	
25°C	0sec	1cycle

Done by Sophie C, Abby T
Supervised by Gabriel

Date:20260729

Title: Make 2% gel(x3)

Purpose: Make gel for electrophoresis

Procedure:

1. Weight 1.2g of agarose powder and 60mL 1X TAE buffer
2. Add the components into a conical flask
3. Wrap a plastic wrap and poke 3 hole
4. Microwave for 30 sec and agitate. Repeat this step until all mixture dissolved
5. Take the flask out and wait for it to cool down
6. When there's no more smoke but still hot, add 3 μL of safety dye and agitate it
7. Pour the melted gel into the assembled mold and tray
8. Wait 20mins until the gel solidify in room temp

Done by Sophie C

Date:20260729

Title:Transformation

Purpose: Transfer bacterial liquid to agar plates

Procedure:

A:(STK) 1A:PtagC+RBS-Ho1:pcayA-L3S2P21

L:(STK)Exp1:3P-EEBD-RBS0-spacer-GanaA-Gala-spacer-B0015

1. Turn on the water bath at 42°C
2. Prepare an ice bath
3. Throw competent cells on ice until 50% thawed
4. Add 10 µL of plasmid
5. Incubate on ice for 20 min
6. Heat shock in the water bath for 1 min
7. Recover on ice for 3 min
8. Add 1 mL of SOC medium
9. Incubate at 37°C for 1 hr
10. Centrifuge at 5000 g for 5 min
11. Remove 1 mL of supernatant
12. Resuspend the pellet in remaining supernatant and spread on LBA plate
13. Incubate at 37°C for 16 hours

Done by SophieC

2026/7/30 

Date: 20260730

Title: Miniprep

Purpose: Collect the concentrated DNA plasmid

Procedure:

D7: pSTK-0A-Phyperspank

D10: pSTK-0AB-PtagC+RBS

D11: pSTK-0B-3P promoter

E8: pSTK0-B-RBS2

G12: pSTK-0D-L3S2P21

H2: pSTK-0D-bkdB Cold double terminator

1. Transfer 1000 μ l liquid culture into 1.5mL tube
2. Put in centrifuge and centrifuge for 30 seconds at 23444 g
3. Discard the supernatant
4. Repeat step 1-3 until all medium are collected
5. Add 200 μ l of Mx1 buffer
6. Resuspending by vortexing or pipetting
7. Add 250 μ l of Mx2 buffer
8. Invert slowly 8-10 times and incubate at RT for 2 mins
9. Add 350 μ l of Mx3 buffer and invert slowly 4-5 times
10. Centrifuge at 16000 g for 5 min
11. Transfer around 800 μ l supernatant into the column
12. Centrifuge at 7000 g for 30 min
13. Discard the flow through
14. Add 500 μ l of WN buffer, centrifuge 7000 g for 30 seconds, discard the flow through
15. Add 700 μ l of WS buffer, centrifuge 7000 g for 30 seconds, discard the flow through
16. Centrifuge 1000 g for 8min
17. Assemble 1.5mL microtube
18. Add 50 μ l of EB buffer on the center of the membrane
19. Incubate in RT for 3 min
20. Centrifuge 10000 g for 3min

Biodrop results:

	concentration(ng/ μ l)	A260/230	A260/280
D7	187.0	1.345	2.055
D10	69.0	1.255	1.974
D11	134.0	1.149	1.925
E8	108.0	1.091	1.929
G12	138.0	1.155	1.443

H2	110.0	0.846	1.914
----	-------	-------	-------

Done by Sarah and Sophie H

Date: 20260730

Title: Make B sub electrocompetent cells(x4)

Purpose: restock for 0730 electroporation

Procedure:

1. Dilute an overnight culture of B sub 16-fold in growth medium(LB medium containing 0.5M sorbitol) by adding 312.5 µl of B sub in 5mL GM and grow at 37°C to OD600 of 0.944
2. Cool the cells on ice and put on shaker for 10 minutes
3. Centrifuge at 4°C at 5000 g for 5 minutes to harvest cells
4. Add 1mL of ice cold electroporation medium to resuspend cells in the 1.5mL microtube by washing it along the side of the microtubes. Discard the medium. Repeat this step for 2 times to wash the cells
5. Pipette out the residue medium
6. Suspend the cells in 1/40 of the culture volume of the electroporation solution(132.8 µl for 3 tubes and 82.8 µl to the testing tube) with a cell concentration of $1-1.3 \times 10^{10}$ cfu/mL
7. Store the competent cells at -80°C

Done by Abigail and Ryan W

Supervised by Gabriel

Date: 20260730

Title: Colony PCR

Purpose: To verify bacterial colonies after a transformation step in molecular cloning

Procedure:

A: PtagC+RBS-Ho1:pcyA-L3S2P21

K: Pveg-SamyQ-TAT-HA-AmilCP(B.sub)

L:3P-GanA-Gala-B0015

1. Prepare the following reaction(master mix) and distribute solutions into PCR tubes

A: 13 reactions

Forward primer: Biobrick prefix	13 µl
Reverse primer: STK reverse(seq)	13 µl
Q-amp master mix	130 µl
ddH ₂ O	104 µl

K: 9 reactions

Forward primer: SamyQ F	9 μ l
Reverse primer: AmilCP R	9 μ l
Q-amp master mix	90 μ l
ddH ₂ O	72 μ l

L: 13 reactions

Forward primer: EXP fw	13 μ l
Reverse primer: EXP-seq4 rv	13 μ l
Q-amp master mix	130 μ l
ddH ₂ O	104 μ l

2. Circle 12 colonies each(from A and L) and 8 from K
3. Pick up those colonies with tip and dispense in PCR tubes
4. Take out tips and close the lid after incubating in RT for 5 min
5. Place into the PCR machine for the thermal cycle below:

95°C	5 min	x1
95°C	20 sec	x35
65°C	20 sec	
72°C	3 min	
72°C	5 min	x1

Done by Sophie C, Ryan Y, and Sophia

Date: 20260730

Title: Make 1% agarose gel

Purpose: Do electrophoresis for colony PCR

Procedure:

1. Assemble mold, tray, and comb
2. Weigh 0.6g of agarose powder and 60 mL of 1xTAE buffer
3. Add the components into a conical flask
4. Cover with plastic wrap and poke holes
5. Place into microwave for 30 secs
6. Agitate
7. Put into the microwave for 5 sec then agitate. Repeat until no sediments
8. When it stops smoking, add 3 μ L of safety dye, agitate

9. Pour the melted gel into the assembled mold
10. Let solidify at RT

Done by Chloe and Sophie C

Date: 20260730

Title: Make LB, LBA with glycerol, and LBE with glycerol 5µg/mL

Purpose: Restock plates

Procedure:

1. Measure and add the following into flasks:

LBA with glycerol

	Mass(g)
Broth powder	6.25g
Agar powder	3.75g
ddH ₂ O	250mL
Glycerol	2mL
Antibiotic(ampicillin)	500 µl

LBE with glycerol

	Mass(g)
Broth powder	6.25g
Agar powder	3.75g
ddH ₂ O	250mL
Glycerol	2mL
Antibiotic(Erythromycin)	1250 µl

LB

	Mass(g)
Broth powder	6.25g
Agar powder	3.75g
ddH ₂ O	250mL

2. Add one magnetic stirring bar to the flask
3. Place the flask containing the mixture on the hot plate and stir until no sediments
4. Cover the flask with aluminum foil and apply an indicator tape
5. Autoclave for 20 minutes at 127°C.
6. Place the new petri dishes in BSC and equipments
7. Set up water bath at 55°C
8. Spray the working surface of the BSC and equipments with 75% ethanol
9. Add 500 µl of Ampicillin/1250 µl of Erythromycin into flask
10. Use pipette to add 20 mL of melting LB into the petri dish
11. Let the plates dry on the BSC for 1 hour.

Done by Sarah and Sophie H

Date: 20260730

Title: Gel electrophoresis

Purpose: check bps of 0730 colony pcr

Procedure:

1. Load the samples in the following order:

Gel 1: M K1 K2 K3 K4 K5 K6 K7 K8
 5 5 5 5 5 5 5 5 5


M: [DM2300] ExcelBand™ 100 bp+3K DNA ladder

K: Pveg-SamyQ-TAT-HA-AmilCP

Forward primer: SamyQ F

Reverse primer: AmilCP R

bp: 842

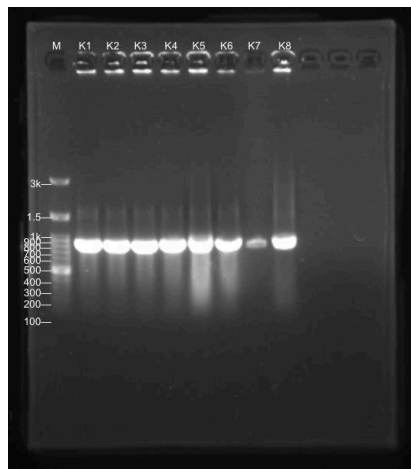

Gel 2: M A1 A2 A3 A4 A5 A6 A7 A8 A9 A10 A11 A12 L1 L2 L3 L4
 5 5 5 5 5 5 5 5 5 5 5 5 5 5 5 5
L5 L6 L7 L8 L9 L10 L11 L12
 5 5 5 5 5 5 5 5

M: 1kb DNA ladder


A: PtagC+RBS-Ho1:pcyA-L3S2P21

- Forward primer: Biobrick prefix
- Reverse primer: STK-seq rv
- bp: 1833

L: 3P-GanA-Gala-B0015

- Forward primer: EXP-seq fw
- Reverse primer: EXP-seq4 rv
- bp: 3626

A1, A2, A3 Liquid culture 7/30 -> sequencing L all fail

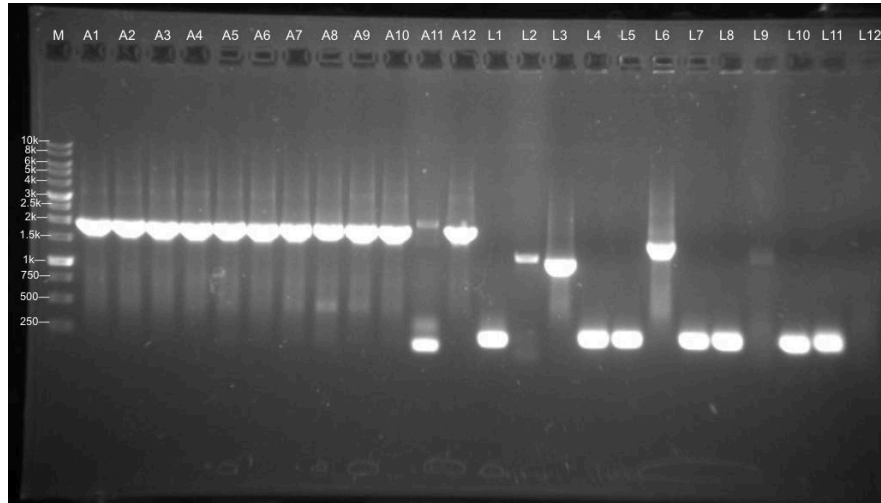

2. Resolve the gel at 100V for 30 minutes

Done by Chloe and Sophie C

Date: 20260730

Title: Transformation

Purpose: Transform liquid bacteria on plates

Procedure:

J: Pveg-csn:ACCD

1. Turn on the water bath and set it to 42°C and prepare ice
2. Thaw competent cells on ice until 50% thawed
3. Add 10 µl of plasmid
4. Incubate on ice for 20minutes
5. Heat shock in the water bath for 1 minute
6. Recover on ice for 3 min
7. Spread on LBA plate
8. Incubate at 37°C for 16 hours

Done by Sarah and Sophie H

Supervised by Chloe and Gabriel

Date: 20260730

Title: Electroporation

Purpose: Use electric pulse to make temporary holes in a cell membrane

Procedure:

J: Pveg-EEBD-RBS0-csn:ACCD-B0015

1. Prewarm soc medium at 37°C
2. Precold 0.2 cm electroporation cuvette
3. Mix 2.2 µl of plasmid(concentration 410 ng/ µl) with 86 µl of electrocompetent cell in a 1.5mL microtube
4. Set up electroporator and set up voltage resistance
 - a. Voltage: 23kv
 - b. Resistance: 300 Ω
 - c. Capacitance: 35F
5. Transfer mixture of cells and DNA to electroporation cuvette and tap suspension to bottom
6. Place in shock pad and close lid
7. Wait until beep
8. Remove from shock pad and immediately add 1 mL soc medium
9. Transfer to a polypropylene tube and incubate at 37°C for 3 hours. Shaking at 200 rpm
10. Check and record the parameter
11. Centrifuge tubes at 5000 g for 5 minutes
12. Spread the LB plate with erythromycin

TC(ms): 6.9

VOLts(V): 2285

Done by Chloe and Sophie C

Date: 20260730

Title: Golden Gate Assembly(G, I, M)

Purpose:

Procedure:

- 1.

L: PtagC+RBS-Spacer CX-GanA-Gala-Histag CZ-L3S2P21

Part	Dilution
Vector (shuttle vector): pSTK-EXP-1-sfGFP	1.39 µL (dilute 1/5)
Insert promoter(Pveg)	1.1 µL (dilute 1/10)
Insert RBS(EEBD-RBS0)	2 µL (Dilute 1/10)

Insert CDS(GanA-Gala)	1.53 (dilute 1/5) μ L
Insert Terminator (B0015)	1.94 μ L (dilute 1/10)
Spacer (pSTK-OCX-spacer)	1.07 μ L (dilute 1/10)
Histag (pSTK-OCZ-HisGX)	1.16 μ L (dilute 1/10)
T4 DNA ligase buffer	2 μ L
Sterile water	5.81 μ L
T4 DNA ligase	1 μ L
BsaI-HFv2	1 μ L
Total	20 μ L

G: (STK) Module 4-PcpcG2-veg - EEBD- RBS0 - CSN::LEA - B0015

Vector (1D) - concentration(ng/μL): 132 - bp: 3120	1.5 μ L (dilute 1/10)
Insert promoter(PcpcG2-veg) - concentration(ng/μL): 168 - bp: 235	1.175 (dilute 1/50)
Insert RBS(EEBD-RBS0) - concentration(ng/μL): 150 - bp: 2392	1.04 μ L (Dilute 1/5)
Insert CDS(CSN::LEA) 77 ng/ μ L 634 bp	1.08 μ L (dilute 1/20)
Insert Terminator (B0015) - concentration(ng/μL): 107 - bp: 2213	1.34 μ L (dilute 1/5)

T4 DNA ligase buffer	2 μ L
Sterile water	10.05 μ L
T4 DNA ligase	1 μ L
Bsal-HFv2	1 μ L
Total	20 μ L

M: (STK) Module 4-PcpcG2-veg - EEBD- RBS0 - spacer CX - GanAGala - histag CZ - dkdB cold

Vector (1D) - concentration(ng/μL): 132 - bp: 3120	1.5 μ L (dilute 1/10)
Insert promoter(PcpcG2-veg) - concentration(ng/μL): 65 - bp: 235	1.175 (dilute 1/50)
Insert RBS(EEBD-RBS0) - concentration(ng/μL): 150 - bp: 2392	1.04 μ L (Dilute 1/5)
Insert CDS(CSN::LEA) 77 ng/ μ L 634 bp	1.08 μ L (dilute 1/20)
Insert Terminator (B0015) - concentration(ng/μL): 107 - bp: 2213	1.34 μ L (dilute 1/5)
T4 DNA ligase buffer	2 μ L

Sterile water	10.05 µL
T4 DNA ligase	1 µL
Bsal-HFv2	1 µL
Total	20 µL

Date: 20260730

Title: Liquid culture

Procedure:

A:Ptag+RBS-Ho1:PcyA-L3S2P21

K:pExp1-Pveg-SamyQ-TAT-HA-AmilCP

1. Label 7 liquid culture tube(A1,A2,A3,K1,K2,K3,K8)
2. Turn on Bunsen Burner
3. Add 3mL of TB medium and 3µL of Ampicillin
4. Pip a pipette tip into the colonies
5. Dispense tips in their respective tubes
6. Place the tubes in the incubator, slant around 45°C

Done by Ryan Wei

Supervised by Gabriel

Date: 20260730

Title: Plant seeds for set C and vernalize in hydroponic system

Purpose: heat stress testing

Procedure:

1. Create lids with 3 mm diameter holes, and autoclave for 20 min at 127°C
2. Eppendorf
3. Add 0.8% agar into 48 lids (270 µL each)
4. Wait until agar solidifies
5. Add 2.2 mL of ½ MS medium into 48 tube racks
6. Put lids onto the tube rack (hole facing up)
7. Use an autoclaved toothpick to pick one seed, and add the seed into the hole of the lid
8. Repeat step 6 until all lids have seeds
9. Close the lid of the tube rack
10. Wrap the tube rack with aluminium foil
11. Place the tube rack into the refrigerator at 4°C

Done by: Alex

Date: 20260730

Title: Prepare 50 mL 0.8% ½ MS agar

Purpose: add agar into lids for new set of seeds

Procedure:

The following protocol is for 50 mL ½ MS agar (0.8%)

1. ½ MS powder (w/o sucrose): 0.115 g
2. MES buffer: 0.025 g
3. Distilled water: 50 mL
4. Sucrose: 0.5 g
5. Agar: 0.5 g 0.4 g
6. Add ½ MS powder, MES buffer, sucrose, and distilled water into a flask
7. Shake the flask with hands until no sediments are in the solution
8. Use the pH probe to measure the pH of the solution, ensure the pH is between 5.7–5.8
9. If pH is too low, add 0.1 M KOH slowly to increase it
10. Add agar into the flask
11. Wrap the beaker with aluminium foil and stick a piece of indicator tape on top
12. Autoclave for 20 min at 127°C

Done by: Alex

Date: 20260730

Title: Prepare 2 1500 mL ½ MS medium

Purpose: prepare medium for changing medium and salinity test's medium

Procedure:

The following protocol is for 1500 mL ½ MS medium

1. ½ MS powder (w/o sucrose): 3.45 g
 - a. MES buffer: 0.75 g
 - b. Distilled water: 1500 mL
2. Add the contents into a beaker
3. Use hot plate to stir (using a magnetic stirring bar) until no sediments are left in the solution
4. Use the pH probe and measure the pH of the solution, ensure the pH is between 5.7 and 5.8
5. If pH is too low, add 1 M KOH to slowly increase it
6. Wrap the beaker with aluminium foil and stick a piece of autoclave indicator tape on top
7. Repeat steps 1–6 to prepare the second 1500 mL ½ MS medium
8. Autoclave for 20 min at 127°C

Done by: Alex

2026/7/31 

Date: 20260731

Title: Colony PCR of J

Purpose: To verify bacterial colonies after a transformation step in molecular cloning

Procedure:

J:Pveg-EEBD-RBS0-csn:AccD-B0015(B.sub)

1. Prepare the following reaction (master mix) and distribute solutions to PCR tubes
JX9

Forward primer Csn fw	9µL
Reverse primer AccD R	9µL
Q-amp master mix	90µL
ddH ₂ O	72µL

2. Circle 8 colonies from plate J
3. Pick up those colonies with tips and dispense in PCR tubes
4. Take out the tips and close the lid after incubating in RT for 5 min
5. Place into PCR machines

95°C	5 min	x1
95°C	20sec	x35
55°C	20sec	
72°C	2min	
72°C	5min	x1

Done by Sophie C, Ryan Y

Date: 20260731

Title: Make 1% agarose gel

Purpose: for electrophoresis of colony PCR

Procedure:

1. Assemble tray, madd, comb
2. Weigh 0.2g of agarose powder and measure 20mL 1xTAE buffer
3. Add all into a conical flask
4. Cover with plastic wrap&poke holes
5. Microwave for 30sec agitate
6. Repeat until no sediments
7. When no smoke , add 1µL safety dye agitate
8. Pour melted gel into mold
9. Let it solidify in RT

Done by Ryan Y, Sophie C

Date: 20260731

Title: Gel electrophoresis

Purpose: Check base pair of colony PCR product


Procedure:

J:Pveg-EEBD-RBS0-csn:AccD-B0015

1. Load the samples in the following order:

<u>M</u>	<u>J1</u>	<u>J2</u>	<u>J3</u>	<u>J4</u>	<u>J5</u>	<u>J6</u>	<u>J7</u>	<u>J8</u>
5	5	5	5	5	5	5	5	5
2. Resolve the gel at 100vv for 30min

20260731 Colony PCR (J in B. Sub)

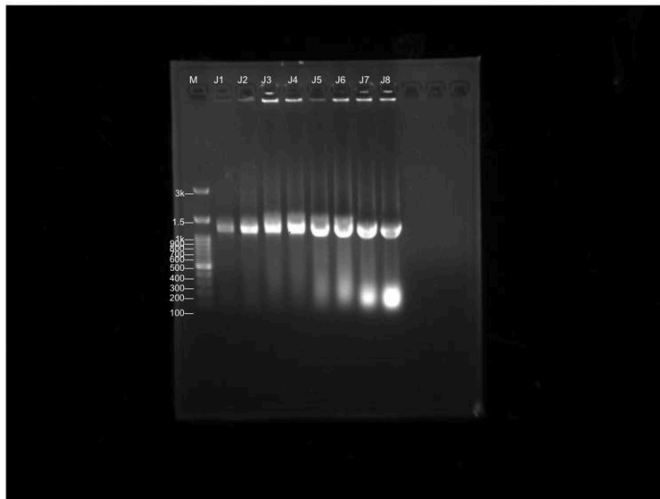

M: DNA ladder

J: Pveg-EEBD-RBS0-Csn:ACCD-B0015

- Forward primer: Csn F
- Reverse primer: ACCD R

bp: 1145(?)

Date: 20260731

Title: Transformation

Purpose: Transform and cultivate engineered bacteria on agar plate

Procedure:

G:PcpcG2-veg-EEBD

1. Prepare ice bucket and set the water bath to 42°C
2. Get LBA plates
3. Thaw competent cells until 50% thawed and add 10µL plasmid
4. Incubate on ice for 20 min
5. Heat shock at 42°C for 45 sec-1min
6. Recover on ice for 3min
7. Add 1mL SOC medium
8. Incubate at 37°C for 1 hr
9. Centrifuge at 5000g for 5min
10. Take out 1000µL supernatant&resuspend pellet
11. Spread on LBA plate and incubate at 37°C overnight

Done by Chloe